

This Page Is Inserted by IFW Operations
and is not a part of the Official Record

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images may include (but are not limited to):

- BLACK BORDERS
- TEXT CUT OFF AT TOP, BOTTOM OR SIDES
- FADED TEXT
- ILLEGIBLE TEXT
- SKEWED/SLANTED IMAGES
- COLORED PHOTOS
- BLACK OR VERY BLACK AND WHITE DARK PHOTOS
- GRAY SCALE DOCUMENTS

IMAGES ARE BEST AVAILABLE COPY.

**As rescanning documents *will not* correct images,
please do not report the images to the
Image Problem Mailbox.**

THIS PAGE BLANK (USPTO)




INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(51) International Patent Classification ⁷ : G06F 19/00		A2	(11) International Publication Number: WO 00/28460
			(43) International Publication Date: 18 May 2000 (18.05.00)
(21) International Application Number: PCT/US99/26521		(74) Agents: MEHRMAN, Michael, J. et al.; Jones & Askew, LLP, 2400 Monarch Tower, 3424 Peachtree Road, Atlanta, GA 30326 (US).	
(22) International Filing Date: 8 November 1999 (08.11.99)			
(30) Priority Data: 60/107,704 9 November 1998 (09.11.98) US 60/144,705 20 July 1999 (20.07.99) US		(81) Designated States: AE, AL, AM, AT, AU, AZ, BA, BB, BG, BR, BY, CA, CH, CN, CR, CU, CZ, DE, DK, DM, EE, ES, FI, GB, GD, GE, GH, GM, HR, HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MA, MD, MG, MK, MN, MW, MX, NO, NZ, PL, PT, RO, RU, SD, SE, SG, SI, SK, SL, TJ, TM, TR, TT, TZ, UA, UG, US, UZ, VN, YU, ZA, ZW, ARIPO patent (GH, GM, KE, LS, MW, SD, SL, SZ, TZ, UG, ZW), Eurasian patent (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM), European patent (AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE), OAPI patent (BF, BJ, CF, CG, CI, CM, GA, GN, GW, ML, MR, NE, SN, TD, TG).	
(71) Applicant (for all designated States except US): LIFESTREAM TECHNOLOGIES, INC. [US/US]; 510 Clearwater Loop, Suite 101, Post Falls, ID 83854 (US).			
(72) Inventors; and (75) Inventors/Applicants (for US only): MAUS, Christopher, T. [US/US]; 3277 Jewel Lake, Sagle, ID 83860 (US). CONNOLLY, Jackson, B. [US/US]; 4677 East Alpine Drive, Post Falls, ID 83854 (US). COAD, Craig, A. [US/US]; 20600 Canary Road, Cataldo, ID 83810 (US). COAD, Noah, M. [US/US]; 20600 Canary Road, Cataldo, ID 83810 (US). MOODY, James, L. [US/US]; 32501 SW Juliette Drive, Wilsonville, OR 97070 (US). NESBITT, Kenn, A. [US/US]; 2719 West Summit Boulevard, Spokane, WA 99201 (US). CLEGG, Kenneth, D. [US/US]; 14534 North Fairview Drive, Mead, WA 99021 (US).		Published Without international search report and to be republished upon receipt of that report.	

(54) Title: HEALTH MONITORING AND DIAGNOSTIC DEVICE AND NETWORK-BASED HEALTH ASSESSMENT AND MEDICAL RECORDS MAINTENANCE SYSTEM

(57) Abstract

A health monitoring and diagnostic device (LIFESTREAM cholesterol meter) configured as a self-contained testing and diagnostic unit in a clam-shell type case. One side of the case includes a spring-loaded finger stick and a compartment for carrying one or more packages of disposable items including a test strip, a needle for the finger stick, and an alcohol swipe. The other half of the case includes a test strip reader, a key pad, and a liquid crystal display. The meter reads a test strip carrying a droplet of blood and receives additional diagnostic information from the patient, such as age, gender, weight, and family history of heart disease. Within minutes, the

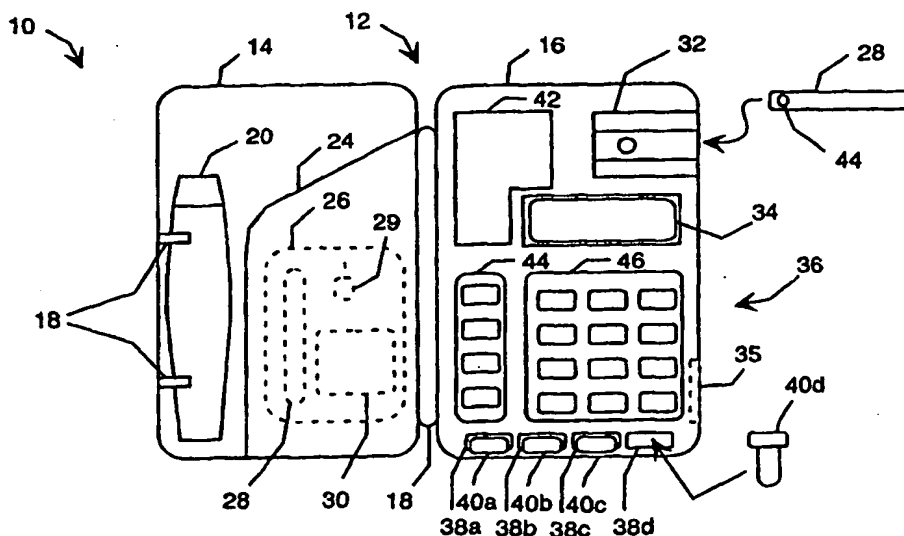

meter displays test results, including total cholesterol levels. The meter also displays additional diagnostic results, such as the patient's "cardiac age", recommended weight loss, and a cardiac risk assessment. The meter also works in connection with a network-based comprehensive health analysis and reporting system. The meter writes patient data to a smartcard. This patient data typically includes patient identification information, the test results, the diagnostic information, and the diagnostic results. A computer station reads the smartcard and establishes a network connection with a health report server over the Internet. The computer then downloads the patient data to the health report server, which prepares a comprehensive health report. Within minutes, this report is transmitted back to the computer station, where it is printed out and delivered to the patient.

FOR THE PURPOSES OF INFORMATION ONLY

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AL	Albania	ES	Spain	LS	Lesotho	SI	Slovenia
AM	Armenia	FI	Finland	LT	Lithuania	SK	Slovakia
AT	Austria	FR	France	LU	Luxembourg	SN	Senegal
AU	Australia	GA	Gabon	LV	Latvia	SZ	Swaziland
AZ	Azerbaijan	GB	United Kingdom	MC	Monaco	TD	Chad
BA	Bosnia and Herzegovina	GE	Georgia	MD	Republic of Moldova	TG	Togo
BB	Barbados	GH	Ghana	MG	Madagascar	TJ	Tajikistan
BE	Belgium	GN	Guinea	MK	The former Yugoslav Republic of Macedonia	TM	Turkmenistan
BF	Burkina Faso	GR	Greece	ML	Mali	TR	Turkey
BG	Bulgaria	HU	Hungary	MN	Mongolia	TT	Trinidad and Tobago
BJ	Benin	IE	Ireland	MR	Mauritania	UA	Ukraine
BR	Brazil	IL	Israel	MW	Malawi	UG	Uganda
BY	Belarus	IS	Iceland	MX	Mexico	US	United States of America
CA	Canada	IT	Italy	NE	Niger	UZ	Uzbekistan
CF	Central African Republic	JP	Japan	NL	Netherlands	VN	Viet Nam
CG	Congo	KE	Kenya	NO	Norway	YU	Yugoslavia
CH	Switzerland	KG	Kyrgyzstan	NZ	New Zealand	ZW	Zimbabwe
CI	Côte d'Ivoire	KP	Democratic People's Republic of Korea	PL	Poland		
CM	Cameroon	KR	Republic of Korea	PT	Portugal		
CN	China	KZ	Kazakhstan	RO	Romania		
CU	Cuba	LC	Saint Lucia	RU	Russian Federation		
CZ	Czech Republic	LI	Liechtenstein	SD	Sudan		
DE	Germany	LK	Sri Lanka	SE	Sweden		
DK	Denmark	LR	Liberia	SG	Singapore		
EE	Estonia						

5

10 **HEALTH MONITORING AND DIAGNOSTIC DEVICE AND**
 NETWORK-BASED HEALTH ASSESSMENT AND MEDICAL
 RECORDS MAINTENANCE SYSTEM

REFERENCE TO RELATED APPLICATION

15 This application claims priority to commonly-owned U.S. Provisional Patent Application No. 60/107,704, filed November 9, 1998; and commonly-owned U.S. Provisional Patent Application No. 60/144,705, filed July 20, 1999.

TECHNICAL FIELD

20 This invention relates to health monitoring and diagnostic devices and, more particularly, relates to a hand-held device operable for determining blood lipid levels from test-strip analyses, obtaining additional diagnostic information from a user, displaying corresponding diagnostic results, and storing this data on a secure patient-held data carrier, such as a smartcard. The invention also relates to a
25 secure network-based health assessment and medical records maintenance system that receives medical information from the health monitoring and diagnostic devices, produces health assessments based on the received medical information, and stores the received medical information in a secure medical records maintenance system.

30

BACKGROUND OF THE INVENTION

American health care is undergoing a revolution. By the year 2000, more than two-thirds of all American workers with health insurance will be enrolled in

some kind of managed care plan, where the emphasis is on early detection of disease and preventive care.

Fueling this revolution is the skyrocketing cost of health care, combined with new medical research showing lifestyle is important to good health. In fact, in its 1982 report on "Health and Behavior," the National Academy of Sciences concluded that half of the ten leading causes of death in the United States are primarily related to lifestyle. Dietary patterns are identified as key lifestyle choices.

Cholesterol levels are particularly important in the United States. For this reason, the American Heart Association, the American Medical Association and the health related agencies of the U.S. government have embarked on national education campaigns to inform the public about the importance of making lifestyle changes to lower blood cholesterol and prevent heart disease. Although the dangers of high cholesterol have been widely publicized, many people fail to make effective use of this information because they do not know their own blood cholesterol levels. In other words, a great many of the people with high cholesterol levels fail to heed the advice to lower their cholesterol levels simply because they are unaware of their own cholesterol levels.

This situation persists because of the high cost and inconvenience presently involved in obtaining cholesterol information. To obtain this information, most people go to a physician's office, have blood drawn, and wait for the return of the blood chemistry analysis. Often, obtaining the results involves a second trip to the physician's office. This is expensive and time consuming; the average cost is about \$83 for each office cholesterol consultation, and the average wait for the results is several days.

The cost and inconvenience involved in obtaining cholesterol tests inhibits many people from testing their cholesterol frequently enough to provide effective positive feedback. As a result, many people who begin corrective exercise, diet, or drug therapy programs in response to high cholesterol tests often give up their corrective programs because they do not monitor their cholesterol frequently enough to remain aware of the benefits of their programs. Moreover, blood cholesterol numbers by themselves are often poor motivators for patients who feel and look fine, and do not immediately feel or look differently when they take their prescriptions. In fact, studies have shown that 80% of the patients prescribed

cholesterol-lowering drug therapies stop taking their prescriptions within a few months. And the attrition rates for exercise and diet programs may be even higher.

In addition, there is a need for a medical records maintenance system, not only for blood cholesterol tests but for many types of medical information that can be obtained outside of the hospital environment. This need will increase with increases in the availability of remote health monitoring devices in the future, such as blood pressure measuring devices, blood sugar testing devices, blood cholesterol testing devices, AIDS testing devices, heart monitoring devices, sleep respiration monitoring devices, reproductive cycle and pregnancy monitoring devices, epileptic and other types of seizure monitoring devices, and a wide range of other remote health monitoring devices that may be developed in the future. As the availability of the remote health monitoring devices increases, users will have an increasing need for securely storing the tests results in electronic format. The current system of hard-copy and electronic medical records maintained in doctors' offices will become increasingly obsolete and inconvenient as the availability of electronically-stored medical data increases. Because a patient's medical records are highly confidential, there is a need for a highly secure and permanent medical records maintenance system under the control of individual patients and their doctors.

Thus, there is a general need in the art for a less expensive and more convenient approach to providing cholesterol tests. There is a further need for making motivational information regarding cholesterol levels more readily available and more effective. And there is yet another need for a highly secure and permanent medical records maintenance system under the control of individual patients and their doctors.

SUMMARY OF THE INVENTION

The present invention meets the needs described above in a health monitoring and diagnostic device referred to as a LIFESTREAM cholesterol meter. This meter is configured as a self-contained testing and diagnostic unit in a clam-shell type case. One side of the case includes a biological sample gathering device, such as spring-loaded finger stick, and a compartment for carrying one or

more packages of disposable items, typically including a test strip, a needle for the finger stick, and an alcohol swipe. The other half of the case includes a test strip reader, a user input device such as a key pad, and a display device such as a liquid crystal display. The meter reads a test strip carrying a biological sample, such as a droplet of blood, and within minutes displays test results, such as total cholesterol levels, on the meter's display.

The hand-held LIFESTREAM cholesterol meter drastically reduces the costs and inconvenience associated with obtaining cholesterol tests by performing total cholesterol tests in virtually any location, including a physician's office, a pharmacy, a clinic, or in the privacy of the patient's home. The meter produces the test results within minutes using on-board circuitry and programming. The meter also includes an on-board diagnostic program that prompts for additional diagnostic information, such as the patient's age, gender, weight, family history of heart disease, blood pressure, and so forth.

The meter then translates this diagnostic information, along with the test results, into diagnostic results that may be more meaningful to the user than the test results alone. For example, the meter may use a well-known methodology, such as the Framingham Medical Study, to produce diagnostic results including the user's cardiac age (as compared to chronological age), recommended weight loss, 5-year risk of heart attack, 10-year risk of heart attack, an assessment of stroke risk, and other results that will be easily and immediately understood by the patient. Like the test results themselves, these more meaningful diagnostic results are displayed on the meter within minutes.

Producing diagnostic results like "cardiac age" and "5-year risk of heart attack" rather than total cholesterol levels alone may motivate more people to change their lifestyles and reduce their cholesterol levels. Moreover, producing these diagnostic results instantaneously, inexpensively, and in a convenient location encourages frequent testing and provides patients with the positive feedback necessary to encourage continued compliance with drug therapies and lifestyle changes. Ultimately, widespread use of the LIFESTREAM cholesterol meter can be expected to improve cardiac health nationwide, shift the focus of cardiac treatment from corrective to preventative, improve the cardiac health of the population in general, and reduce medical costs and health insurance rates.

The benefits of the LIFESTREAM cholesterol meter may be improved over time and extended to other health problems because the meter is programmable and configured to perform multiple types of tests. That is, although the meter will be initially configured to perform total cholesterol tests using test strips and human blood samples, it is also configured to perform multiple types of tests using different types of test strips or other test media carrying other types of biological fluid or tissue samples. For example, the meter may also produce other types of blood lipid test results, such as HDL cholesterol, triglycerides, LDL cholesterol, etc. The meter may also perform other types of tests, such as blood glucose tests, AIDS tests, cancer tests, and virtually any other type of test that can be performed using a test strip or another suitable test medium carrying a sample of biological fluid or tissue. To accommodate multiple tests, the meter typically includes four romkey sockets that allow the meter to carry and read four different romkeys.

The LIFESTREAM cholesterol meter also works in connection with a network-based comprehensive health analysis and reporting system. The meter includes a data drive that writes patient data stored within the meter to a patient-held data storage device, such as a smartcard. This patient data typically includes patient identification information, the test results, the diagnostic information, and the diagnostic results. A computer station, such as a typical desktop or laptop personal computer, can then read the smartcard and establish a network connection with a health report server, typically over the Internet. The computer then downloads the patient data to the health report server, which prepares a comprehensive health report. This report is then transmitted back to the computer station, where it is printed out and delivered to the patient.

The health report server typically works in concert with the patient's physician or pharmacist, who may provide additional diagnostic information to the server, such as a newly-prescribed drug therapy, other currently-prescribed drugs for the patient, exercise and dietary recommendations, and so forth. Within minutes, the health report server assembles a comprehensive health report including a data sheet for the newly-prescribed drug, cross-reaction information for the newly-prescribed drug and the other currently-prescribed drugs, weight and total cholesterol goals, exercise and dietary recommendations, any food or

activity warnings associated with the overall therapy package, and recommendations for on-going monitoring using the meter. This provides a complete written record of the patient's current condition, the therapy prescribed by the physician and filled by the pharmacist, and a roadmap for monitoring the patient's progress during the ensuing therapy.

The comprehensive health report may also include additional patient-specific information, such as the diagnostic information and results compiled by the meter, and additional diagnostic and health assessment information compiled by the server. For example, the report may include a trend analysis showing how cholesterol, blood glucose, and weight levels have changed over multiple readings. The report may also include generally-applicable educational information, such as coronary risk factors, dietary guidelines for reducing cholesterol levels, diabetes information, cancer information, and the like. At present, a patient may have to undergo a physical examination, pay thousands of dollars, and wait weeks to obtain a similar comprehensive health report. The network-based comprehensive health analysis and reporting system, working in concert with the LIFESTREAM cholesterol meter, allows the patient to obtain the report within minutes at a fraction of the cost.

The meter also includes a number of advantageous security features. For example, the meter cannot be activated until a user enters a proper activation code. This typically requires that the user call the manufacturer, which provides an opportunity to verify the meter's authenticity, set up a data file for the meter in the health report server, and tell the user how to update the meter software, if necessary. If a software update is indicated, the user may be instructed to activate the meter, initialize a smartcard, load the smartcard into a computer station, and establish a network connection with the health report server. The server can then download the new software (e.g., new version of an existing software module or a new software module) to the smartcard, which, in turn, can be placed back in the meter. The new software can then be uploaded to the meter.

The meter may also require validation of all test strips. Validation is important for some types of tests because readings obtained from each test strip will have to be interpreted correctly to obtain correct test results, and the calibration data used to interpret the readings from different lots of test strips may

vary significantly. To allow proper calibration, each lot of test strips has a corresponding memory device, such as a romkey, that must be placed into the meter. The romkey includes a code number, an expiration date, and the calibration data for interpreting readings from the corresponding test strips. A test strip identification number that is mathematically derived from the code number is printed on the test strips or their packaging. The user must enter the proper test strip identification number into the meter, which the meter verifies with reference to the code number and the expiration date read from the romkey. This allows the meter to prevent the use of expired test strips and to also prevent test strips from being used in combination with incorrect romkeys.

Test strip validation is also an important aspect of one business model for deploying the meters. That is, the meters themselves may be provided for use at little or no charge to individual patients, whereas proprietary test strips will be sold to generate revenue from use of the meter. This may be a desirable business model for deploying the devices because it minimizes the initial cost that an individual patient must pay to begin using the device. Having to sell each device at its full cost, on the other hand, would undermine the economic feasibility of using the device in many contexts. For this business model, the meter should only activate for use with proprietary test strips after validation of the test strips.

The meter may also require each smartcard to be initialized with a personal identification number (PIN). Patient-specific PINs allow multiple patients to use the same meter, and also allows each patient's data to be secure to that patient. That is, only the patient or someone authorized by the patient (i.e., knowing the patient's PIN) can read the medical data stored on the smartcard. In this manner, each patient controls his or her own medical data, which can be a particularly important attribute for highly sensitive medical data, such as AIDS tests, cancer tests, and the like.

Generally described, the invention provides a test strip for use with a health monitoring device or meter. The test strip, when carrying a sample of biological fluid or tissue, may be read by the meter to obtain test results based on the sample and calibration data specific to the test strip. The test strip also corresponds to a memory device that stores a code number and the calibration data, which may also be read by the meter. The test strip has an associated test strip identification

number that is mathematically derived from the code number and printed on the test strips, the packaging for the test strips, or a tag packaged with the test strips.

To verify test strips, the meter reads the code number from the memory device, mathematically derives a test strip identification number corresponding to
5 the code number, compares the received test strip identification number to the derived test strip identification number, and activates the meter for use with the test strip only if the received test strip identification number corresponds to the derived test strip identification number.

The memory device may also store an expiration date for the test strip,
10 which may be read by the meter. In this case, the meter may activate for use with the test strip only if the expiration date is prior to a current date read by the meter from an internal clock. The memory device may be a romkey that is inserted into a socket housed within the meter. The romkey is typically packaged with an associated group of the test strips, and the test strip identification number is
15 typically printed on the test strips, printed on packaging for the test strips, or printed on a tag packaged with the test strip.

The invention also provides a hand-held health monitoring device or meter that includes an enclosure for housing a disposable test strip for use with the meter. The meter also includes a holder for removably supporting a device for
20 gathering a sample of biological fluid or tissue, such as a finger stick. The meter also includes a test strip reader operable for reading the test strip carrying the sample of biological fluid or tissue and obtaining test results based on the sample and calibration data specific to the test strip. A memory reading device (e.g., romkey socket) functionally connected to the test strip reader reads the calibration
25 data from a memory device (e.g., romkey). A user input device, such as a key pad, receives user input commands and a display device, such as a liquid crystal display, displays information on the meter.

The meter also includes a processor that is functionally connected to the test strip reader, the user input device, and the display device. The processor contains
30 a program module that obtains the test results from the test strip reader and causes the display device to display the test results. A data drive functionally connected to the processor writes the test results to a removable memory storage device, such as a smartcard. The meter may be packaged in a clam-shell case that

opens to reveal first and second compartments. The first compartment may contain the enclosure for housing the disposable test strip and the holder for removably supporting the biological fluid or tissue gathering device, and the second compartment may contain the test strip reader, the memory reading device, the display device, the processor, and the data drive.

To provide activation verification, the meter may receive an activation code through the user input device, compute an activation code based on the current date and instructions contained in an activation routine stored within the meter, and activate the meter only if the computed activation code corresponds to the received activation code. In addition, to provide security to a patient's medical data, the meter may determine whether a PIN has been previously stored on the removable memory storage device. If a PIN has not been previously stored on the removable memory storage device, the meter prompts the user to enter a PIN and stores the received PIN on the removable memory storage device. Alternatively, if a PIN has been previously stored on the removable memory storage device, the meter prompts the user to enter a PIN, compares the stored PIN to the received PIN, and writes the test results to the removable memory storage device only if the stored PIN corresponds to the received PIN.

The test strip reader may also be operable for reading a second type of test strip carrying a second sample of biological fluid or tissue and obtaining health-related test results based on the second sample of biological tissue or fluid and calibration data specific to the second type of test strip. In this case, the meter may include a second memory reading device (e.g., romkey socket) functionally connected to the test strip reader and operable for reading calibration data from a second memory device (e.g., romkey) corresponding to the second type of test strip. For example, the meter may read both blood lipid test strips and blood glucose test strips. As noted previously, the meter typically includes four romkey sockets that allow the meter to carry and read four different romkeys.

The meter may also prompt the user to enter diagnostic information using the user input device, such as gender, ethnicity, family history of heart disease, personal history of heart disease, personal history of diabetes, personal history of smoking, height, weight, age, blood pressure, and fitness level. The meter may then perform a diagnostic analysis and produce diagnostic results based on the test

results and diagnostic information, and display diagnostic results. For example, the diagnostic results may include a medical risk index, a recommended weight loss, a five-year risk of heart attack, a ten-year risk of heart attack, a cardiac age, an extended age, and a risk of stroke.

5 The invention also provides a system for remotely producing health reports. This system includes a health monitoring device or meter, as described above, a computer station, and a health report server connected with the computer station through a network, such as the Internet. The meter writes health-related test results to a memory storage device. The computer station reads the test results
10 from the memory storage device, establishes a network connection with the health report server, receives additional diagnostic information from a user, and transmits the test results and the additional diagnostic information to the health report server. The server, in turn, compiles a health report based on the test results and the additional diagnostic information and transmits the health report to the
15 computer station, where the report may be printed and delivered to the patient.

 The health report may include a trend analysis with test results compiled for a number of samples, such as total cholesterol level and blood glucose level trend reports. The additional diagnostic information may include a newly-prescribed drug and other currently-prescribed drugs, and the health report may include a
20 data sheet for the newly-prescribed drug and information relating to cross-reactions between the newly-prescribed drug and the other currently-prescribed drugs. The health report may also include a target weight and total cholesterol levels, a schedule for future testing using the meter, health assessment summary, a coronary risk assessment, dietary guidelines to lower cholesterol, and other
25 educational information.

 The business model described above is largely dependent on the sale of proprietary test strips for the collection of revenue from end users. That is, the health monitoring device itself may be made available to individual patients at little or no cost, with the sale of proprietary test strips providing a major source of
30 revenue for the proprietor of the health monitoring device. As noted previously, this may be a desirable business model for deploying the devices because it minimizes the initial cost that an individual patient must pay to begin using the

device. Having to sell each device at its full cost, on the other hand, would undermine the economic feasibility of using the device in many contexts.

Nevertheless, it may also be desirable to provide a health monitoring device that does not rely on the sale of proprietary test strips as a major source of revenue. For example, the health monitoring device may be adapted to read non-proprietary test strips, or may incorporate a reusable and/or non-invasive testing device, such as an electrode, blood pressure monitoring device, sonic testing device, thermometer, saliva testing device, optical testing device, and the like. Of course, a non-invasive multi-use testing device may be used many times without affording the proprietor of the health monitoring device an opportunity collect revenue associated with each use of the device.

To provide an opportunity for the proprietor of the health monitoring device to collect revenue based on use of the device, the removable memory storage device may be utilized as a type of "debit card" or payment source for use with the health monitoring device. That is, the removable memory storage device may be purchased with a monetary value, or it may have a monetary value that is replenishable over the Internet using a bank credit or debit card or other conventional payment source. The health monitoring device may then deduct the cost of performing particular services from the monetary value represented by the monetary balance stored on the removable memory storage device. In other words, the health monitoring device may be configured to activate for the performance of a service upon deducting a charge for the service from a monetary value stored on a removable memory storage device inserted into the device.

This business model includes a health monitoring device operable for obtaining medical data associated with a patient and reading an initial monetary balance stored on a removable memory storage device. The health monitoring device determines whether the initial monetary balance is sufficient to pay a monetary value assigned to performance of a test involving the medical data to be performed by the testing device. If the initial monetary balance is sufficient to pay for the test, the health monitoring device computes a revised monetary balance by deducting the monetary value assigned to performance of the test from the initial monetary balance, replaces the initial monetary balance with the revised monetary

balance on the removable memory storage device, and activates the health monitoring device for performance of the specified service.

5 The business model also includes a system that includes one or more of the health monitoring devices described above, one or more removable memory storage devices, and a network-based server operable for remotely charging a cost to a payment source and crediting the cost to an initial balance stored on the removable memory storage device. The network-based server may also remotely store the monetary value assigned to performance of the test on the removable memory storage device. In this case, the health monitoring device reads the
10 monetary value assigned to performance of the test from the removable memory storage device. Thus, rate schedules for various services to be performed by the health monitoring device may be changed from time to time, based on quantity discounts or other considerations.

The invention also includes a secure medical records maintenance system.
15 Although this system is specifically adapted for use with the health monitoring device described above, it may be used to store any type of electronic data including a wide variety of medical records, and is particularly convenient for storing a wide range of electronic medical data generated remotely from the hospital or doctor's office environment. The secure medical records maintenance
20 system includes a number of removable memory storage devices, which are each operable for storing medical data for an associated patient. Each removable memory storage device also stores a patient-specified personal identification number (PIN), a medical records identification number secured by the PIN, and a patient identification number secured by the PIN.

25 The data stored on the removable memory storage device is downloadable to a two-server system including a first remote server that stores patient identification information indexed by patient identification numbers, and a second remote server that stores patient medical data indexed by the medical records identification number. For security purposes, the medical data maintained in the
30 second remote server cannot be correlated to the associated patient identification information maintained in the first remote server based on the information contained in the first and second remote servers.

To allow correlation of the data stored in the two servers, the secure medical records maintenance system includes a correlation table uniquely associating each medical records identification number with a particular one of the patient identification numbers. The correlation table for a particular patient typically resides on the patient's removable memory storage device. The correlation table for a practitioner's patients may also reside on the practitioner's computer, such as a doctor's or pharmacist's computer, that is associated with a licensed medical practitioner having an assigned professional registration number. For further security, the first and second remote servers are accessed by the practitioner's computer through encrypted communications secured by an application procedure that includes validation of the practitioner's registration number. The application procedure may be further secured by receipt and validation of a practitioner-supplied PIN. Moreover, the application procedure typically includes issuance of a client certificate insuring that access to the first and second remote servers occurs from the same practitioner's computer and browser that initiated the application procedure.

Because the data on the servers is separate and secure from each other, access may be granted to either server without identifying any particular patient's medical data. For example, access may be granted to the first remote server, but not to the second server, for the purpose of generating a mailing list of patients without divulging any medical data associated with the patients. Similarly, access may be granted to the second remote server, but not to the first server, for the purpose of conducting investigative analyses involving the medical data without divulging any patient identification information associated with the patients.

For further data security and because each removable memory storage device only has a limited data storage capability, the medical data stored on each removable memory storage device may be automatically erased from the memory storage device after the data is entered into the second remote server. To obtain the medical data, the removable memory storage device is receivable within a hand-held health monitoring device operable for storing the medical data on the removable memory storage device. And to download the medical data to the medical records maintenance system, the removable memory storage device is

receivable within a computer operable for reading the medical data and transmitting it to the second remote server over the Internet.

That the invention improves over the drawbacks of health monitoring and diagnostic systems and accomplishes the advantages described above will become
5 apparent from the following detailed description of the exemplary embodiments and the appended drawings and claims.

BRIEF DESCRIPTION OF THE DRAWINGS

10 FIG. 1A is a front view of a hand-held health monitoring and diagnostic device in an open position.

FIG. 1B is a rear view of the hand-held health monitoring and diagnostic device of FIG. 1 in an open position.

FIG. 2 is a block diagram illustrating a system for remotely producing health reports.

15 FIG. 3 is a functional block diagram of a health monitoring and diagnostic device.

FIG. 4 is a logic flow diagram illustrating a routine for activating a health monitoring and diagnostic device.

20 FIG. 5 is a logic flow diagram illustrating a routine for computing an activation code for a health monitoring and diagnostic device.

FIG. 6 is a logic flow diagram illustrating a routine for verifying a test strip for a health monitoring and diagnostic device.

FIG. 7 is a logic flow diagram illustrating a routine for computing a test strip identification number for a health monitoring and diagnostic device.

25 FIG. 8 is a logic flow diagram illustrating a routine for entering diagnostic program modules into a health monitoring and diagnostic device.

FIG. 9 is a logic flow diagram illustrating a routine for computing immediate diagnostic results in a health monitoring and diagnostic device, and for remotely producing health reports.

30 FIG. 10 is a logic flow diagram illustrating a routine for obtaining cholesterol-related diagnostic information for a health monitoring and diagnostic device.

FIG. 11 is a logic flow diagram illustrating a routine for computing immediate cholesterol-related diagnostic results for a health monitoring and diagnostic device.

FIG. 12 is a logic flow diagram illustrating a routine for remotely producing health reports.

FIG. 13 is a logic flow diagram illustrating a routine for saving medical data to a PIN-secured removable memory storage device for a health monitoring and diagnostic device.

FIG. 14 is a functional block diagram of a system for using a health monitoring and diagnostic device in connection with a secure medical records maintenance system.

FIG. 15 is software architecture diagram illustrating a system for conducting secure communications between a health monitoring and diagnostic device and a secure medical records maintenance system.

FIG. 16 a functional block diagram illustrating security aspects of a secure medical records maintenance system.

FIG. 17 is a logic flow diagram illustrating a process for a communicating with a secure medical records maintenance system.

FIG. 18 is a logic flow diagram illustrating a process for applying for access to a secure medical records maintenance system.

FIG. 19 is a logic flow diagram illustrating a process for logging into a secure medical records maintenance system.

FIG. 20 is an illustration of a "switchboard" user interface in a secure medical records maintenance system.

FIG. 21 is an illustration of an "address" user interface in a secure medical records maintenance system.

FIG. 22 is an illustration of a "billing information" user interface in a secure medical records maintenance system.

FIG. 23 is an illustration of a "cover letter" user interface in a secure medical records maintenance system.

FIG. 24 is an illustration of a "patient selection" user interface in a secure medical records maintenance system.

FIG. 25 is an illustration of a "patient information" user interface in a secure medical records maintenance system.

FIG. 26 is an illustration of a "questionnaire data" user interface in a secure medical records maintenance system.

5 FIG. 27 is an illustration of a "generate reports" user interface in a secure medical records maintenance system.

FIG. 28 is an illustration of typical health assessment charts generated by a secure medical records maintenance system.

10 FIG. 29 is an illustration of an additional health assessment chart generated by a secure medical records maintenance system.

FIG. 30 is a logic flow diagram illustrating a routine for adding a monetary value to a smartcard for use with a health monitoring device.

FIG. 31 is a logic flow diagram illustrating a routine for using a smartcard to pay for a service provided by a health monitoring device.

15

DETAILED DESCRIPTION OF EMBODIMENTS OF THE INVENTION

Hand-Held Health Monitoring and Diagnostic Device

20 Turning now to the figures, in which like numerals refer to like elements through the several figures, FIG. 1A is a front view of a hand-held health monitoring and diagnostic device **10**, which is also referred to as a meter or a LIFESTREAM cholesterol meter. The meter **10** is housed in a clam-shell case **12** including a first compartment **14** and a second compartment **16**. The case **12** may be opened, as shown in FIG. 1A, or closed about a hinge **18**. This allows a
25 patient to close the meter **10** for transportation or storage, and then easily open it for use. When in use, the patient may place the meter **10** in the open position on a flat surface, such as a table or seat, or hold the meter by hand.

Although the meter **10** is shown in a hinged clam-shell, hand-held configuration, it could alternatively be embodied in other configurations, such
30 wall-mounted, built into a movable cart, built into a desktop computer, built into a fixed podium, and so forth. In addition, the hinged clam-shell case could be replaced by a non-hinged case, a separable multi-piece case, a case with a pull-out drawer, a case with a flat cover, a meter that fits into a separate zippered case, and

other types of single- or multi-piece configurations. Many other variations of the meter case configuration will be apparent to those skilled in the art.

The first compartment **14** includes a holder **18** for removably supporting a biological sample gathering device, in this instance a conventional spring-loaded finger stick **20**. Although the holder **18** is shown as a clip with two arms that fit snugly against the finger stick **20**, the holder may have any other configuration suitable for removably supporting the finger stick, such as a channel into which a pencil-like finger stick is inserted, an openable enclosure, snaps, a VELCRO fastener, and the like. If samples other than blood are to be gathered, the first compartment **14** could alternatively house other types of biological sample gathering devices, such as a skin sample collector, a saliva collector, a stool sample collector, and so forth. In addition, the meter **10** may include other types of instruments for gathering test data, for example the meter may be adapted to read non-proprietary test strips, or may incorporate a reusable and/or non-invasive testing device, such as an electrode, blood pressure monitoring device, sonic testing device, thermometer, saliva testing device, optical testing device, and the like.

The first compartment **14** also includes an openable enclosure **24** for storing one or more packages **26** of disposable items. Specifically, each package may contain a test strip **28**, a needle **29** for the finger stick **20**, and an alcohol swipe **30**. These disposable items are tailored for one-time use with the finger stick **20**. If the meter **10** includes biological sample gathering devices other than the finger stick **20**, other types of disposable items may be stored in the enclosure **24**. In addition, the enclosure **24** may have other configurations suitable for storing or holding disposable items, such as a drawer, a tilting channel, a clip, and so forth.

The second compartment **16** houses the electronic components of the meter **10**, including a test strip reader **32**, a display device **34**, a user input device **36**, and one or more memory reading devices **38a-d**. Each of these memory reading devices is configured to receive a corresponding memory device **40a-d**. The second compartment **16** also includes an instructional label **42** located adjacent to the display device **34**. Internally, the second compartment **16** houses a motherboard, an analyzer board, and a data drive that control the functionality of the meter **10**. These internal components are described with reference to FIG. 3,

and the functionality of the meter **10** is described with reference to FIGS. 4-13. Additional functionality of the meter **10** for use with a debit-card type payment system is described with reference to FIGS. 30-31.

The test strip reader **32** may be a GLUCOTREND Basic TIM (test
5 instrument module) assembly No. 1739905-741 manufactured by Boehringer Mannheim, Roche Diagnostics GmbH. This is a commercially-available optical test strip reader suitable for reading chemical test strips carrying human blood samples and producing either blood glucose readings, total cholesterol readings, or both. Alternatively, other suitable types of test strip readers may be included in
10 the meter **10**, and multiple test strip readers may be included in the meter, if appropriate. This may be desirable, for instance, if biological samples other than blood are to be analyzed by the meter. The meter **10** may be configured with additional reusable and/or non-invasive testing devices, such as an electrode, blood pressure monitoring device, sonic testing device, thermometer, saliva testing
15 device, optical testing device, and the like.

The display device **34** may be a conventional liquid crystal display (LCD) configured to display at least two lines of text including at least 14 characters per line. This visual display works in concert with a speaker **35** that beeps to convey audible messages. The speaker may also produce other types of audible messages,
20 such as tones, recorded messages, a simulated human voice, and the like. The meter **10** may also include other types of visual display devices, such as an electronic capacitive matrix, a small video display, or other types of suitable visual display devices. The meter **10** may also include a jack for connecting the meter to external display devices, such as a computer monitor or video display.

25 The user input device **36** may be a keypad with a first key section **44** and a second key section **46**. The first key section **44** includes four keys, a "scroll" key, a "yes" key, a "no" key and an "enter" key. The second key section **46** includes twelve keys, including ten numerical keys, a "clear" key, and an "on/off" key. The user input device **36** may include other key patterns and other types of
30 user input devices, such as a touch-sensitive screen, a voice-recognition device, or other input devices. The meter **10** may also include a jack for connecting the meter to external input devices, such as a keyboard or joystick.

The memory reading devices **38a-d** may be romkey sockets, and the memory devices **40a-d** may be romkeys that removably insert into the sockets. As shown in FIG. 1A, the meter **10** preferably includes four romkey sockets. Nevertheless, the meter **10** could also be configured with only one socket because
5 the romkeys themselves are removable. The romkeys, which store identification, expiration, and calibration data for a corresponding lot of test strips **28**, are desirable because they are small, may be easily packaged with the corresponding test strips, and have an adequate amount of computer-readable memory. But the romkey sockets may be replaced by a magnetic card reader, an optical reader, or
10 another reader suitable for use with a memory storage device that can be easily shipped with a corresponding lot of test strips **28** and has an adequate amount of computer-readable memory.

The instructional label **42** located adjacent to the display device **34** typically includes instructions for entering diagnostic information into the meter **10**. This
15 label may also include instructions for using the meter **10** in concert with a remote health report system, which is described with reference to FIG. 2. As this type of information may be understood best when explained by a physician or pharmacist, the instructional label **42** may be included on meters provided to physicians and pharmacists, but may be excluded from meters provided for home use by
20 individual patients. The programmed functionality of the meter may also be adjusted accordingly.

To use the meter **10**, a patient first opens the meter and removes the finger stick **20** and the package of disposable items **26**. The patient then opens the package, installs the needle **29** in the finger stick **20**, and wipes the alcohol swipe
25 **30** on the area of the finger to be stuck with the needle. The patient also inserts the correct romkey for the test strip **28**, represented by the romkey **40d**, into a corresponding romkey socket **38d** and manipulates the keypad **44**, **46** to indicate to the meter **10** which romkey socket contains the correct romkey. The patient then sticks the selected finger with the finger stick **20**, places a droplet of blood **44**
30 on an indicated area of the test strip **28**, and inserts the test strip **28** into the test strip reader **32**.

The user then manipulates the input device **36** by following prompts displayed on the display device **34** to complete the test. Within minutes, the meter

10 completes the test and displays the test results, such as total cholesterol levels, blood glucose levels or another testing service provided by the meter, on the display device **34**. If appropriate, the user may also manipulate the input device **36** to enter additional diagnostic information into the meter, such as gender, ethnicity, family history of heart disease, personal history of heart disease, personal history of diabetes, personal history of smoking, height, weight, age, blood pressure, and fitness level. Within minutes, the meter **10** performs a diagnostic analysis and produces diagnostic results, such as a medical risk index, a recommended weight loss, a five-year risk of heart attack, a ten-year risk of heart attack, a cardiac age, an extended age, and a risk of stroke. These diagnostic results are also immediately displayed on the display device **34**.

For a blood lipid or cholesterol test, the well-known Framingham Medical Study may provide the methodology used by the meter **10** to produce the diagnostic results from the test results and the diagnostic information. Other methodologies, such as those sanctioned by the National Cholesterol Education Program, the American Heart Association, the American Medical Association, or another appropriate organization may also be used. In fact, the meter **10** may allow the user to select among several alternative diagnostic program modules stored within the meter. These diagnostic program modules may be updated from time to time, and new diagnostic program modules may be added to the meter **10** through the data drive, which is described below.

FIG. 1B is a rear view of the meter **10**, which shows the outside of the meter. The outside of the first compartment **14** includes the manufacturer's name, Lifestream Technologies, Inc., and the meter's trademark, LIFESTREAM. The outside of the second compartment **16** includes a data drive **50**, which includes an opening **52** for receiving a removable memory storage device **54**. For example, the data drive **50** may be a smartcard drive, such as an STM IC: MC33560ADW manufactured by Motorola, and the removable memory storage device **54** may be a smartcard usable with this drive. This smartcard typically includes an electrical contact **56** for reading and writing data and a small microprocessor **58**, which typically controls security aspects of the smartcard. Specifically, the smartcard includes a PIN-protected secure memory and an

unsecure memory. The microprocessor 58 controls the PIN and any other functionality resident on the smartcard.

The smartcard is an advantageous memory storage device because of its small size, its on-card PIN security feature, and its on-board programmable processing unit. Moreover, it is expected that smartcards will become increasingly popular in the near future, and most personal computers will come with factory-installed smartcard drives. Nevertheless, the meter 10 could include other types of data drives, such as a floppy disk drive, an optical disk drive, a removable RAM chip, or any other suitable type of removable memory storage device. Furthermore, the meter 10 could also include a wire-line data port for connecting a cable or another computer station in addition to or as an alternative to the data drive. Similarly, the wire-line data port could be replaced by a wireless communication device, such as a radio-frequency link, a laser link, an infra-red link, and so forth.

The second compartment 16 also includes a battery enclosure 60 housing a removable battery 62 for powering the meter 10. The battery 60, which may be a disposable 9 Volt battery, may be replaced or augmented by an A/C power cord and an appropriate power inverter. The battery 60 also be replaced by a rechargeable battery augmented with a battery charger that may be located within the meter 10 or in a separate enclosure, such as a storage container for housing the meter when it is not in use.

Remote Health Report System

FIG. 2 is a block diagram illustrating a system 100 for remotely producing health reports. The PIN-securable smartcard described above allows multiple patients to use the same meter, and also allows each patient to control access to his or her medical data. The smartcard is typically used as a temporary storage location for one or several test results. From time to time, each patient is expected to download the medical data contained on his or her smartcard to a health report server 102 for permanent storage. The smartcard is then erased (except for the PIN and any other security information contained in the secure memory), which frees the memory space for new medical data. The smartcard may also be used to receive new program modules or new versions of program

modules from the health report server **102** and upload these program modules to the meter **10**.

The system **100** includes the meter **10**, which downloads medical data to the removable memory storage device **54**. This medical data typically includes
5 patient identification information, test results generated by the meter, diagnostic information entered into the meter **10** by the patient, and the diagnostic results **104** computed by the meter and immediately displayed on the meter at the time of the test. If the removable memory storage device **54** is a smartcard, it may then be placed in a converter, such as a conventional smartcard-to-floppy-disk
10 converter **106**, which can be directly inserted into a computer station **108**. The converter **106** will be unnecessary, of course, if the computer station **108** includes a smartcard drive, or if another communication mechanism is employed to transmit information between the computer station **108** and the meter **10**.

Once the medical data arrives at the computer station **108**, it establishes a
15 network connection with the health report server **102**, typically over the Internet **110**. The medical data is then transmitted to the health report server **102**, which may also prompt the user for additional diagnostic and health report information. Specifically, the health report server **102** typically works in concert with the patient's physician or pharmacist, who may provide additional diagnostic
20 information to the server, such as a newly-prescribed drug therapy, other currently-prescribed drugs for the patient, exercise and dietary recommendations, and so forth. Within minutes, the health report server **102** assembles a comprehensive health report **112** that is transmitted back to the computer station **108**, where it may be printed on a local printer **114**. User access procedures and
25 a menu-driven user interface system for generating the health reports is described with reference to FIGS. 17-29.

The comprehensive health report **112** typically includes a data sheet for the newly-prescribed drug, cross-reaction information for the newly-prescribed drug and the other currently-prescribed drugs, weight and total cholesterol goals,
30 exercise and dietary recommendations, any food or activity warnings associated with the overall therapy package, and recommendations for on-going monitoring using the meter. This provides a complete written record of the patient's current

condition, the therapy prescribed by the physician and filled by the pharmacist, and a roadmap for monitoring the patient's progress during the ensuing therapy.

The comprehensive health report may also include additional patient-specific information, such as the diagnostic information and results compiled by the meter, and additional diagnostic and health assessment information compiled by the server. For example, the report may include a trend analysis showing how blood lipid, blood glucose, and weight levels have changed over multiple readings. The report may also include generally-applicable educational information, such as coronary risk factors, dietary guidelines for reducing cholesterol levels, diabetes information, cancer information, and the like. At present, a patient may have to undergo a physical examination, pay thousands of dollars, and wait weeks to obtain a similar comprehensive health report. The network-based comprehensive health analysis and reporting system, working in concert with the LIFESTREAM cholesterol meter, allows the patient to obtain the report within minutes at a fraction of the cost.

The Internet **10** allows a wide variety of users to access the health report server **102**, which allows meters to be deployed in a variety of settings. For example, the accessing computer stations may include a physician's computer station **116**, a pharmacist's computer station **118**, an individual's computer station **120**, and many others. This will allow meters to be effectively deployed for multi-patient use in clinics, physicians' offices and pharmacies, as well as for individual patient or family use in the privacy of their own homes.

Functional Operation of the Meter

FIG. 3 is a functional block diagram of the health monitoring and diagnostic device **10**, which is also referred to as a meter. The meter includes an analyzer board **150**, a mother board **152**, a memory reading device **154** including the romkey sockets **38a-d** and corresponding romkeys **40a-d**, and a user interface **156** including the display device **34**, the speaker **35**, and the input device **36**. The memory reading device **154** and the user interface **156** were described with reference to FIG. 1A. Although the analyzer board **150** and the mother board **152** may be configured as two separate integrated circuit boards, alternatively

they may be combined into a single integrated circuit board, or deployed on more than two integrated circuit boards.

The analyzer board **150** may be part of the GLUCOTREND Basic TIM (test instrument module) assembly No. 1739905-741 manufactured by Boehringer Mannheim, Roche Diagnostics GmbH. This board includes the test strip reader **32**, which was described with reference to FIG. 1A, a test instrument module **158**, and a romkey driver **160**. The test strip reader **32** reads the test strip **28** carrying the blood sample **44**. The romkey driver **160** reads the calibration data for the test strip **28** from a corresponding romkey, such as the romkey **40d**, and the test instrument module **158** computes the test results from the test strip reading and the calibration data. These test results are then passed to the motherboard **152**.

The motherboard **152** includes a non-interruptible battery **162**, such as a lithium battery. The non-interruptible battery **162**, which powers the on-board clock **164**, is distinct from the power supply battery **62** (shown on FIG. 1B). The additional non-interruptible battery **162** allows the clock to continue functioning even when the power supply battery **62** runs down or is removed from the meter **10**. The motherboard **152** also includes a processor **166** and a memory **168** that control the functionality of the meter **10**. Specifically, the memory **168** stores and the processor **166** implements a meter activation routine, test strip verification routine, diagnostic routines, and a read/write security routine. These program modules are described with reference to FIGS. 4-13. The motherboard **152** also includes a data drive **170** that reads data from and writes data to the removable data storage device **54**, such as a smartcard.

The mother board processor **166** may be a 40-pin DIP CPU Model No. MC68HC705C9ACP manufactured by Motorola. The data drive **170** may be a STM IC Model No. MC33560ADW ISO read/control card manufactured by Motorola, supported by an ISO Card Socket, Model No. 145206-3 physical interface manufactured by AMP. However, any of these specific devices may be replaced by equivalent devices capable of performing the functionality of the meter **10** described in this specification.

To verify test strips, the romkey driver **160** reads a code number from a romkey, such as the romkey **40d**, installed in the meter **10**. The code number typically includes the lot number for the corresponding test strips and a test type

ID stored on the romkey by the manufacturer of the meter **10**. The test strip **28** has an associated test strip identification number **172** that is mathematically derived from the code number and printed on the test strip itself, the packaging for the test strip, or a tag packaged with the test strip. The meter **10** prompts the user to enter the test strip identification number **172** into the meter using the keypad **46**.

The meter **10** also reads the code number from the memory device, mathematically derives a test strip identification number corresponding to the code number, compares the received test strip identification number to the derived test strip identification number, and activates the meter for use with the test strip only if the received test strip identification number corresponds to the derived test strip identification number. The romkey **40d** also stores an expiration date for the corresponding test strip **28**. The meter **10** reads the expiration date and activates for use with the test strip **28** and the romkey **40d** only if the expiration date is prior to a current date read by the meter from the internal clock **164**.

FIGS. 4-14 are logic flow diagrams illustrating examples of the functionality that may be implemented by the meter **10** and the system **100** for remotely generating health reports. The following description of these logic flow diagrams will also refer to the elements shown on FIGS. 2 and 3. It should be understood that these examples illustrate the use of these components in producing total cholesterol tests and related health reports. In particular, the specific diagnostic information gathered and the specific diagnostic results described are those associated with the well-known Framingham Medical Study, which is incorporated herein by reference. Although this particular program module illustrates the operation of an illustrative embodiment of the invention, those skilled in the art will understand that the meter **10** and the system **100** could be programmed with additional and different program modules.

In addition, as the meter **10** is configured with multiple romkey sockets and programmed to accept additional program modules in the future, it will be appreciated that similar functionality may be implemented in the future for blood glucose tests and diabetes-related health reports, AIDS tests and related health reports, cancer tests and related health reports, and so forth. Additional

functionality of the meter **10** for use with a debit-card type payment system is described with reference to FIGS. 30-31.

Meter Activation

5 FIG. 4 is a logic flow diagram illustrating a routine **400** for activating the meter **10**. This routine requires that a user enter a proper activation code to activate the meter **10**. This typically requires that the user call the manufacturer, which provides an opportunity to verify the meter's authenticity, set up a data file for the meter in the health report server, and tell the user how to update the meter
10 software, if necessary. If a software update is indicated, the user may be instructed to activate the meter, initialize a smartcard, load the smartcard into a computer station, and establish a network connection with the health report server. The server can then download the new software (e.g., new version of an existing software module or a new software module) to the smartcard **54**, which,
15 in turn, can be placed back in the meter **10**. The new software can then be uploaded to the meter.

 In step **402**, the meter **10** is programmed at the factory by setting the internal clock **164** to the correct date and storing an activation code algorithm within the meter. Step **402** is followed by routine **404**, in which the activation
20 site, typically the manufacturer, computes an activation code for the current day. That is, each day the activation site computes a new activation code that is valid only for that day. The activation code is computed using the same algorithm that was stored in the meter **10** in step **402**. This algorithm is described below with reference to FIG. 5.

25 Step **404** is followed by routine **406**, in which the activation site receives an activation communication for the purpose of activating a meter **10**. The user of the meter typically places this communication by placing a telephone call to a telephone number (e.g., a "1-800" or other toll-free telephone number) on the meter or on the packaging or documentation provided with the meter. Step **406**
30 is followed by routine **408**, in which the activation site delivers the activation code for the current date to the calling user. Step **408** is followed by routine **410**, in which the calling user enters the activation code into the meter **10** by manipulating the user input device **36**.

Step **410** is followed by step **412**, in which the meter **10** verifies the received activation code by computing an activation code using its on-board activation code algorithm and the current date read from its internal clock **164**. In doing so, the meter **10** uses the same algorithm that was used by the activation site to compute the activation code that was delivered to the user and entered into the meter (i.e., the algorithm described with reference to FIG. 5). Step **412** is followed by routine **414**, in which the meter **10** compares the received activation code to the computed activation code. Step **414** is followed by routine **416**, in which the meter **10** determines whether the activation is verified by determining whether the received activation corresponds to the computed activation code.

If the activation is not verified, the "NO" branch is followed from step **416** to step **418**, in which the meter **10** determines whether a timeout condition or an allowed number of tries has been reached. If a timeout condition or an allowed number of tries has not been reached, the "NO" branch loops from step **418** to step **410**, and the meter **10** displays an error and prompts the user to reenter the activation code. If a timeout condition or an allowed number of tries has been reached, the "YES" branch is followed from step **418** to the "END" step, and the meter **10** is not activated.

Referring again to step **416**, if the activation is verified, the "YES" branch is followed from step **416** to step **420**, in which the meter **10** is activated. Step **412** is followed by the "END" step, in this case with the meter **10** activated.

FIG. 5 is a logic flow diagram illustrating a routine **404** for computing an activation code for the meter **10** or the activation site, referred to collectively below as the "computing entity." Routine **404** begins following step **402** shown on FIG. 4. In step **502**, the computing entity gets the day, month, and year (six digit decimal) from the internal clock **164**. Step **502** is followed by step **504**, in which the computing entity combines the six digits defining the date into a hex value (six digit hex). Step **504** is followed by step **506**, in which the computing entity applies a predetermined mathematical operation to this hex value to compute a new hex value (new six digit hex).

Step **506** is followed by step **508**, in which the computing entity applies a mask to the new six digit hex value to obtain a four digit hex value (new four digit hex). In other words, the computing entity selects a predetermined four of the six

digits for further processing. Step **508** is followed by step **510**, in which the computing entity converts this new four digit hex value to a decimal value (four digit decimal value). Step **510** is followed by step **512**, in which the computing entity relocates one or more of the digits and adds one or more null numbers at predefined locations (new four digit decimal value). Step **512** is followed by step **514**, in which the computing entity uses the resulting new four digit decimal number as the activation code. Step **514** is followed by the "RETURN" step, which goes to step **406** on FIG. 4.

The mathematical operation applied in step **506** may be any of a variety of algorithms designed to quasi-randomize the result. For example, the digits may be grouped into subsets that, in turn, are used in one or more linear mathematical operations, such as addition, subtraction, multiplication, division, raising to a power, raising to a fraction, and the like. For example, a polynomial formula may be applied to the digits or subsets of the digits. The digit shuffling operation applied in step **512** is also applied to quasi-randomize the result. Many other types of quasi-randomizing methodologies will be apparent to those skilled in the art.

Test Strip Validation

FIG. 6 is a logic flow diagram illustrating a routine **600** for verifying a test strip for the meter **10**. In step **602**, the meter **10** is programmed at the factory with a test strip validation algorithm. Step **602** is followed by step **604**, in which a romkey with a corresponding lot of test strips is programmed with a test type ID. This test type ID, together with a lot number for the test strips placed on the romkey by the test strip manufacturer, forms a four-digit code number that is resident on the romkey when it leaves the factory.

Step **604** is followed by routine **606**, in which a test strip ID is mathematically derived from the romkey code number. Routine **606** is described below with reference to FIG. 7. Routine **606** is followed by step **608**, in which the test strip ID is printed on the lot of test strips corresponding to the romkey, on the packaging for the test strips, or tags that are packaged with the test strips. The test strips are then packaged with corresponding romkeys and distributed for use with various meters. It should be understood that many identical romkeys may be produced for each lot of test strips because one romkey is included in each

distribution-sized package, and a production-sized lot of test strips may be many times larger than a distribution-sized package.

Step 608 is followed by step 610, in which the user of a meter 10 puts a romkey into the meter and enters the corresponding test strip ID into the meter using the user input device 36. This is the test strip ID that was printed on the test strip, its packaging, or a tag packaged with the test strip in step 608. Step 610 is followed by step 612, in which the meter 10 gets the current date from the internal clock 164, reads the code number and expiration date from the romkey, and computes a test strip ID number based on the code number. The test strip ID number is derived from the code number using the same algorithm that was used to compute the test strip ID number at the factory in routine 606.

Step 612 is followed by step 614, in which the meter 10 determines whether the current date is prior to the expiration date read from the romkey. If the current date is not prior to the expiration date read from the romkey, the "NO" branch is followed to the "END" step, and the meter is not activated for use with the instant romkey. If the current date is prior to the expiration date read from the romkey, the "YES" branch is followed to step 616, in which the meter 10 compares the received test strip (input by the user) to the computed test strip ID (derived from the code number read from the romkey). Step 616 is followed by step 618, in which the meter 10 verifies the test strip if the received test strip corresponds to the computed test strip ID.

If the meter 10 verifies the test strip, the "YES" branch is followed from step 618 to step 622, in which the meter 10 activates for use with the instant test strip and romkey. Step 622 is followed by the "END" step. If the meter 10 does not verify the test strip, the "NO" branch is followed from step 618 to step 620, in which the meter 10 determines whether a timeout condition or an allowed number of tries has been reached. If a timeout condition or an allowed number of tries has not been reached, the "NO" branch loops from step 620 to step 610, and the meter 10 displays an error and prompts the user to reenter the test strip ID and/or insert the correct romkey. If a timeout condition or an allowed number of tries has been reached, the "YES" branch is followed from step 620 to the "END" step, and the meter 10 is not activated for the instant test strip and romkey.

Test Strip ID Assignment

FIG. 7 is a logic flow diagram illustrating routine **606** for computing a test strip ID for the meter **10**. Routine **606** begins following step **604** shown on FIG. 6. In step **702**, the meter **10** reads the code number (i.e., lot number plus test type ID) from the romkey. Step **702** is followed by step **704**, in which the meter **10** converts this four digit decimal value to 16-bit binary value. Step **704** is followed by step **706**, in which the meter **10** applies a predetermined mathematical operation to this 16-bit binary value.

Step **706** is followed by step **708**, in which the meter **10** checks a value at a predetermined bit of the 16-bit binary value. Step **708** is followed by step **710**, in which the meter **10** performs a conditional bit swap. That is, a first type of bit swap is performed if the value of the checked bit is a "1" and a second type of bit swap is performed if the value of the checked bit is a "0." Step **710** is followed by step **712**, in which the meter **10** converts the resulting number to a four digit decimal value. Step **712** is followed by step **714**, in which the meter **10** uses the resulting four digit decimal number as the test strip ID. Step **714** is followed by the "RETURN" step, which goes to step **608** on FIG. 6.

The mathematical operation applied in step **706** may be any of a variety of algorithms designed to quasi-randomize the result. For example, the digits may be grouped into subsets that, in turn, are used in one or more linear mathematical operations, such as addition, subtraction, multiplication, division, raising to a power, raising to a fraction, and the like. For example, a polynomial formula may be applied to the digits or subsets of the digits. The digit shuffling operation applied in steps **708** and **710** is also applied to quasi-randomize the result. Many other types of quasi-randomizing methodologies will be apparent to those skilled in the art.

Diagnostic Data Entry

FIG. 8 is a logic flow diagram illustrating a routine **800** for entering diagnostic program modules into the meter **10**. In step **802**, the meter **10** is programmed at the factory with one or more initial diagnostic algorithms. In addition, a computer maintained at the programming site, typically the

manufacturer, is programmed with a cross reference table indicating the serial number for the meter and the version number of each diagnostic program module installed in the meter **10**.

At some later point, a user downloads medical data from the meter **10** to a smartcard in the course of normal meter use. At this time, the serial number for the meter is also stored on the smartcard. The smartcard is then read by a computer station, which establishes a network connection with a programming server, which may be the same as, or coordinated with, the health report server **102** shown on FIG. 2. At step **804**, the programming server receives the data that was stored on the smartcard. Step **804** is followed by step **806**, in which the programming server instructs the computer station to erase the unsecured data stored on the smartcard. That is, the computer station erases the medical data stored in the unsecure memory but does not erase the PIN or any data stored in the secure memory.

Step **806** is followed by step **808**, in which the programming server gets the serial number for the meter **10** from the received data and looks up the version numbers for the diagnostic program modules installed on the meter. Step **808** is followed by step **810**, in which the programming server determines whether the meter **10** includes the latest version of all of the program modules that should be installed on the meter. If the meter **10** includes the latest version of all of the program modules that should be installed on the meter, the "YES" branch is followed to the "END" step, and the meter software is not updated.

If, on the other hand, the meter **10** does not include the latest version of all of the program modules that should be installed on the meter, the "NO" branch is followed to step **812**, in which the programming server loads new diagnostic program modules, new versions of diagnostic program modules, or updates for existing program modules on the smartcard. At some later point in step **814**, the smartcard is placed back in the meter **10**. At this point, step **814** is followed by step **816**, in which the new diagnostic program modules, new versions of diagnostic program modules, or updates for existing program modules stored on the smartcard are uploaded to the meter **10**. Step **816** is followed by the "END" step.

Diagnostic Analysis

FIG. 9 is a logic flow diagram illustrating a routine **900** for computing immediate diagnostic results in the meter **10**, and for remotely producing health reports. In step **902**, the meter **10** prompts the user to select a romkey and insert
5 a corresponding test strip, with a blood sample, into the meter. Step **902** is followed by step **904**, in which the meter **10** prompts the user to select a desired diagnostic program module corresponding to the type of test strip inserted into the meter. Step **904** is followed by step **906**, in which the meter **10** performs the test strip verification algorithm shown on FIG. 6 and, if the test strip is verified,
10 obtains test results.

Step **906** is followed by routine **908**, in which the meter **10** prompts the user for additional diagnostic information. Routine **908** is described below with reference to FIG. 10. Routine **908** is followed by routine **910**, in which the meter **10** computes immediate diagnostic results. Routine **910** is described below with
15 reference to FIG. 11. Routine **910** is followed by step **912**, in which the meter **10** displays the diagnostic results on the display device **34**. Step **912** is followed by step **914**, in which the meter **10** stores the test results, the diagnostic information, and the diagnostic results (and also stores the meter's serial number) on the smartcard.

20 At some later point, in step **916** the user reads the smartcard with a computer station. Step **916** is followed by step **918**, in which the computer station establishes a network connection with the health report server **102**. Step **918** is followed by step **920**, in which the computer station downloads the data from the smartcard to the health report server **102**. Step **920** is followed by
25 routine **922**, in which health report server **102** receives additional diagnostic and health report information from the user of the computer station and compiles a personal health report based on the data received from the computer station. Routine **922** is described below with reference to FIG. 12. Routine **922** is followed by step **924**, in which the health report server **102** transmits the health
30 report to the computer station, where the health report is printed and delivered to the patient. Step **924** is followed by the "END" step.

FIG. 10 is a logic flow diagram illustrating routine **908** for obtaining cholesterol-related diagnostic information for the meter **10**. Routine **908** begins

following step **906** shown on FIG. 9. In step **1002**, the meter **10** prompts for and receives the patient's gender. Step **1002** is followed by step **1004**, in which the meter **10** prompts for and receives the patient's ethnicity. Step **1004** is followed by step **1006**, in which the meter **10** prompts for and receives an indication of the patient's family history of heart disease. Step **1006** is followed by step **1008**, in which the meter **10** prompts for and receives an indication of the patient's personal history of heart disease.

Step **1008** is followed by step **1010**, in which the meter **10** prompts for and receives an indication of whether the patient is a type-1 diabetic. Step **1010** is followed by step **1012**, in which the meter **10** prompts for and receives an indication of whether the patient is a type-2 diabetic. Step **1012** is followed by step **1014**, in which the meter **10** prompts for and receives an indication of whether the patient is a smoker.

Step **1014** is followed by step **1016**, in which the meter **10** prompts for and receives an indication of the patient's height. Step **1016** is followed by step **1018**, in which the meter **10** prompts for and receives an indication of the patient's weight. Step **1018** is followed by step **1020**, in which the meter **10** prompts for and receives an indication of the patient's age. Step **1020** is followed by step **1022**, in which the meter **10** prompts for and receives an indication of the patient's blood pressure. Step **1022** is followed by step **1024**, in which the meter **10** prompts for and receives an indication of the patient's fitness. Step **1024** is followed by the "RETURN" step, which goes to routine **910** shown on FIG 9. It will be appreciated that the preceding list of diagnostic information is only illustrative of the type of information that may be gathered, and that less data, different data, or more data could be gathered, if desired.

FIG. 11 is a logic flow diagram illustrating routine **910** for computing immediate cholesterol-related diagnostic results for the meter **10**. Routine **910** begins following routine **908** shown on FIG. 9. In step **1102**, the meter **10** prompts for and receives a selection of a diagnostic algorithm resident on the meter. Step **1102** is followed by step **1104**, in which the meter **10** gets total cholesterol test results from the test instrument module **158** or, if they are not available, prompts the user to input the test results manually.

Step **1104** is followed by step **1106**, in which the meter **10** calculates and displays a medical risk index associated with heart disease or heart attack, such as “very high,” “high,” “moderate,” “low,” or “very low.” Step **1106** is followed by step **1108**, in which the meter **10** calculates and displays a recommended weight loss, if appropriate. Step **1108** is followed by step **1110**, in which the meter **10** calculates and displays a 5-year cardiac risk (e.g., risk of cardiac arrest in five years is 10%). Step **1110** is followed by step **1112**, in which the meter **10** calculates and displays a 10-year cardiac risk (e.g., risk of cardiac arrest in ten years is 20%).

Step **1112** is followed by step **1114**, in which the meter **10** calculates and displays a cardiac age (to compare against the patient’s chronological age). Step **1112** is followed by step **1114**, in which the meter **10** calculates and displays an extended cardiac age (e.g., cardiac age compared to chronological age for five, ten, fifteen, etc. years into the future). Step **1116** is followed by step **1118**, in which the meter **10** calculates and displays a medical risk index associated with stroke, such as “very high,” “high,” “moderate,” “low,” or “very low.” Step **1118** is followed by the “RETURN” step, which goes to step **912** on FIG. 9. It will be appreciated that the preceding list of diagnostic results is only illustrative of the type of information that may be generated, and that less data, different data, or more data could be generated, if desired.

Remote Health Report Generation

FIG. 12 is a logic flow diagram illustrating routine **922** for remotely producing health reports. Routine **922** begins following step **920** shown on FIG. 9. In step **1202**, in which the health report server **102** prompts for and receives a new drug therapy, such as a cholesterol-lowering prescription. Step **1202** is followed by step **1204**, in which the health report server **102** gets a prescription data sheet for the drug therapy. This data sheet typically includes instructions for taking the prescription, such as dosage, times to take each dose, whether to take with food or liquid, whether to avoid driving, pregnancy-related instructions, foods to avoid, and so forth.

Step **1204** is followed by step **1206**, in which the health report server **102** prompts for and receives information regarding any other prescription drugs that

the patient is currently taking. Step 1206 is followed by step 1208, in which the health report server 102 gets cross-reaction information regarding the new drug therapy and the other current drug prescriptions. Step 1208 is followed by step 1210, in which the health report server 102 prompts for and receives a specific description of exercise therapy for the patient, or a selection of standard exercise sections for inclusion in the health report. Step 1210 is followed by step 1212, in which the health report server 102 prompts for and receives diet therapy for the patient, or a selection of a standard diet section for inclusion in the health report. Step 1212 is followed by step 1214, in which the health report server 102 prompts for and receives indications of additional standard sections for inclusion in the health report.

Step 1214 is followed by step 1216, in which the health report server 102 assembles the preceding information for inclusion in a health report. Step 1216 is followed by step 1218, in which the health report server 102 prompts for and receives trend analysis selections. Step 1218 is followed by step 1220, in which the health report server 102 prepares the selected trend analysis, such as total cholesterol and blood glucose levels over a series of tests. The trend analysis may be provided alone or as part of the health report. Step 1220 is followed by step 1222, in which the health report server 102 assembles the preceding information into a health report and/or trend analysis report. Step 1222 is followed by the "RETURN" step, which goes to step 924 shown on FIG. 9. It will be appreciated that the preceding list of health report information is only illustrative of the type of information that may be compiled, and that less data, different data, or more data could be compiled, if desired.

Smartcard PIN Security

FIG. 13 is a logic flow diagram illustrating a routine 1300 for saving medical data to the PIN-secured removable memory storage device 54, represented by a smartcard. In step 1302, a user places the smartcard in the meter 10. Step 1302 is followed by step 1304, in which the meter 10 prompts for and receives a PIN from the user. Step 1304 is followed by step 1306, in which the meter 10 reads the PIN field on the smartcard. Step 1306 is followed by step

1308, in which the meter 10 determines whether the PIN field is occupied (i.e., whether a PIN has been previously stored on the smartcard).

If the PIN field is not occupied, the "NO" branch is followed from step 1308 to step 1310, in which the meter 10 stores the PIN received from the user in the PIN field on the smartcard. Step 1310 is followed by step 1314, in which the meter 10 activates for use with the current smartcard. Step 1314 is followed by the "END" step.

Referring again to step 1308, if the PIN field is occupied, the "YES" branch is followed to step 1312, in which the meter 10 determines whether the PIN received from the user matches the PIN stored in the PIN field on the smartcard. If the PIN received from the user matches the PIN stored in the PIN field on the smartcard, the "YES" branch is followed to step 1314, in which the meter 10 activates for use with the current smartcard. Step 1314 is followed by the "END" step.

Referring again to step 1312, if the PIN received from the user does not match the PIN stored in the PIN field on the smartcard, the "NO" branch is followed to step 1314, in which the meter 10 determines whether a timeout condition or an allowed number of tries has been reached. If a timeout condition or an allowed number of tries has not been reached, the "NO" branch is followed to step 1318, in which the meter 10 displays an error and prompts the user to reenter the activation code. From step 1318, routine 1300 loops to step 1312. If a timeout condition or an allowed number of tries has been reached, the "YES" branch is followed from step 1316 to the "END" step, and the meter 10 is not activated for use with the instant smartcard.

Secure Medical Records Maintenance System

FIG. 14 is a functional block diagram of a system for using the meter 10 in connection with a secure medical records maintenance system 1400. The meter 10 stores a patient's test results and diagnostic information on a removable memory storage device, such as the smartcard 54. A conventional interface 106 may then be used to interface the smartcard 54 with a conventional desktop, laptop or other type of computer 108, which in turn communicates with other computers over a network-based computer system, such as the Internet 110. As

discussed previously, this Internet link may be used to produce a printed health report booklet **112** including a health assessment based on a patient's test results and diagnostic information, which a health-care provider typically provides to the patient.

5 Although this secure medical records maintenance system **1400** is specifically adapted for use with the meter **10**, it may be used to store any type of electronic medical records, and is particularly convenient for storing a wide range of electronic medical data generated remotely from the hospital or doctor's office environment. In addition, the secure medical records maintenance system **1400**
10 may read medical data stored on memory storage devices other than the smartcard **54**, and may operate over computer networks other than the Internet **100**.

 The secure medical records maintenance system **1400**, which is presently known as the "PRIVALINK" system, operates in connection with a software
15 module **1402** installed on the accessing computer **108**. This software is presently known as the "PRIVALINK" user site software. The server-side "PRIVALINK" system **1400** includes a first remote server **1404** that stores patient identification information, a second remote server **1406** that stores patient medical data, an encryption/decryption module **1408** that implements encryption and other
20 security-related functions, and a booklet generation module **1410**, which produces printed health report booklets **112** based on the information stored in the servers **1404**, **1406**. The patient identification information and medical data are maintained in separate, secure servers **1404**, **1406** to prevent correlation of a specific patient's medical data with the associated patient identification
25 information.

 Because the data on the servers **1404**, **1406** is separate and secure from each other, access may be granted to either server without identifying any particular patient's medical data. For example, access may be granted to the first remote server **1404**, but not to the second server **1406**, for the purpose of
30 generating a mailing list of patients without divulging any medical data associated with the patients. Similarly, access may be granted to the second remote server **1406**, but not to the first server **1404**, for the purpose of conducting investigative

analyses involving patient medical data without divulging any patient identification information associated with the medical data.

For further data security and because each smartcard **54** only has a limited data storage capability, the medical data stored on each smartcard may be automatically erased from the smartcard after the data is entered into the second remote server **1406**. To obtain the medical data, the smartcard **54** is received within the meter **10**, which stores the medical data on the smartcard. And to download the medical data to the medical records maintenance system **1400**, the removable memory storage device is receivable within the interface **106** for communication with the computer **108**, which is operable for transmitting the medical data to the second remote server **1406** over the Internet **110**.

FIG. 15 is software architecture diagram illustrating a system for conducting secure communications between the computer **108** and the secure medical records maintenance system **1400**. The "PRIVALINK" user site software **1402** includes an encryption/decryption module **1412** that implements encryption/decryption services on the client computer **108**. A corresponding encryption/decryption module **1408** on the secure medical records maintenance system **1400** implements complimentary encryption/decryption services on the server side, typically using a well-known encryption/decryption package, such as that known as "BLOWFISH" or "DES." The server-side encryption/decryption module **1408** also maintains access to a table of valid professional registration numbers **1411**, typically DEA numbers, received from the registering agency. This table is used to validate the professional registration numbers of practitioners attempting to access the secure medical records maintenance system **1400**.

The server-side encryption/decryption module **1408** also maintains a record of all transactions with accessing computers for future analysis. In addition, all data transfers from the secure medical records maintenance system **1400** to the user site **108**, such as health report booklets **112**, are encrypted/decrypted through an Internet secure sockets layer connection **1500**, which is well known to those skilled in the art. These encryption/decryption services prevent theft or inadvertent loss of patient medical data through unintended transmission or extraction of communications occurring on the Internet **110**.

The "PRIVALINK" user site software **1402** also implements client application, certification and login processes for accessing the secure medical records maintenance system **1400**. The application, certification and login process is described below with reference to FIGS. 17-19. Upon successful login to the secure medical records maintenance system **1400**, the "PRIVALINK" user site software **1402** implements a menu-driven user interface system **1414** for conducting communications between the practitioner computer **108** and the secure medical records maintenance system **1400**. This interface system includes a health care provider data entry screen **2100** shown in FIG. 21, a patient address data entry screen **2500** shown in FIG. 25, a patient test data entry screen **2600** shown in FIG. 26, a booklet selection screen **2700** shown in FIG. 27, and a booklet preview/printing screen **2800**. The menu-driven user interface system **1412** also includes other related user interface screens. The operation of the menu-driven user interface system **1412** is described below with reference to FIG. 17.

FIG. 16 a functional block diagram illustrating security aspects of the secure medical records maintenance system **1400**. The secure medical records maintenance system **1400** preferably includes a large number of smartcards **54a-n**, which operate in concert with a large number of meters **10a-n**. Although each smartcard **54a-n** is preferably used to store medical data for an associated patient, each card could be used to store medical data for multiple patients. Each smartcard **54** also stores a patient-specified personal identification number (PIN), and may store PINs for multiple patients if the smartcard is configured to store medical data for multiple patients. Each PIN is used to gain access to a secure storage area on the smartcard **54**, which stores an associated patient identification number and medical records identification number, which are assigned by the secure medical records maintenance system **1400**.

The first remote server **1404** of the secure medical records maintenance system **1400** stores patient identification information indexed by patient identification numbers, and the second remote server **1406** stores patient medical data indexed by the medical records identification numbers. Typically, each patient identification number and medical identification number is a unique number assigned by the operator of the secure medical records maintenance

system **1400** upon entry of each patient into the system. For example, the patient identification numbers and medical identification numbers may be social security numbers or sequential registration numbers. Alternatively, the patient identification numbers and medical identification numbers may be unique numbers
5 computed by a non-repeating pseudo-random number generator. The patient identification number and the medical records identification number may be 16-bit hexadecimal or another suitable value. In addition, either or both of the patient identification number and the medical records identification number may be a global user identification number (GUID), which is a secure communication key
10 generated by well-known secure encryption systems, such as currently available private-key and public-key encryption systems including "BLOWFISH," "DES" and others. Many other schemes for assigning unique patient identification numbers will become evident to those skilled in the art.

For security purposes, the medical data maintained in the second remote
15 server **1406** cannot be correlated to the associated patient identification information maintained in the first remote server **1408** based on the information contained in the first and second remote servers. To allow correlation of the data stored in the two servers **1404**, **1406** the secure medical records maintenance system **1400** includes a correlation table **1600** uniquely associating each medical
20 records identification number with a particular one of the patient identification numbers. The correlation table **1600** for a particular patient typically resides in the PIN-secured storage area on the patient's smartcard **54**. The correlation table **1600** for a practitioner's patients may also reside on the practitioner's computer **108**, such as a doctor's or pharmacist's computer, that is associated with a
25 licensed medical practitioner having an assigned professional registration number (DEA number). Each practitioner's correlation table **1600** is preferably encrypted and maintained in a secure file. The proprietor of the secure medical records maintenance system **1400** may also maintain a complete back-up correlation table **1600**, typically in a secure encrypted file located on a separate file
30 server.

For further security, the first and second remote servers **1404**, **1406** are accessed by the practitioners' computers **108a-n** through encrypted communications secured by an application procedure that includes validation of

the practitioner's registration number (DEA number). This access will be limited to medical records and patient identification information associated with the accessing practitioner. In other words, each practitioner will only have access to his or her patients' medical records and patient identification information.

5 Similar access procedure may be implemented for individual patients, except that access will be limited to that particular patient's medical records and patient identification information. For example, individual patients may register in advance with the proprietor of the system **1400**, which will issue each patient a unique registration number. In this case, the first and second remote servers
10 **1404**, **1406** may accessed by computers operated by the individual patients through encrypted communications secured by an application procedure that includes validation of the patient's registration number.

For both practitioners and individual patients, the application procedure may be further secured by receipt and validation of a client-supplied PIN.
15 Moreover, the application procedure typically includes issuance of a client certificate insuring that access to the first and second remote servers occurs from the client computer that initiated the application and certification process.

Those skilled in the art will appreciate that the data distribution system implemented by the secure medical records maintenance system **1400** includes
20 many aspects of data security. For example, a patient's medical records identification number cannot be obtained from the patient's smartcard **54** without access to the patient-assigned PIN. In addition, while the patient's medical data is indexed by the patient's medical records identification number in the second server **1406**, the patient's name and other identification information cannot be
25 retrieved from the first server **1404** using this data. Similarly, a hacker obtaining access to one or both of the servers **1404**, **1406** cannot correlate patient identification information with patient medical data. In addition, the correlation data for the entire secure medical records maintenance system **1400** is distributed among the various smartcards **54a-n** and practitioner computers **108a-n**
30 registered for use with the system. Thus, a hacker cannot obtain the correlation data for any single patient, much less the entire database, through access to the centrally-maintained data servers **1404**, **1406**.

For further security, the secure medical records maintenance system **1400** cannot be accessed from the a medical practitioner's computer **108** without knowledge of the proper practitioner-assigned PIN. And all communication between the practitioner's computer **108** and the secure medical records maintenance system **1400** are encrypted for transmission security. Furthermore, each correlation table **1600**, which provides the link between patient identification numbers and medical records identification numbers for a particular practitioner, may itself be encrypted, with the key to this encryption stored in a separate location. For example, this encryption key may be a practitioner PIN or GUID stored on a PIN-secured area the practitioner's computer **108** or on a smartcard **54** assigned to the practitioner.

FIG. 17 is a logic flow diagram illustrating an illustrative process **1700** for communicating with the secure medical records maintenance system **1400**. In routine **1702**, a medical practitioner conducts a client application procedure to obtain a secure client certificate, which is also known as a "CA" or certificate-authentication. Routine **1702** is described in greater detail with reference to FIG. 18. Routine **1702** is followed by a client login procedure **1704**, which is described in greater detail with reference to FIG. 19. These procedures ensure that access to the secure medical records maintenance system **1400** is limited to registered medical practitioners using the same computer and browser that the practitioner used to obtain a secure client certificate through the application procedure. That is, any attempt to access the secure medical records maintenance system **1400** by a person without a valid DEA number, or from a non-certified computer or browser, will be rejected. As noted previously, similar procedures may be implemented for allowing individual patients to access their records on the secure medical records maintenance system **1400**.

FIG. 18 is a logic flow diagram illustrating routine **1702** used by a medical practitioner to apply for access to the secure medical records maintenance system **1400**. Routine **1702** begins at the start of FIG 17. In step **1802**, the client-side PRIVALINK software **1402** issues a client application request to the server-side encryption/decryption module **1408**, which initiates a transaction document record for the communication. That is, all communications between the client-side PRIVALINK software **1402** and the server-side encryption/decryption module

1408 are recorded for future analysis by the server-side encryption/decryption module **1408**. The client application request is typically received by an Internet link to a specified web page associated with the server-side encryption/decryption module **1408**. The medical practitioner can obtain the proper Internet address by placing a telephone call to the operator of the system **1400**. The phone number, and possibly the Internet address, will be printed on each meter **10** and may also be available through other forms of notification, such as advertisement and direct-mail notification to licensed practitioners (e.g., physicians and pharmacists).

Step **1802** is followed by step **1804**, in which the PRIVALINK software **1402** displays an application screen and receives client input including the practitioner's professional registration number, typically a DEA number. Step **1804** is followed by step **1806**, in which the PRIVALINK software **1402** receives a certification request including selection of an encryption type. This is typically associated with completion of the data-entry fields of the application screen and selection of a "submit" control item. Step **1806** is followed by step **1808**, in which the PRIVALINK software **1402** downloads an encryption program module from the server-side encryption/decryption module **1408**, such as an "ACTIVE X" control or applet, to the client's computer **108**. This encryption program module typically includes a "key" or nugget of information to be stored in the accessing browser. The server-side encryption/decryption module **1408** can later check an accessing computer for the presence of the "key" to ensure that the accessing computer and browser is the same as the one going through the application procedure.

Step **1808** is followed by step **1810**, in which the server-side encryption/decryption module **1408** validates the registration number (e.g., DEA number) entered by the applicant, typically by comparing the received registration number to the table of valid registration numbers **1411** received from the registration authority (e.g., table of valid DEA numbers). If the registration number is properly validated, step **1810** is followed by step **1812**, in which the server-side encryption/decryption module **1408** transmits an e-mail message to the PRIVALINK software **1402** including a URL for accessing the secure servers **1404**, **1406**. Step **1812** is followed by step **1814**, in which the PRIVALINK software **1402** transmits a client link to the designated URL. Upon receipt of the

link, the server-side encryption/decryption module **1408** checks for the presence of the “key” to ensure that the accessing computer and browser is the same as the one that went through the application procedure.

If the “key” is properly validated, step **1814** is followed by step **1816**, in which the server-side encryption/decryption module **1408** transmits a “server root CA” and a “client certificate” to the PRIVALINK software **1402** on the client’s computer **108**. Step **1816** is followed by step **1818**, in which the PRIVALINK software **1402** links the applicant to a secure area of the encryption/decryption module **1408**, which validates the presence of the server root CA. If the server root CA is properly validated, step **1818** is followed by step **1820**, in which the PRIVALINK software **1402** prompts the user to identify the client certificate for use in the transaction. If the client certificate is properly validated, step **1820** is followed by step **1822**, in which the PRIVALINK software **1402** links the applicant to a login screen. Step **1822** is followed by step **1824**, in which the encryption/decryption module **1408** saves the transaction documentation for the client’s application procedure. Step **1824** is followed by the “CONTINUE” step, which returns to routine **1704** shown in FIG. 7.

FIG. 19 is a logic flow diagram illustrating a routine **1704** for logging into the secure medical records maintenance system **1400**. Routine **1704** begins following routine **1702** shown in FIG. 7. In step **1902**, the PRIVALINK software **1402** displays a login screen on the client’s computer **108**. Step **1902** is followed by step **1904**, in which the PRIVALINK software **1402** loads the client’s registration number, which was obtained during the application procedure described with reference to FIG. 18, into the login screen. That is, an accessing client does not have an opportunity to enter the login procedure without having first having gone through the application procedure described with reference to FIG. 18, which requires the applicant to provide a valid registration number, which typically may be a DEA number for a licensed medical practitioner or a registration number issued by the proprietor of the secure medical records maintenance system **1400** for an individual patient.

Step **1904** is followed by step **1906**, in which the PRIVALINK software **1402** receives the practitioner’s personal identification number (PIN). Step **1906** is followed by step **1908**, in which the encryption/decryption module **1408**

validates the practitioner's personal identification number (PIN) for use with the received professional registration number. Step **1908** is followed by step **1910**, in which the encryption/decryption module **1408** determines whether the received PIN is a new PIN for use in connection with the received professional registration number. If the received PIN is a new PIN for use in connection with the received professional registration number, the "YES" branch is followed from step **1910** to step **1912**, in which the PRIVALINK software **1402** prompts the user to complete the interface screens **2000-2007** (FIGS. 20-27). Step **1912** and the "NO" branch from step **1910** are followed by step **1914**, in which the client is linked to the interface screens **2000-2007** (FIGS. 20-27). Step **1914** is followed by "CONTINUE" step, which returns to step **1706** shown on FIG. 7.

User Interface Design

Referring again to FIG. 17, once the medical practitioner has successfully logged into the secure medical records maintenance system **1400**, routine **1704** is followed by step **1706**, in which the PRIVALINK software **1402** displays a "switchboard" user interface **2000**. FIG. 20 is an illustration of a typical "switchboard" user interface **2000**, which includes a number of selection items corresponding to functions available on the server. For example, the interface **2000** typically includes a "provider" selection item **2002** and a "patients" selection item **2004**. To illustrate the operation of the user interface it is assumed that the practitioner initially selects the "provider" selection item **2002**. Thus, step **1706** is followed by step **1708**, in which the PRIVALINK software **1402** receives a "provider" selection from the "switchboard" user interface **2000**.

Step **1708** is followed by step **1710**, in which the PRIVALINK software **1402** displays an "address" user interface **2100**. FIG. 21 is an illustration of a typical "address" user interface **2100**. The interface **2100** includes a first field **2102** for entering a practitioner registration number, typically a DEA number, and a second field **2104** for entering a practitioner-assigned PIN. The interface **2100** also includes a number of other fields for entering the practitioner's contact information, such as address, phone number, and so forth. The interface **2100** also includes an "address" tab **2106**, a "billing info" tab **2108**, and a "cover letter" tab **2110** displayed adjacent to the user interface **2100**. These tabs allow

the user to toggle among corresponding user interface screens. As noted above, the interface **2100** initially appears in a “default” mode with the “address” tab **2106** selected.

For example, step **1710** may be followed by step **1712**, in which the
5 PRIVALINK software **1402** displays a “billing info” user interface in response to user selection of the “billing info” tab **2108**. FIG. 22 is an illustration of a typical “billing information” user interface **2108**. The interface **2100** includes a number of fields **2202** for entering payment authorization, such as a bank credit or debit card number. The interface **2100** may include other types of payment options,
10 such as a bank account number for wire transfers, authorization to include the charges on a telephone bill or Internet service provider bill associated with the Internet link to the server-based system **1400**, reference to an authorized account for billing at a later date, and the like.

Similarly, step **1712** may be followed by step **1714**, in which the
15 PRIVALINK software **1402** displays a “cover letter” user interface **2300** in response to user selection of the “cover letter” tab **2110**. FIG. 23 is an illustration of a typical “cover letter” user interface **2300**. This interface allows the practitioner to author a cover letter that will be included with a health report booklet to be produced by the booklet generation module **1410** on the server-
20 based system **1400**, as described previously with reference to FIG. 9. This allows the practitioner to customize the cover letter for each health report booklet produced by the system.

To further illustrate the operation of the user interface, it is assumed that the practitioner next selects the “patients” selection item **2004** on the “switchboard”
25 user interface **2000**. Thus, step **1714** is followed by step **1716**, in which the PRIVALINK software **1402** displays a “patients” user interface **2400**. FIG. 24 is an illustration of a typical “patient selection” user interface **2400**. This user interface allows the practitioner to select a preexisting patient and to enter new patients. Upon entry of a new patient, the PRIVALINK software **1402** typically
30 assigns a unique patient identification number and a unique medical records identification number to the patient, which the user site “PRIVALINK” software **1402** stores in a correlation table **1600**, as described previously with reference to FIG. 16.

Step **1716** is followed by step **1718**, in which the PRIVALINK software **1402** displays a “patient information” user interface **2500**. FIG. 25 is an illustration of a typical “patient information” user interface **2500**. This interface includes a number of fields for entering patient identification information, such as address, phone number, and so forth. As described previously with reference to FIG. 16, this patient identification information is typically indexed by the patient identification number and stored in the first server **1404**, whereas the patient’s medical data is typically indexed by the patient’s medical records identification number and stored in the second server **1406**.

Step **1718** is followed by step **1720**, in which the PRIVALINK software **1402** displays a “questionnaire data” user interface **2600**. FIG. 26 is an illustration of a typical “questionnaire data” user interface **2600**. This interface includes a number of fields for entering patient diagnostic information, such as family history, age, weight, body fat, and so forth. This interface may be used to enter patient diagnostic information in addition to that received through the meter **10**. For example, the interface **2600** may be used to enter diagnostic data that the meter would have recorded, but the patient failed to enter the data into the meter **10**. In addition, the interface **2600** may be automatically filled in, either partially or completely, with the appropriate data stored on a corresponding smartcard **54**. This step may require entry of an appropriate patient PIN for the corresponding smartcard **54** into the interface **2600**. The interface **2600** may also be used to change or supplement the data read from the smartcard **54**, or enter additional diagnostic information that the meter **10** is not configured to collect from the patient.

Step **1720** is followed by step **1722**, in which the PRIVALINK software **1402** displays a “generate reports” user interface **2700**. FIG. 27 is an illustration of a typical “generate reports” user interface **2700**. This interface includes a number of fields for selecting items to be included in a patient’s health report booklet, such as cover letter, summary, evaluation, and so forth. The interface includes a number of fields for selecting therapy items to be prescribed for the patient and reflected in the patient’s health report booklet, such as lifestyle therapy, lipid drug prescription, blood pressure drug prescription, and so forth.

FIG. 28 is an illustration of typical health report charts generated by the secure medical records maintenance system **1400** for inclusion in a patient's health report booklet. These charts include a "coronary risk factors" chart **2802**, a "personal health consequences" chart **2804**, and an "extended health assessment" chart **2806**. The "coronary risk factors" chart **2802** includes test results and diagnostic information along with ideal ranges and patient goals for these items. The "personal health consequences" chart **2804** includes interpretive data, such as pounds overweight, cardiac age, and stroke risk. The "extended health assessment" chart **2806** includes a projection of future interpretive data, such as a projected comparison of the patient's chronological age and cardiac age.

FIG. 29 is an illustration of an additional health assessment chart **2900** generated by the secure medical records maintenance system **1400** for inclusion in a patient's health report booklet. This chart includes a pictorial representation of cardiac risk factors, such as gender, smoker, personal history, and so forth. The health assessment chart **2900** typically presents a pictorial assessment of the "coronary risk factors" shown in chart **2802**.

Smartcard Payment System

The system described above is largely dependent on the sale of proprietary test strips **28** for the collection of revenue from end users. That is, the meter **10** may be made available to individual patients at little or no cost, with the sale of proprietary test strips **28** providing a major source of revenue for the proprietor of the health monitoring meter. This may be a desirable business model for deploying the meters **10** because it minimizes the initial cost that an individual patient must pay to begin using the device. Having to sell each meter **10** at its full cost, on the other hand, would undermine the economic feasibility of using the meter in many contexts.

Nevertheless, it may also be desirable to provide a meter **10** that does not rely on the sale of proprietary test strips **28** as a major source of revenue. For example, the meter may be adapted to read non-proprietary test strips, or may incorporate a reusable and/or non-invasive testing device, such as an electrode, blood pressure monitoring device, sonic testing device, thermometer, saliva testing device, optical testing device, and the like. Of course, a non-invasive multi-use

testing device may be used many times without affording the proprietor of the health monitoring device an opportunity collect revenue associated with use of the device.

To solve this problem, the smartcard **54** may be utilized as a type of “debit card” for use with the meter **10**. That is, the smartcard **54** device may be purchased with a monetary value, or it may have a monetary value that is replenishable over the Internet using a bank credit or debit card or other conventional payment source. The meter **10** may then deduct the cost of performing particular services (e.g., blood cholesterol test, health assessment, blood sugar test, AIDS test, etc.) from the monetary value represented by the monetary balance stored on the smartcard **54**. In other words, the meter **10** may be configured to activate for the performance of a variety of services upon deducting a charge for the selected service from a monetary value stored on a smartcard **54** inserted into the device.

FIG. 30 is a logic flow diagram illustrating a routine **3000** for adding a monetary value to a smartcard **54** for use with the meter **10**. Although single-use smartcard could be sold with a monetary value for use with the meter **10**, the cost of the smartcards makes a replenishable debit card system preferable. In step **3002**, the proprietor of the system loads an illustrative smartcard **54** with debit card software and a zero balance. Step **3002** is followed by step **3004**, in which the smartcard **54** is issued for use with the meter **10**. For example, the smartcard **54** may be sold or given away separately or in connection with a meter **10**. Typically, the proprietor of the system will be motivated to make both the meter **10** and the smartcard **54** available for little or no acquisition cost because the cost of these assets will be recovered through use of the debit-card type smartcard **54** with meter **10**.

Step **3004** is followed by step **3006**, in which the user establishes a PIN and activates the smartcard **54**, which can be performed using the meter **10** or a conventional computer **108** having a smartcard interface **106**, as shown in FIG. 2. Step **3006** is followed by step **3008**, in which the user loads the smartcard **54** into a conventional computer **108** having a smartcard interface **106** (if the smartcard is not already loaded) and accesses a predefined debit card Internet site by linking to an associated URL, which may be printed on the smartcard or made available

through advertisement, direct mail, or other communication media. The debit card Internet site may typically be a separate server operated by the proprietor of the secure medical records maintenance system **1400**, or it may be another site operated by a separate proprietor. Step **3008** is followed by step **3010**, in which the debit card Internet site establishes communications with the smartcard **54** and verifies the authenticity of the smartcard.

Step **3010** is followed by step **3012**, in which the debit card Internet site prompts the user for input payment authorization for a monetary value to be added to the smartcard. The debit card Internet site may accept a variety of payment options, such as a bank credit or debit card number, a bank account number for wire transfers, authorization to include the charges on a telephone bill or Internet service provider bill associated with the Internet link, and the like. Step **3012** is followed by step **3014**, in which the debit card Internet site validates the received payment authorization. Step **3014** is followed by step **3016**, in which the debit card Internet site determines whether the received authorized payment source has a sufficient balance or credit limit for the monetary value that the user has requested for addition to the smartcard **54**.

If the authorized payment source has a sufficient balance or credit limit, the “YES” branch is followed to step **3018**, in which the debit card Internet site adds the requested monetary value to the smartcard **54** and charges the cost to the authorized payment source. Optionally, the debit card Internet site may also load a rate schedule for services to be provided by the meter **10** onto the smartcard **54**. In this case, the meter **10** reads the monetary value assigned to performance of the test from the smartcard **54**. Thus, rate schedules for various services to be performed by the meter **10** may be changed from time to time, based on quantity discounts or other considerations. Step **3020** is followed by the “END” step, which concludes routine **3000**.

Referring again to step **3016**, if the payment source is invalid or does not have a sufficient balance or credit limit, the “NO” branch is followed to step **3017**, in which the user is prompted for another payment source. If the user enters another payment source, the “YES” branch loops to step **3014** for validation of the alternate payment source. If the user does not enter another

payment source, the “NO” branch is followed by the “END” step, which concludes routine 3000.

FIG. 31 is a logic flow diagram illustrating a routine 3100 for using a smartcard 54 to pay for a service provided by the meter 10. In step 3102, a user obtains a smartcard 54 with a usable monetary value, typically by purchasing a smartcard having a usable monetary value or by adding a monetary value to a smartcard as described with reference to FIG. 30. Step 3102 is followed by step 3104, in which the meter 10 receives a request to provide a service, such as a blood cholesterol test and/or production of a health report, with the associated cost to be charged against a monetary value stored on the smartcard 54. Step 3104 is followed by step 3106, in which the meter 10 prompts the user to enter a smartcard 54 with a usable monetary value into the meter. Step 3106 is followed by step 3108, in which the meter 10 receives the smartcard 54 in the reader 50. It should be understood that at this point the meter is inactive for performing services that require payment. Note also that the meter 10 may be activated alternatively through proprietary test strip validation, as described with reference to FIG. 6.

Step 3110 is followed by step 3112, in which the meter 10 prompts the user for a PIN for use with the smartcard 54 and validates the PIN (i.e., checks the received PIN against the PIN stored on the smartcard). If the PIN is validated, step 3112 is followed by step 3114, in which the meter 10 determines whether the monetary value stored on the smartcard 54 is sufficient to pay for the requested service. As noted previously, the meter 10 may also read the rate schedule from the smartcard 54. Alternatively, if no rate schedule is present on the smartcard 54, the meter 10 may use a predefined “default” rate schedule for the requested service.

If the monetary value stored on the smartcard 54 is sufficient to pay for the requested service, the “YES” branch is followed to step 3116, in which the meter 10 deducts the cost for the requested service from the monetary value stored on the smartcard 54. That is, the meter 10 computes a revised monetary balance equal to the initial monetary value stored on the smartcard less the cost for the requested service, and replaces the initial monetary value stored on the smartcard with the revised monetary balance. Step 3116 is followed by step 3118, in which

the meter 10 activates for performance of the requested service. Step 3118 is followed by step 3120, in which the meter 10 completes the requested service and returns to the inactive state. Referring again to step 3114, if the monetary value stored on the smartcard 54 is insufficient to pay for the requested service, the
5 “NO” branch is followed to step 3122, in which the meter displays a “balance insufficient” message. Steps 3120 and 3122 are followed by the “END” step, which concludes routine 3100.

The present invention thus provides a health monitoring and diagnostic device (LIFESTREAM cholesterol meter) that works in connection with a
10 network-based comprehensive health analysis and reporting system. The invention also provides a secure medical records maintenance system that works in conjunction with the health monitoring and diagnostic device, and alternative business models for deploying the health monitoring and diagnostic devices. It should be understood that the foregoing pertains only to the preferred
15 embodiments of the present invention, and that numerous changes may be made to the embodiments described herein without departing from the spirit and scope of the invention.

CLAIMS

The invention claimed is:

- 5 1. A test medium for use with a health monitoring device, comprising:
a test strip for receiving and carrying a sample of biological fluid or tissue;
the test strip carrying the sample of biological fluid or tissue being readable
by the health monitoring device to obtain test results based on the sample of
biological tissue or fluid and calibration data specific to the test strip;
10 the test strip corresponding to a memory device readable by the health
monitoring device, the memory device storing a code number and the calibration
data;
the test strip having an associated test strip identification number that is
mathematically derivable from the code number;
15 the test strip identification number being receivable by the health
monitoring device through a user input device; and
the test strip identification number being usable to activate the health
monitoring device for use with the test strip only if the received test strip
identification number corresponds to a test strip identification number
20 mathematically derived by the health monitoring device based on the code
number read from the memory device.
2. The test medium of claim 1, wherein:
the memory device also stores an expiration date for the test strip;
25 the expiration date is readable by the health monitoring device; and
the expiration date is usable to activate the health monitoring device for use
with the test strip only if the expiration date is prior to a current date read by the
health monitoring device from an internal clock.
- 30 3. The test medium of claim 1, wherein the memory device comprises a
romkey insertable into a socket housed within the health monitoring device.

4. The test medium of claim 1, wherein:

a plurality of the test strips are packaged for distribution together with the romkey; and

5 the test strip identification number is printed on the test strips, printed on packaging for the test strips, or printed on a tag packaged with the test strip.

5. A hand-held health monitoring device, comprising:

an enclosure for housing a disposable test strip for use with the health monitoring device;

10 a holder for removably supporting a device for gathering a sample of biological fluid or tissue;

a test strip reader operable for reading the test strip carrying the sample of biological fluid or tissue and obtaining test results based on the sample of biological tissue or fluid and calibration data specific to the test strip;

15 a memory reading device functionally connected to the test strip reader and operable for reading the calibration data from a memory device;

a user input device operable for receiving user input commands;

a display device operable for displaying information;

20 a processor functionally connected to the test strip reader, the user input device, and the display device, the processor containing a program module operable for obtaining the test results from the test strip reader and causing the display device to display the test results; and

a data drive functionally connected to the processor and operable for writing the test results to a removable memory storage device.

25

6. The health monitoring device of claim 5, wherein the processor is further operable for:

determining whether a personal identification number has been previously stored on the removable memory storage device;

5 if the personal identification number has not been previously stored on the removable memory storage device prompting the user to enter a personal identification number, storing the received personal identification number on the removable memory storage device; and

10 if the personal identification number has been previously stored on the removable memory storage device, prompting the user to enter a personal identification number, comparing the stored personal identification number to the received personal identification number, and writing the test results to the removable memory storage device only if the stored personal identification number corresponds to the received personal identification number.

15

7. The hand-held health monitoring device of claim 5, further comprising:

a clam-shell case openable to reveal first and second compartments;

20 the first compartment containing the enclosure for housing the disposable test strip and the holder for removably supporting the biological fluid or tissue gathering device; and

the second compartment containing the test strip reader, the memory reading device, the display device, the processor, and the data drive.

25 8. The hand-held health monitoring device of claim 5, wherein:
the biological fluid or tissue includes a droplet of human blood; and
the test results include total cholesterol levels.

9. The hand-held health monitoring device of claim 8, wherein the test strip reader is operable for reading a second type of test strip carrying a second sample of biological fluid or tissue and obtaining health-related test results based on the second sample of biological tissue or fluid and calibration data specific to the second type of test strip, further comprising:

a second memory reading device functionally connected to the test strip reader and operable for reading calibration data from a second memory device corresponding to the second type of test strip.

10. The hand-held health monitoring device of claim 9, wherein:
the second biological fluid or tissue includes a droplet of human blood; and
the test results include blood glucose levels.

11. A health monitoring device, comprising:

a test strip reader operable for reading the test strip carrying a sample of biological fluid or tissue and obtaining health-related test results based on the sample of biological tissue or fluid and calibration data specific to the test strip;

a memory reading device functionally connected to the test strip reader and operable for reading the calibration data from a memory device;

a user input device operational for receiving user input commands;

a display device operable for displaying information; and

a processor functionally connected to the test strip reader, the user input device, and the display device, the processor containing a diagnostic program module operable for:

obtaining the test results from the test strip reader,
causing the display device to prompt the user to enter diagnostic information using the user input device,

performing a diagnostic analysis to produce diagnostic results based on the test results and the diagnostic information, and

causing the display device to display the diagnostic results.

12. The health monitoring device of claim 11, wherein:

the biological fluid or tissue includes a droplet of human blood;

the test results include blood lipid levels;

the diagnostic information includes one or more of the following data items

5 corresponding to a person providing the droplet of human blood:

gender, ethnicity, family history of heart disease, personal history of heart disease, personal history of diabetes, personal history of smoking, height, weight, age, blood pressure, fitness level; and

the diagnostic results include one or more of the following data items

10 corresponding to the person providing the droplet of human blood:

a medical risk index, a recommended weight loss, a five-year risk of heart attack, a ten-year risk of heart attack, a cardiac age, an extended age, and a risk of stroke.

15 13. The health monitoring device of claim 12, further comprising a data drive functionally connected to the processor and operable for writing the test results, the diagnostic information, and the diagnostic results to a removable memory storage device.

14. The health monitoring device of claim 13, wherein the processor is further operable for:

determining whether a personal identification number has been previously stored on the removable memory storage device;

5 if the personal identification number has not been previously stored on the removable memory storage device prompting the user to enter a personal identification number, storing the received personal identification number on the removable memory storage device; and

10 if the personal identification number has been previously stored on the removable memory storage device, prompting the user to enter a personal identification number, comparing the stored personal identification number to the received personal identification number, and writing the test results, the diagnostic information, and the diagnostic results to the removable memory storage device only if the stored personal identification number corresponds to the received
15 personal identification number.

15. The health monitoring device of claim 14, wherein the test strip reader is operable for reading a second type of test strip carrying a second sample of biological fluid or tissue and obtaining health-related test results based on the
20 second sample of biological tissue or fluid and calibration data specific to the second type of test strip, further comprising:

a second memory reading device functionally connected to the test strip reader and operable for reading calibration data from a second memory device corresponding to the second type of test strip.

25

16. The health monitoring device of claim 15, wherein:
the second biological fluid or tissue includes a droplet of human blood; and
the test results include blood glucose levels.

17. A health monitoring device for use with a disposable test strip, comprising:

a test strip reader operable for reading the test strip carrying a sample of biological fluid or tissue to obtain test results based on the sample of biological tissue or fluid and calibration data specific to the test strip;

a memory device readable by the health monitoring device and storing a code number and the calibration data specific to the test strip;

a user input device operable for entering a test strip identification number into the health monitoring device; and

a processor functionally connected to the test strip reader, the memory device, and the user input device, the processor operable for:

reading the code number from the memory device,

mathematically deriving a test strip identification number corresponding to the code number,

comparing the received test strip identification number to the derived test strip identification number, and

activating the health monitoring device for use with the test strip only if the received test strip identification number corresponds to the derived test strip identification number.

18. The health monitoring device of claim 17, wherein:

the health monitoring device also includes a clock defining a current date;

the memory device also stores an expiration date for the test strip;

the processor is operative to read the expiration date and the current date;

and

the processor is operative to activate the health monitoring device for use with the test strip only if the expiration date is prior to the current date.

19. The health monitoring device of claim 18, wherein the processor is further operative for:

receiving an activation code through the user input device;

5 computing an activation code based on the current date and instructions contained in an activation routine stored within the health monitoring device; and
activating the health monitoring device only if the computed activation code corresponds to the received activation code.

10 20. The health monitoring device of claim 18, wherein the processor is further operable for:

obtaining the test results from the test strip reader,

causing the display device to prompt the user to enter diagnostic information using the user input device,

15 performing a diagnostic analysis to produce diagnostic results based on the test results and the diagnostic information; and
causing the display device to display the diagnostic results.

21. The health monitoring device of claim 20, wherein:

the biological fluid or tissue includes a droplet of human blood;

20 the test results include blood lipid levels;

the diagnostic information includes one or more of the following data items corresponding to a person providing the droplet of human blood:

25 gender, ethnicity, family history of heart disease, personal history of heart disease, personal history of diabetes, personal history of smoking, height, weight, age, blood pressure, fitness level; and

the diagnostic results include one or more of the following data items corresponding to the person providing the droplet of human blood:

30 a medical risk index, a recommended weight loss, a five-year risk of heart attack, a ten-year risk of heart attack, a cardiac age, an extended age, and a risk of stroke.

22. The health monitoring device of claim 21, further comprising a data drive functionally connected to the processor and operable for writing the test results, the diagnostic information, and the diagnostic results to a removable memory storage device.

5

23. The health monitoring device of claim 22, wherein the processor is further operable for:

determining whether a personal identification number has been previously stored on the removable memory storage device;

10 if the personal identification number has not been previously stored on the removable memory storage device prompting the user to enter a personal identification number, storing the received personal identification number on the removable memory storage device; and

15 if the personal identification number has been previously stored on the removable memory storage device, prompting the user to enter a personal identification number, comparing the stored personal identification number to the received personal identification number, and writing the test results, the diagnostic information, and the diagnostic results to the removable memory storage device only if the stored personal identification number corresponds to the received
20 personal identification number.

24. The health monitoring device of claim 23, wherein the test strip reader is operable for reading a second type of test strip carrying a second sample of biological fluid or tissue and obtaining health-related test results based on the
25 second sample of biological tissue or fluid and calibration data specific to the second type of test strip, further comprising:

a second memory reading device functionally connected to the test strip reader and operable for reading calibration data from a second memory device corresponding to the second type of test strip.

30

25. The health monitoring device of claim 24, wherein:
the second biological fluid or tissue includes a droplet of human blood; and
the test results include blood glucose levels.

5 26. A system for remotely producing health reports, comprising:
a health monitoring device including:

a test strip reader operable for reading a test strip carrying a sample
of biological fluid or tissue and obtaining test results based on the sample of
biological tissue or fluid and calibration data specific to the test strip,

10 a memory reading device functionally connected to the test strip
reader and operable for reading the calibration data from a memory device, and

a data drive functionally connected to the test strip reader and
operable for writing the test results to a memory storage device;

a health report server operable for:

15 receiving the test results and additional diagnostic information,
compiling a health report based on the test results and the additional
diagnostic information, and

transmitting the health report; and

a computer station operable for:

20 reading the test results from the memory storage device,
establishing a network connection with the health report server,
receiving the additional diagnostic information,
transmitting the test results and the additional diagnostic information
to the health report server,

25 receiving the health report from the health report server, and
printing the health report.

27. The system for remotely producing health reports of claim 26,
wherein the health report includes a trend analysis including test results compiled
30 for a number of samples of the biological fluid or tissue.

28. The system for remotely producing health reports of claim 26, wherein:

the biological fluid or tissue includes a droplet of human blood;

the test results include blood lipid levels;

5 the additional diagnostic information includes one or more of the following data items corresponding to a person providing the droplet of human blood: gender, ethnicity, family history of heart disease, personal history of heart disease, personal history of diabetes, personal history of smoking, height, weight, age, blood pressure, fitness level.

10 29. The system for remotely producing health reports of claim 26, wherein:

the diagnostic results include one or more of the following data items: a medical risk index, a recommended weight loss, a five-year risk of heart attack, a
15 ten-year risk of heart attack, a cardiac age, an extended age, and a risk of stroke.

30. The system for remotely producing health reports of claim 29, wherein:

the additional diagnostic information includes a newly-prescribed drug and
20 other currently-prescribed drugs; and

the health report includes a data sheet for the newly-prescribed drug and information relating to cross-reactions between the newly-prescribed drug and the other currently-prescribed drugs.

25 31. The system for remotely producing health reports of claim 30, wherein the health report further includes one or more of the following items:

target weight and blood lipid levels;

a schedule for future testing using the health monitoring device,

a health assessment summary;

30 a coronary risk assessment; and

dietary guidelines to lower cholesterol.

32. The system for remotely producing health reports of claim 31, wherein the test strip reader is operable for reading a second type of test strip carrying a second sample of biological fluid or tissue and obtaining health-related test results based on the second sample of biological tissue or fluid and calibration data specific to the second type of test strip, further comprising:

a second memory reading device functionally connected to the test strip reader and operable for reading calibration data from a second memory device corresponding to the second type of test strip.

33. The system for remotely producing health reports of claim 32, wherein:

the test results include blood glucose levels; and

the health report further includes a trend analysis including test results compiled for a number of sample of the biological fluid or tissue.

34. A secure medical records maintenance system, comprising:

a plurality of removable memory storage devices, each operable for storing medical data for an associated patient, a patient-specified personal identification number, and a medical records identification number secured by the patient-specified personal identification number;

a first remote server operable for storing patient identification information indexed patient identification numbers;

a second remote server operable for storing patient medical data indexed by the medical records identification numbers; and

the medical data maintained in the second remote server cannot be correlated to the associated patient identification information maintained in the first remote server based on the information contained in the first and second remote servers.

35. The secure medical records maintenance system of claim 34, wherein each removable memory storage device also stores a patient identification number corresponding to the medical records identification number stored on the removable memory storage.

36. The secure medical records maintenance system of claim 34, further comprising a correlation table uniquely associating each medical records identification number with a particular one of the patient identification numbers.

5 37. The secure medical records maintenance system of claim 36, wherein:
the correlation table resides on a practitioner computer associated with a
licensed medical practitioner having an assigned professional registration number;
and

10 the first and second remote servers can be accessed by the practitioner
computer through encrypted communications secured by an application
procedure comprising validation of the practitioner's registration number.

38. The secure medical records maintenance system of claim 37, wherein
the application procedure further comprises receipt and validation of a client-
15 supplied personal identification number.

39. The secure medical records maintenance system of claim 38, wherein
the application procedure comprises issuance of a client certificate insuring that
access to the first and second remote servers occurs from the practitioner's
20 computer.

40. The secure medical records maintenance system of claim 34, wherein
access is granted to the first remote server, but not to the second server, for the
purpose of generating a mailing list of patients without divulging any medical data
25 associated with the patients.

41. The secure medical records maintenance system of claim 34, wherein
access is granted to the second remote server, but not to the first server, for the
purpose of conducting investigative analyses involving patient medical data
30 without divulging any patient identification information associated with the
patients.

42. The secure medical records maintenance system of claim 34, wherein the medical data stored on each removable memory storage device is automatically erased from the memory storage device after the data is entered into the second remote server.

43. The secure medical records maintenance system of claim 34, wherein the removable memory storage device is receivable within a hand-held health monitoring device operable for storing the medical data on the removable memory storage device.

44. The secure medical records maintenance system of claim 34, wherein the removable memory storage device is receivable within a computer operable for transmitting the medical data to the second remote server over the Internet.

45. A health monitoring device operable for:
obtaining medical data associated with a patient;
reading an initial monetary balance stored on a removable memory storage device;

determining whether the initial monetary balance is sufficient to pay a monetary value assigned to performance of a test involving the medical data to be performed by the testing device; and

if the initial monetary balance is sufficient to pay for the test,
computing a revised monetary balance by deducting the monetary value assigned to performance of the test from the initial monetary balance,
replacing the initial monetary balance with the revised monetary balance on the removable memory storage device, and
activating the health monitoring device for performance of the test.

46. The health monitoring device of claim 45, further operable for reading the monetary value assigned to performance of the test from the removable memory storage device.

47. The health monitoring device of claim 45, further operable for prompting a user of the health monitoring device to replenish the monetary value stored on the removable memory storage device if the initial monetary balance is in sufficient to pay the monetary value assigned to performance of the test.

5

48. A system operable for performing diagnostic health tests, comprising:
a removable memory storage device;

an network-based server operable for remotely charging a cost to a payment source and crediting the cost to an initial balance stored on the removable memory storage device; and

10

a health monitoring device, operable for:

obtaining medical data associated with a patient;

reading an initial monetary balance stored on a removable memory storage device;

15

determining whether the initial monetary balance is sufficient to pay a monetary value assigned to performance of a test involving the medical data to be performed by the testing device; and

if the initial monetary balance is sufficient to pay for the test,

computing a revised monetary balance by deducting the monetary value assigned to performance of the test from the initial monetary balance,

20

replacing the initial monetary balance with the revised monetary balance on the removable memory storage device, and

25



activating the health monitoring device for performance of the test.

49. The system of claim 48, wherein:

network-based server operable for remotely storing the monetary value assigned to performance of the test on the removable memory storage device; and

30

the health monitoring device is further operable for reading the monetary value assigned to performance of the test from the removable memory storage device.

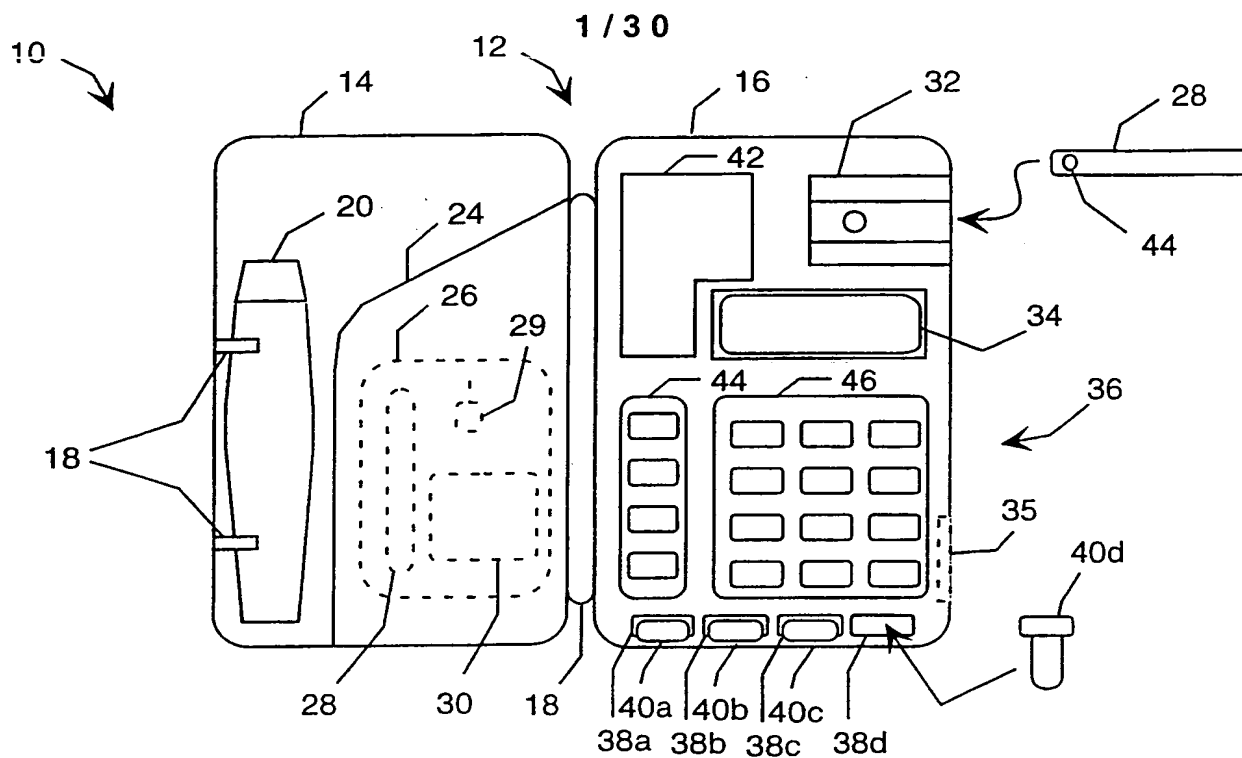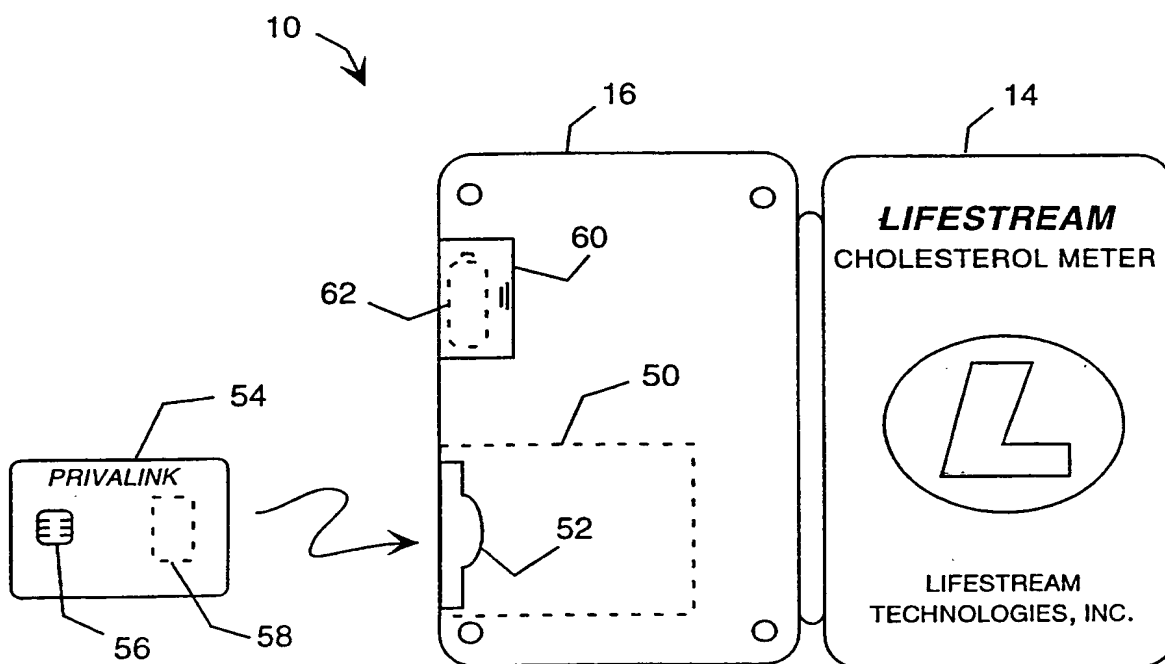**FIG. 1A****FIG. 1B**

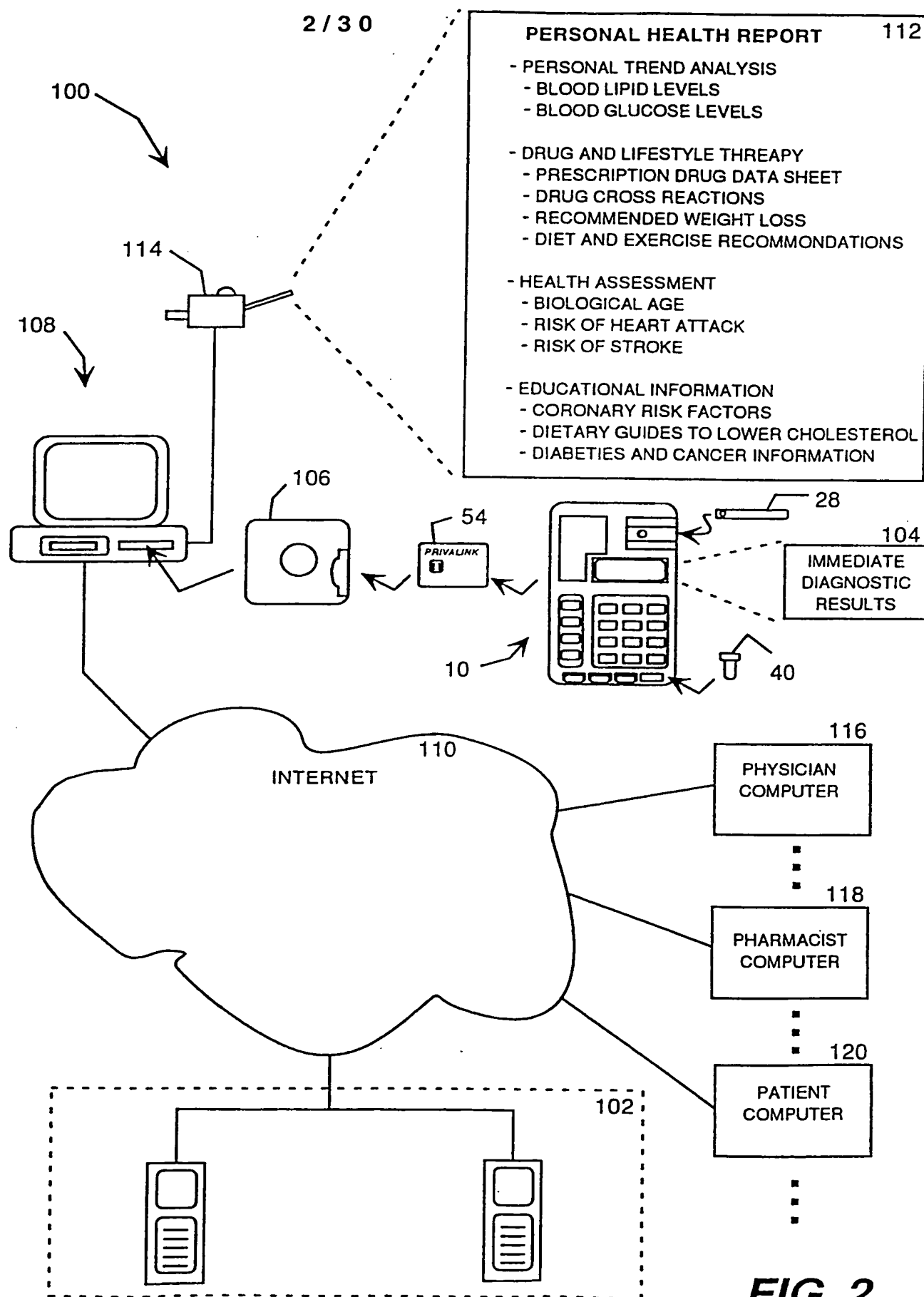
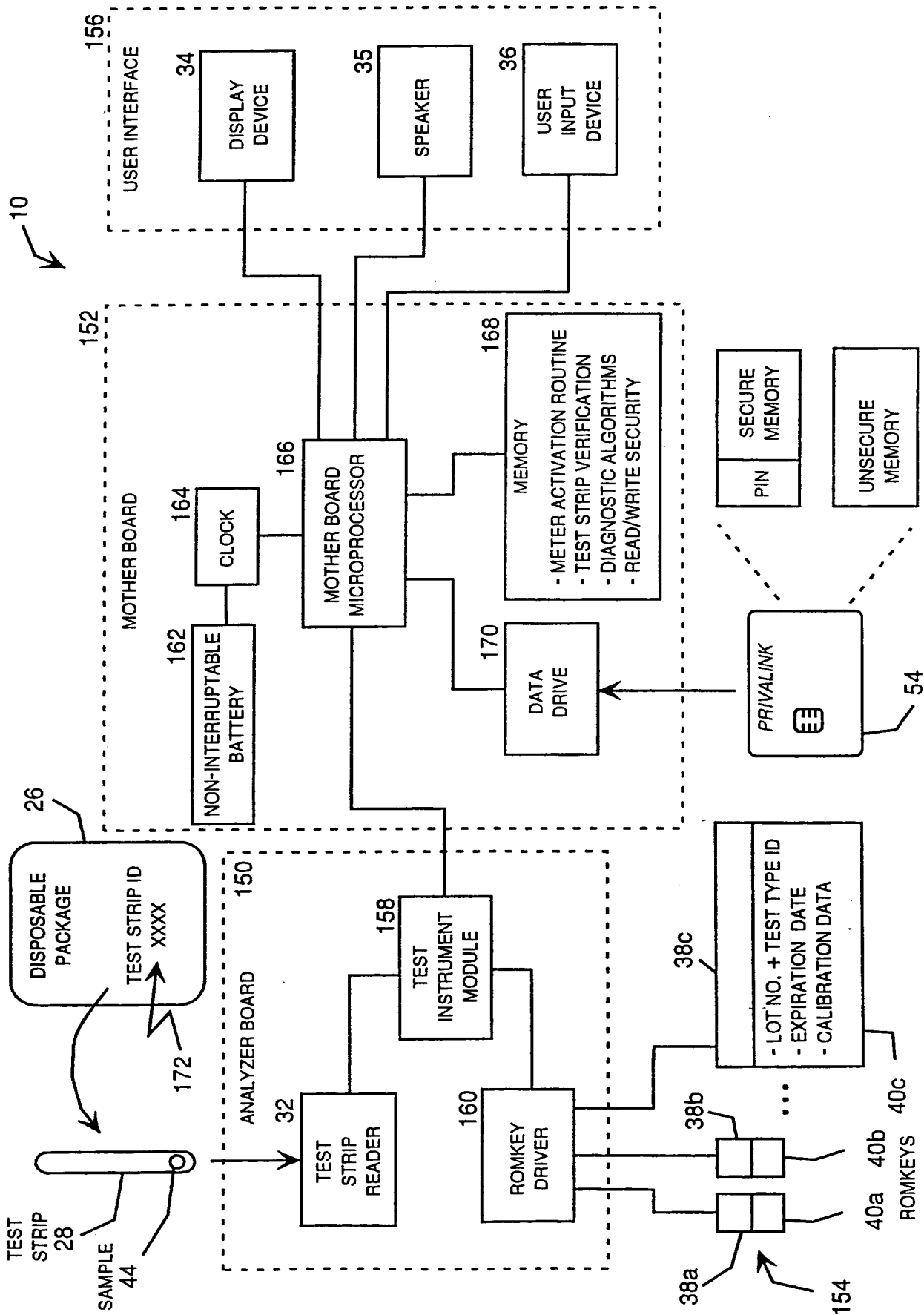
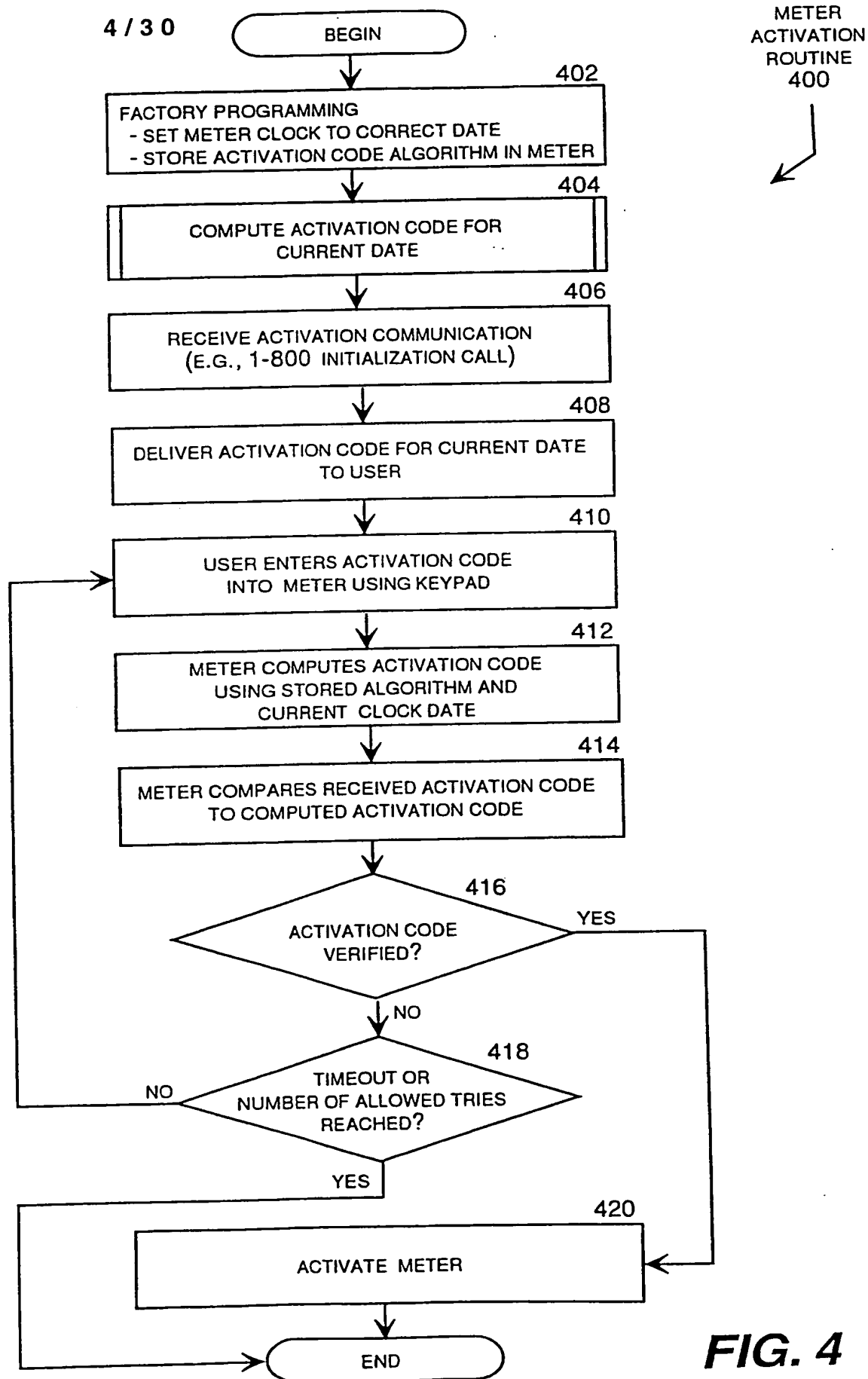


FIG. 3



**FIG. 4**

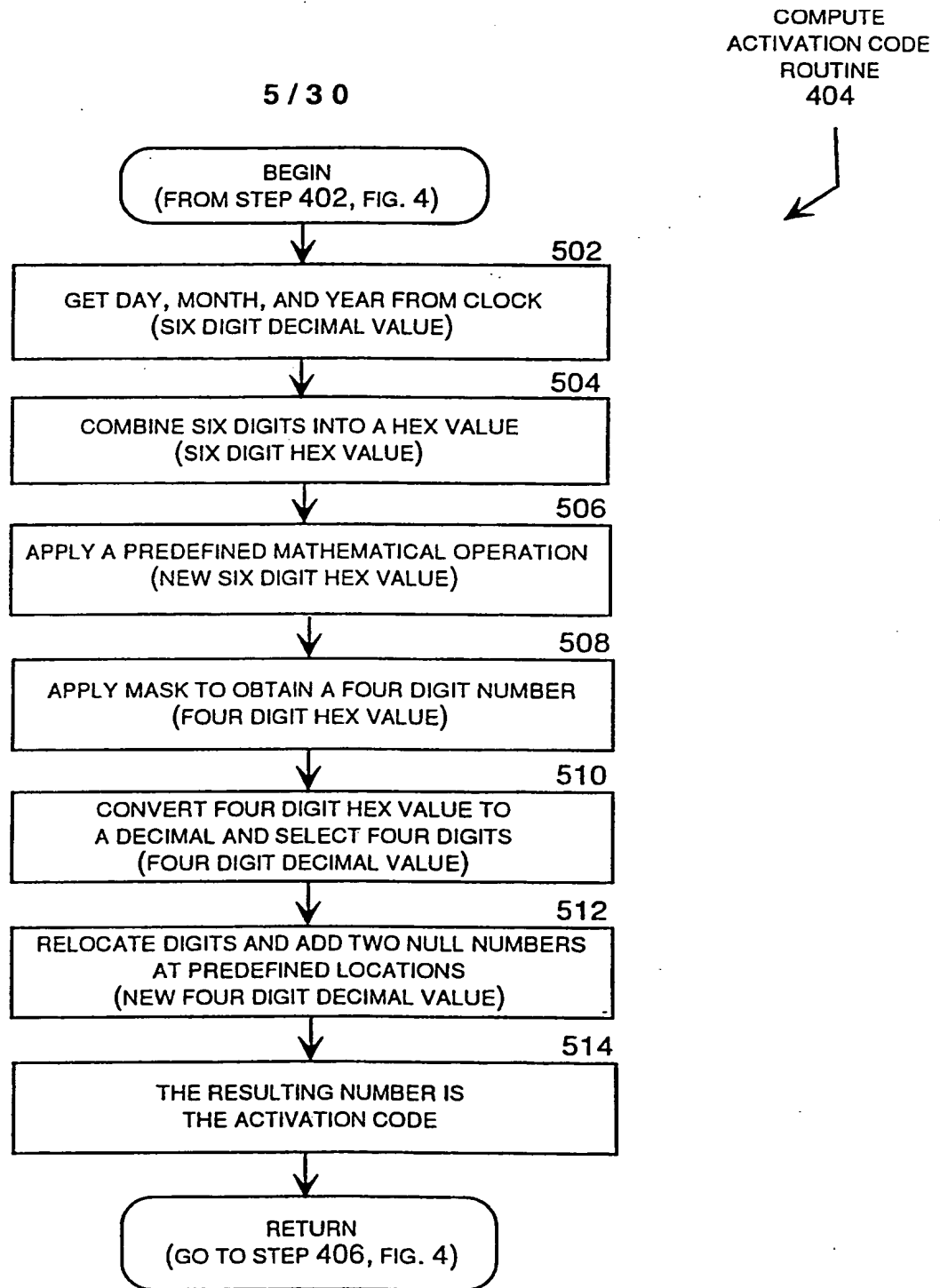
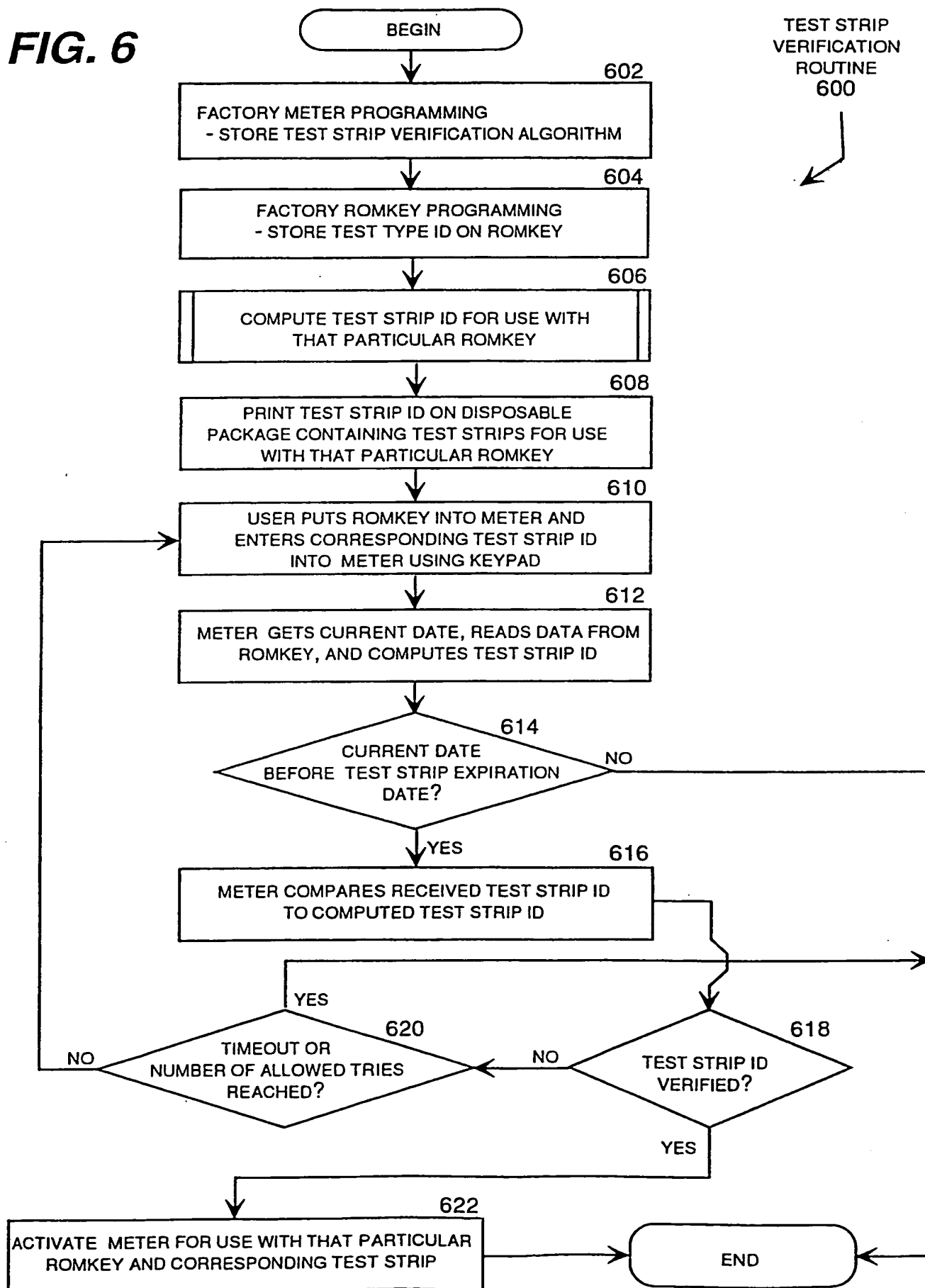
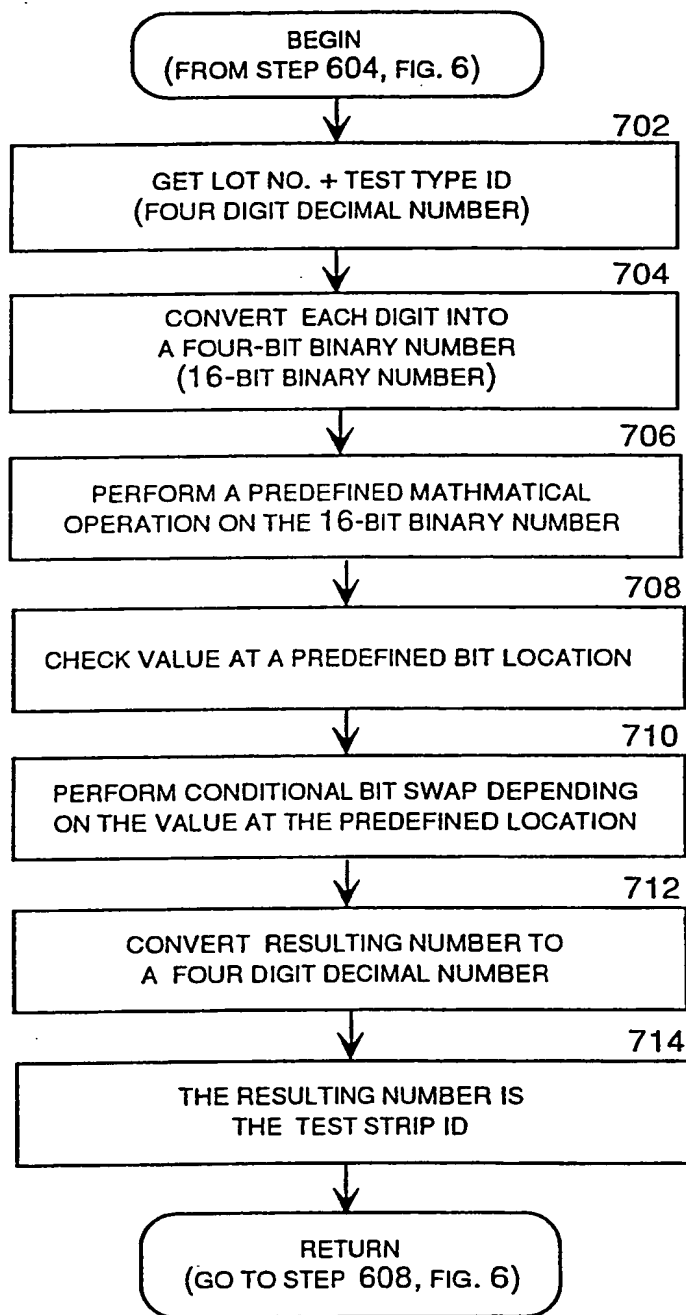
**FIG. 5**

FIG. 6

7 / 3 0

COMPUTE
TEST STRIP ID
ROUTINE
606**FIG. 7**

8 / 3 0



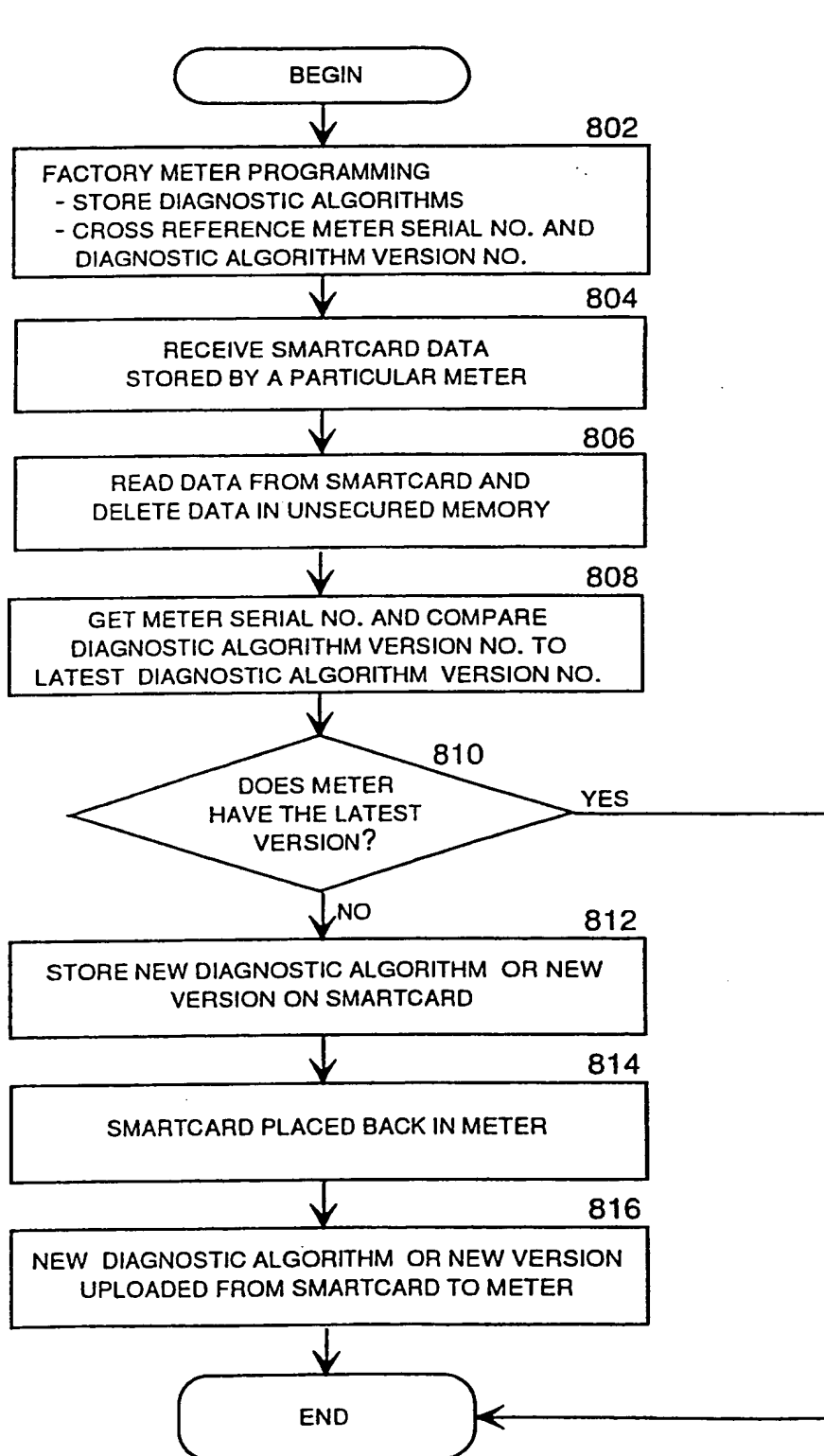
DIAGNOSTIC
PROGRAMMING
ROUTINE
800

FIG. 8


**FIG. 9**

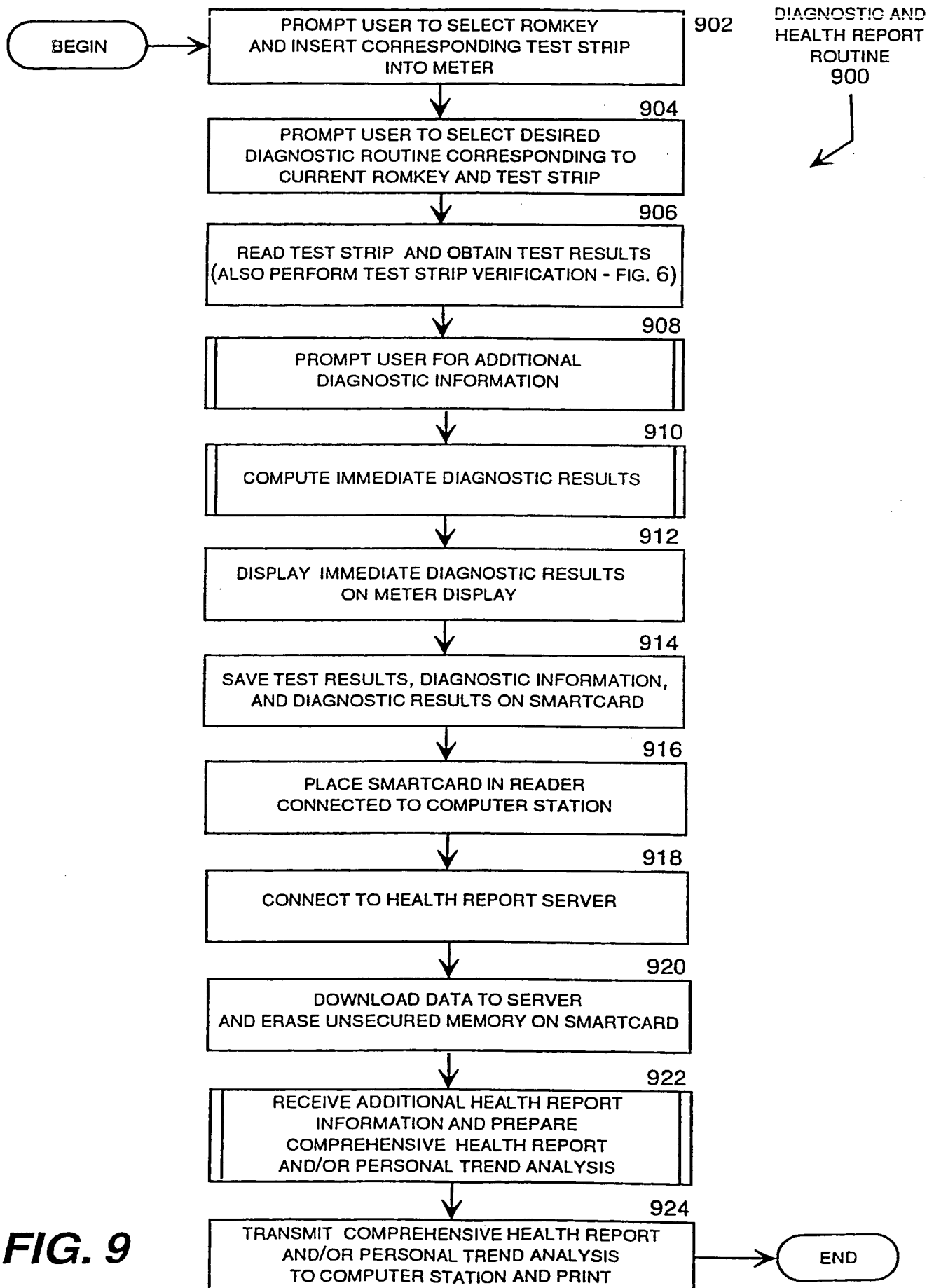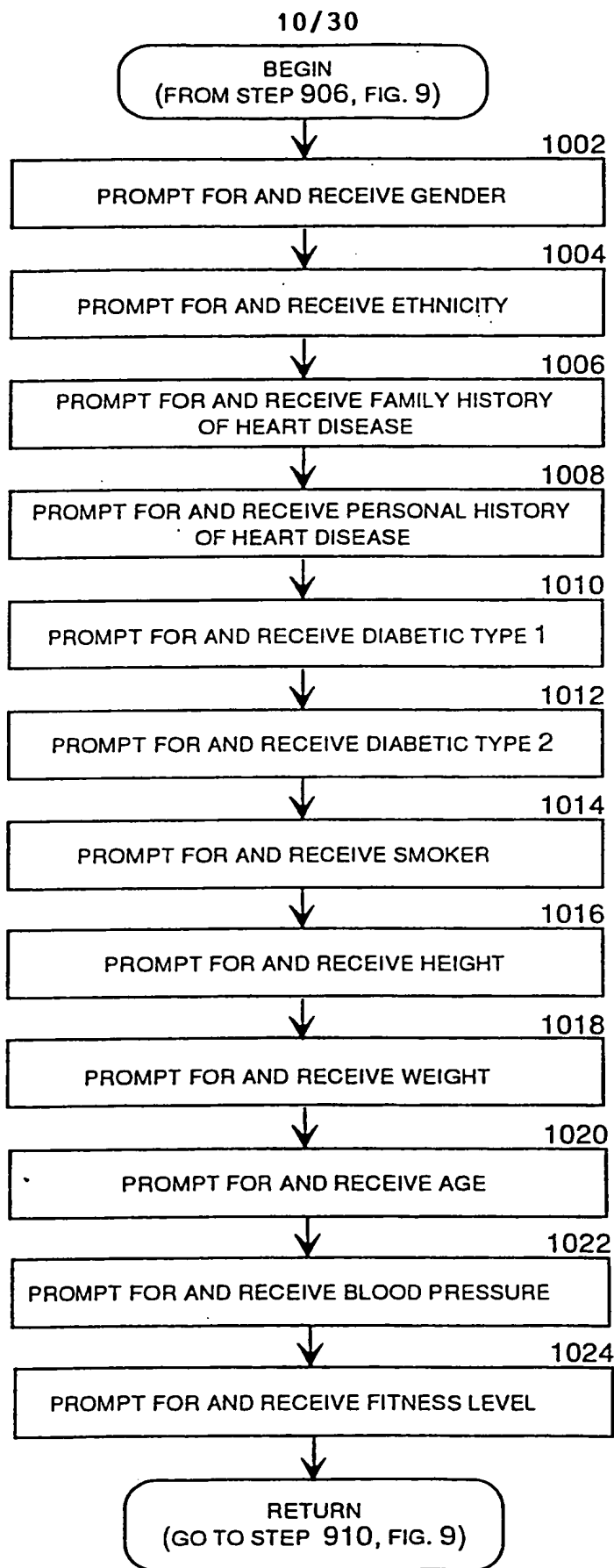

EXAMPLE
CHOLESTEROL
ROUTINE
908

FIG. 10

11/30

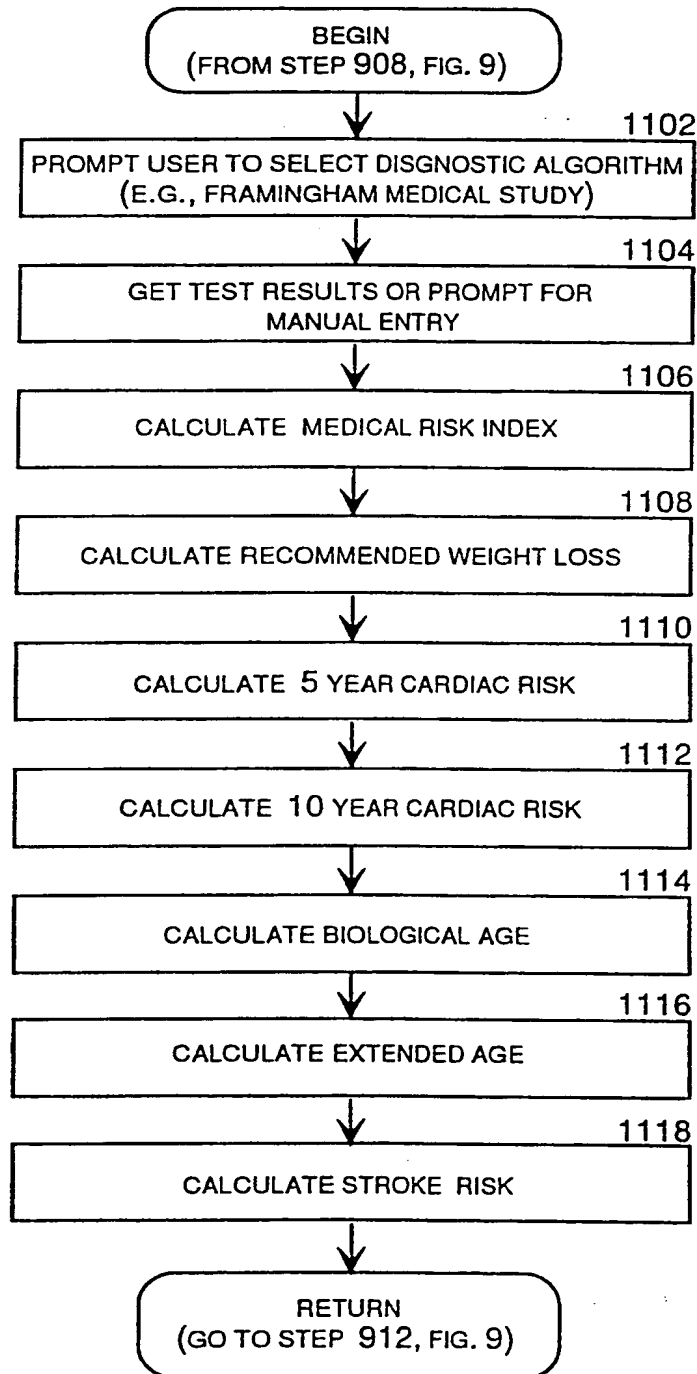EXAMPLE
CHOLESTEROL
ROUTINE
910**FIG. 11**

12/30

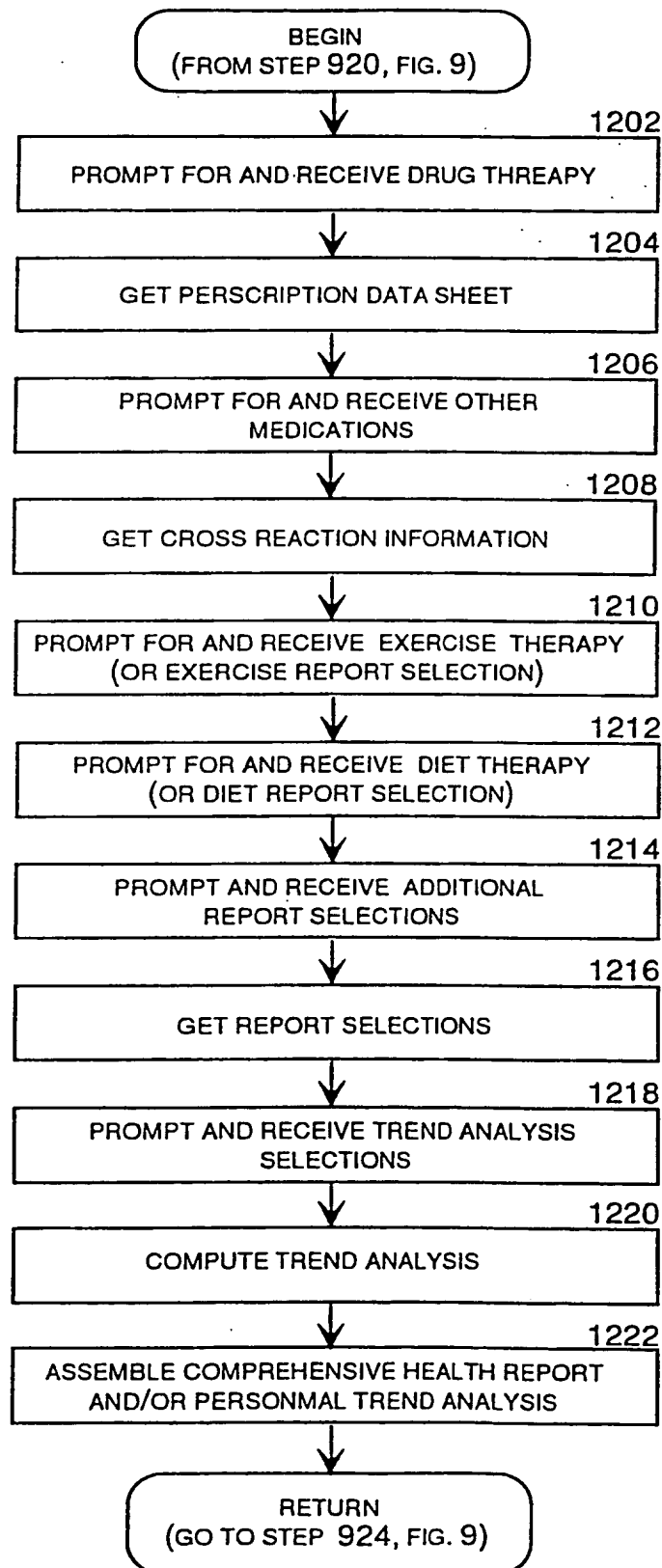
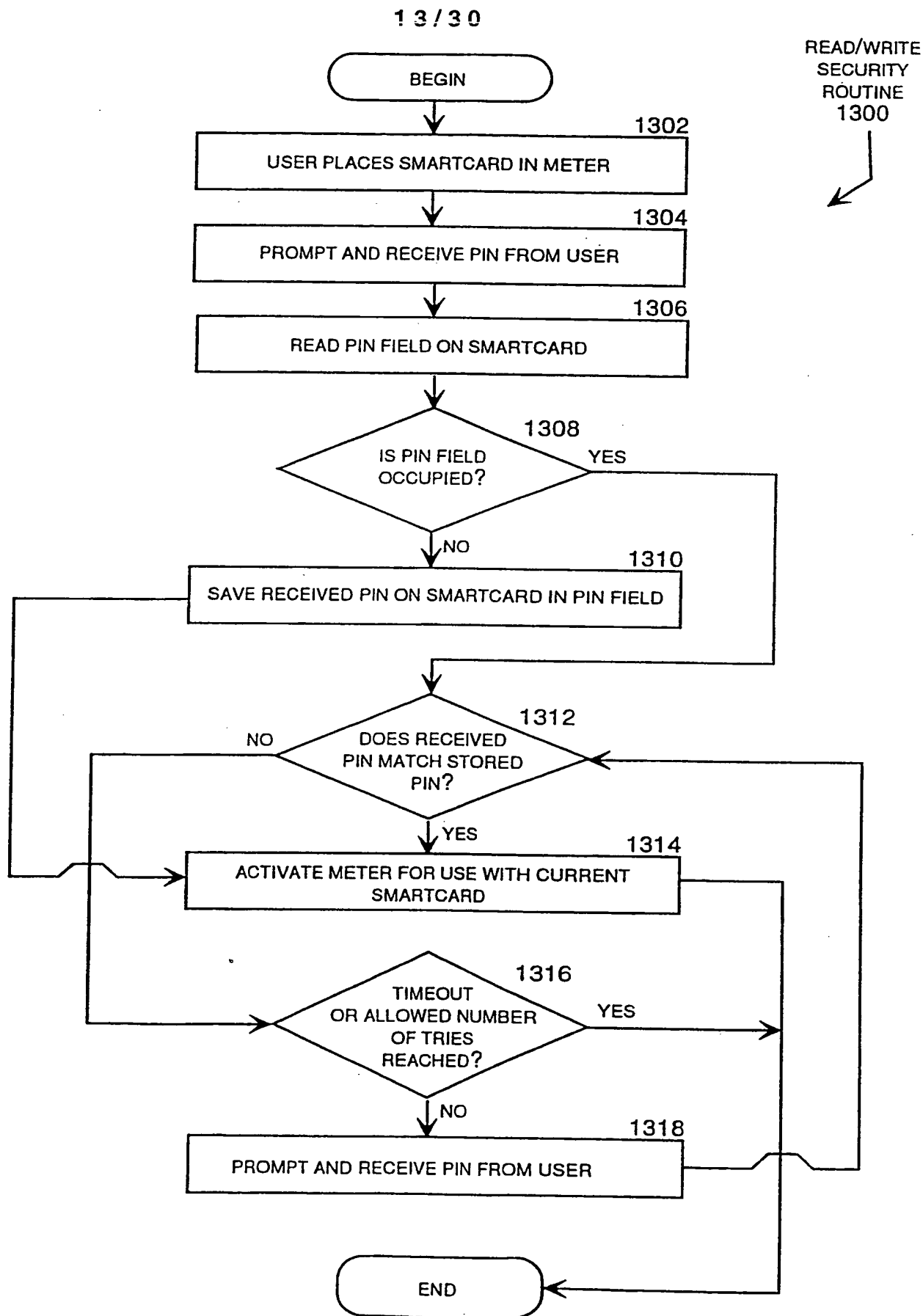
EXAMPLE
CHOLESTEROL
ROUTINE
922

FIG. 12

**FIG. 13**

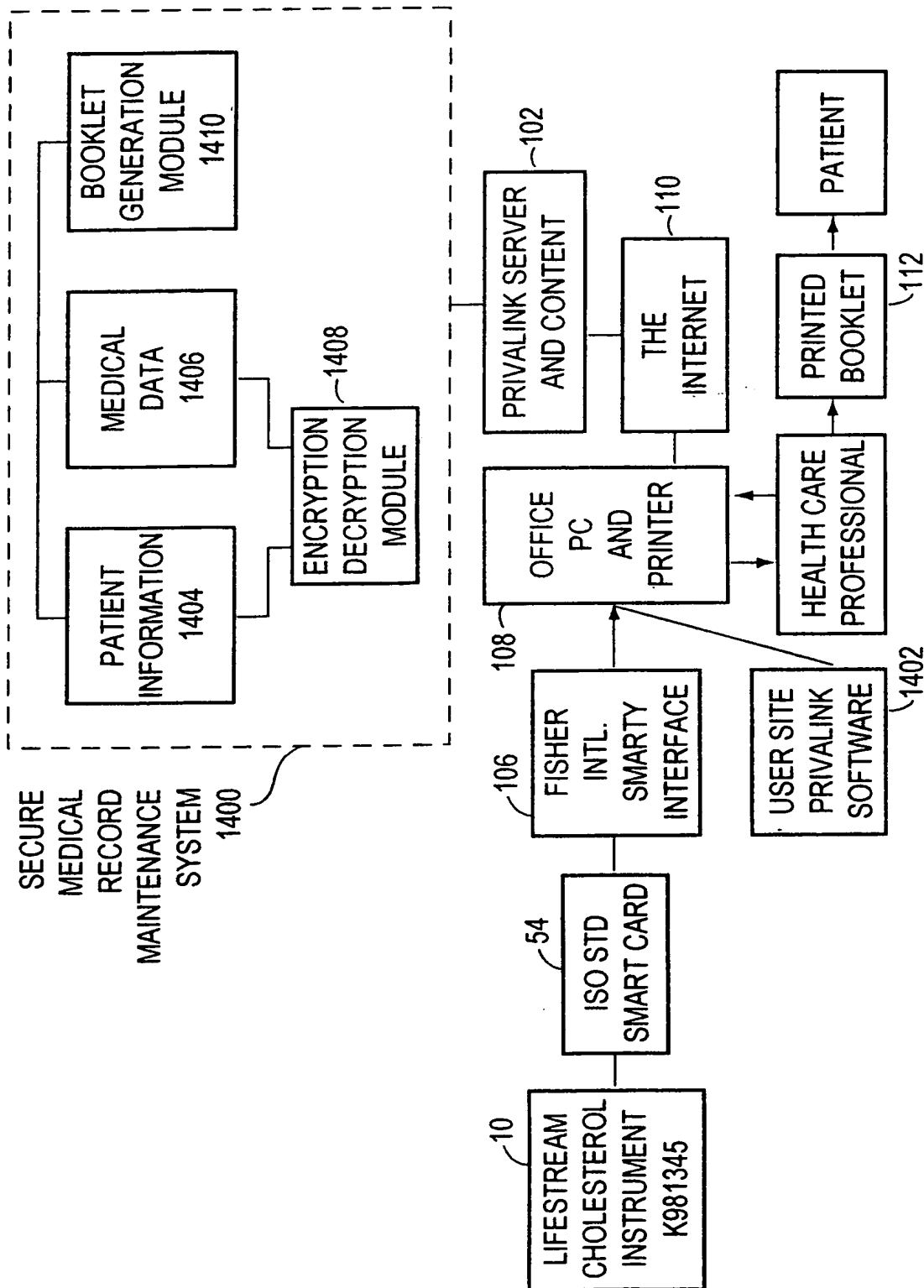


FIG. 14

15/30

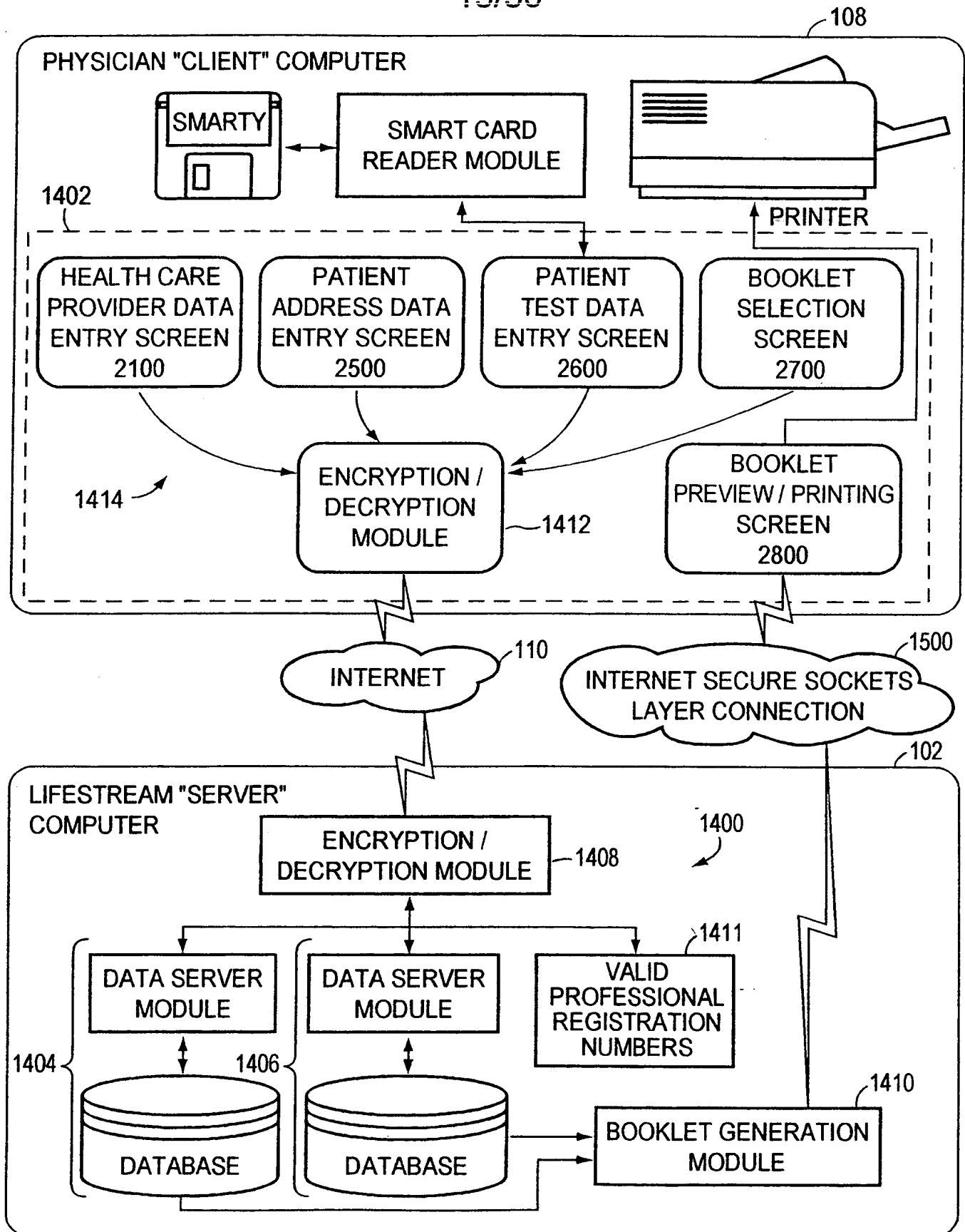


FIG. 15

16/30

GLOBAL SYSTEM HARDWARE PICTORIAL


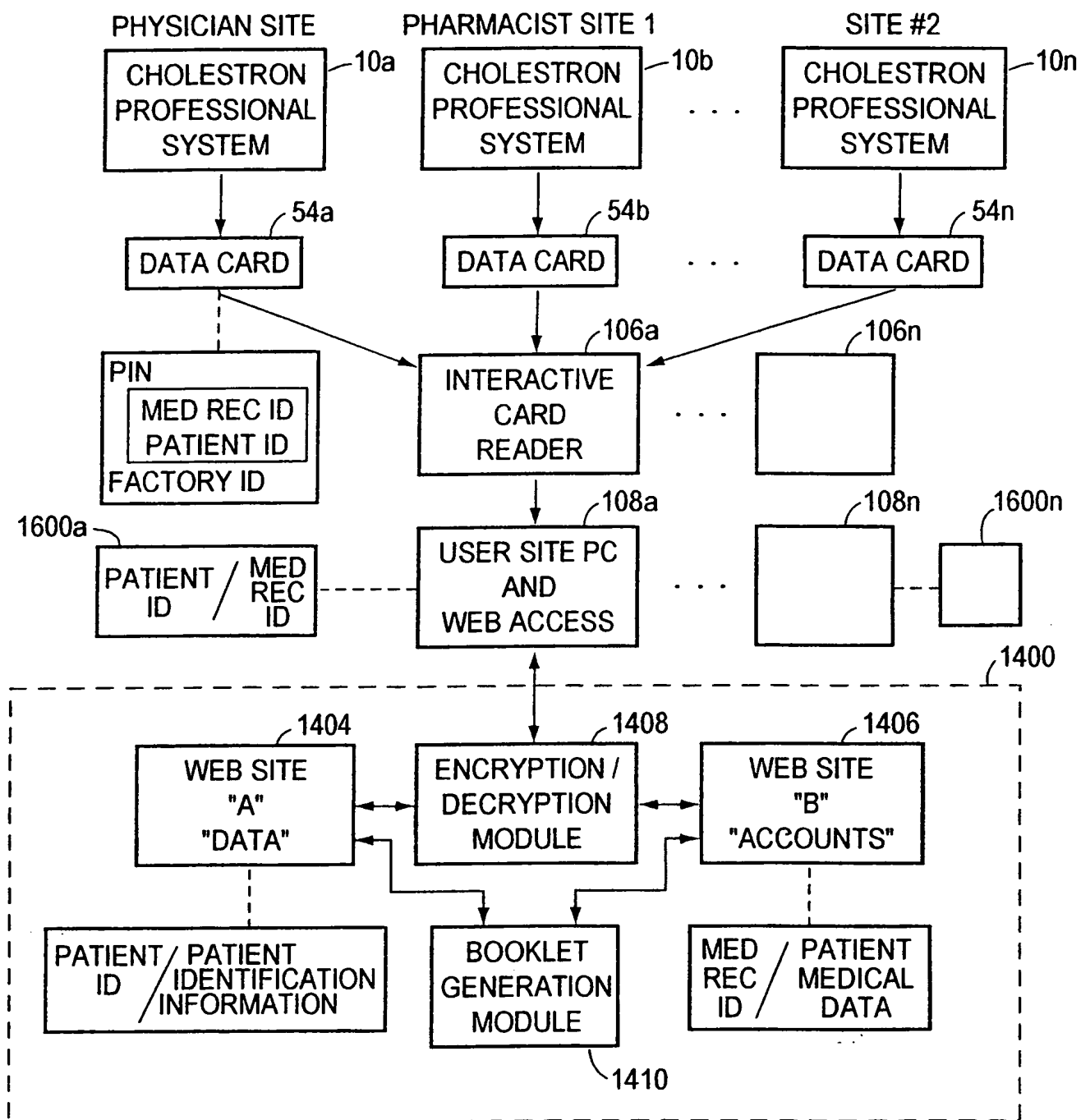


FIG. 16


17/30

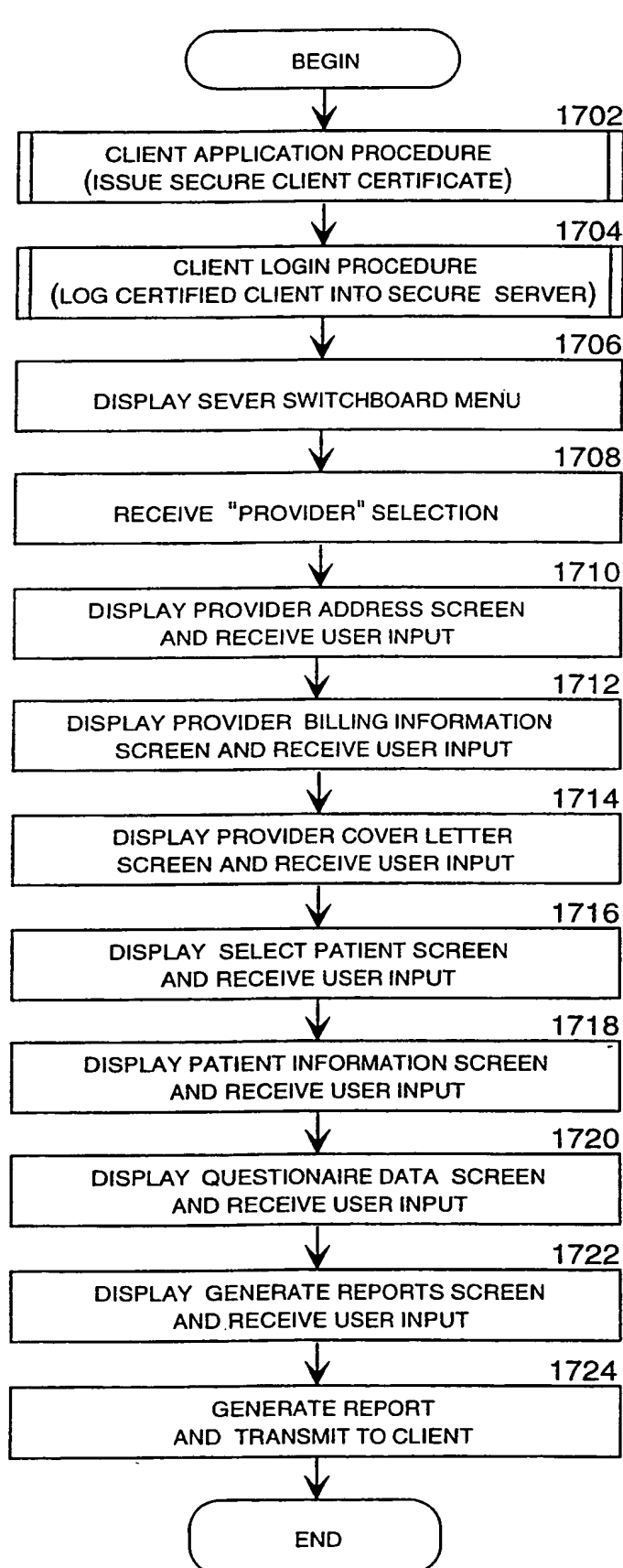


CLIENT
ACCESS
ROUTINE
1700

FIG. 17

18/30

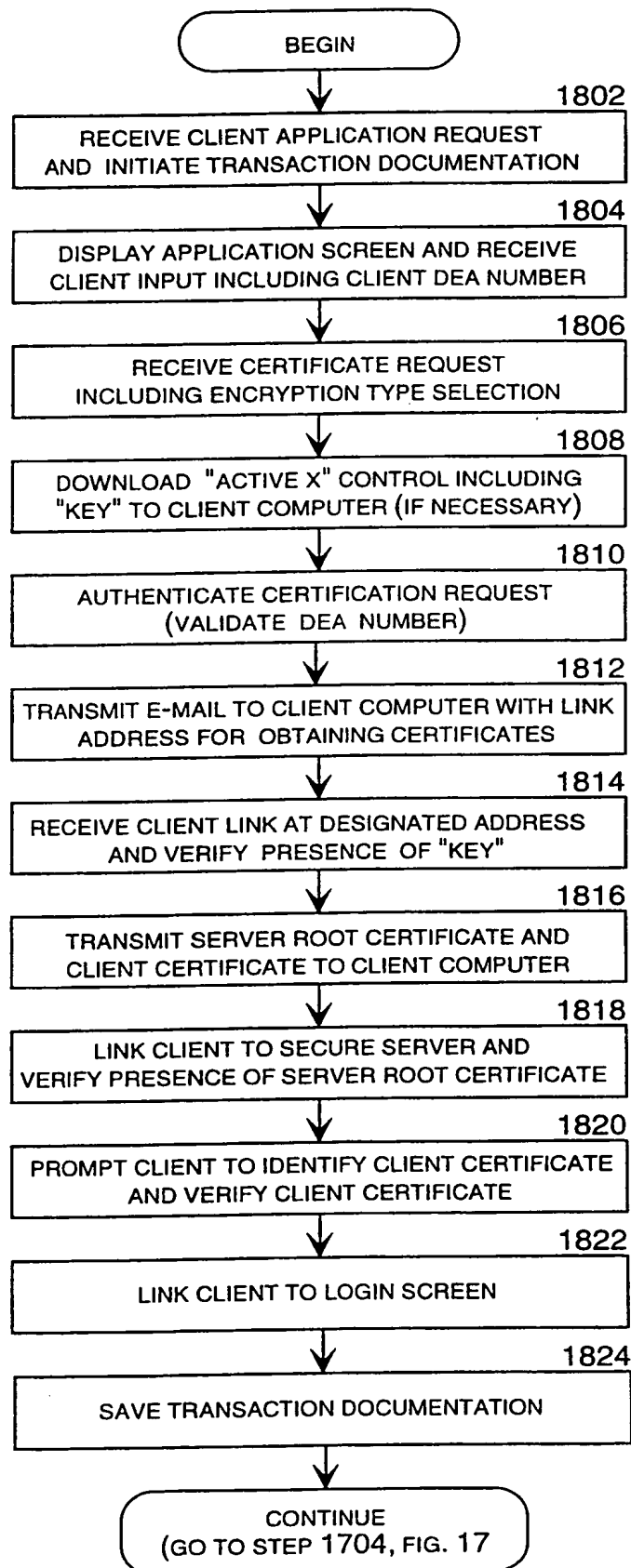CLIENT
APPLICATION
ROUTINE
1702**FIG. 18**

19/30

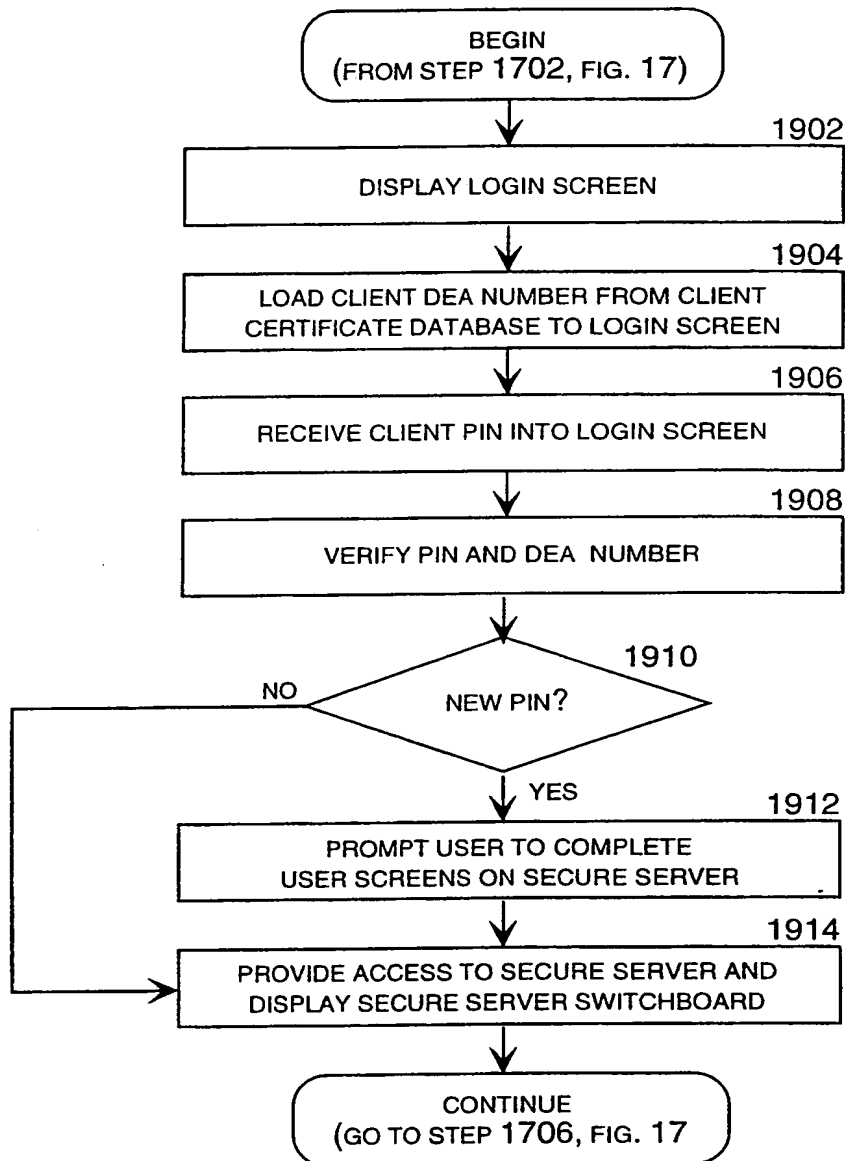CLIENT
LOGIN ROUTINE
1704**FIG. 19**

20/30

INITIAL INPUT SELECTION

2000

2002
Provider

2004
Patients

Programs

Online PDR

Abstracts

FIG. 20

BILLING MODULE

2100

2102

2106

2108

2110

2104

Address

Billing Info

Cover Letter

Account

Provider Number XX1234567

PIN ****

Name

Title	First	Last	Degree
	Lifestream	Technologies	

Address

Street Address 510 Clearwater Loop

Suite 101

City	State/Region	Postal Code
Post Falls	ID	83854

Country United States ▼

Other

Phone (208) 457-9409	Fax (509) 457-9509
Birthdate 1/1/60	E-mail sales@lifestreamtech.com

Associate PIN Numbers

First Name	Last Name
Jane	Worker
Joe	Worker

Add
Modify
Delete

Notes

1.0.3.4

Save

FIG. 21

21/30

2108

2200

2202

Address Billing Info Cover Letter

Billing Method Credit Card

Credit Card Information

Credit Card Number Expiration Date

Name on card

Checking Account Information

Checking Account Routing Number

Save

FIG. 22

22/30

2110

Address Billing Info Cover Letter

From the Office of Lifestream Technologies
510 Clearwater Loop, Suite 101, Post Falls, Idaho 83854
Phone: (208) 457-9409 Fax: (208) 457-9509

~Data~

~Pt.Name~
~Pt.Address~
~Pt.City~, ~Pt.State~ ~Pt.Zip~

Dear ~Pt.FirstName~,



 Editor  Save

FIG. 23


23/30


2400

1.0.253

Select an existing patient or add a new patient


☒ Select Existing Patient





 Read SmartCard


☐ Add New Patient

Administrative

 Backup

 Recreate Card

 Restore

 Quick List


 Next


FIG. 24

24/30

2500

Patient Information

Account		Other	
File Number	Entry Date	Phone	Fax
123	11/24/98	(509) 555-1212	
Name		Birthdate	Gender
First	Last	1/1/50	<input type="radio"/> Male
Jill	Smith		<input checked="" type="radio"/> Female
Address		E-mail	
Street Address			
110 Main			
City	State/Region	Postal Code	
Spokane	WA	99201	
Country	Language		
United States	English		
		<input checked="" type="checkbox"/> Include on Mailing List	
		Notes	

 Save


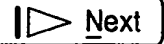
 Previous  Next

FIG. 25

25/30

2600

Questionnaire Data

Patient

File Number	First	Last	Test Time Stamp
123	Jill	Smith	12/8/98 10:16:01AM

Risk Components

Family History	<input checked="" type="checkbox"/>	Height	68 <input type="text"/> inches	Glucose	0 <input type="text"/> mg/dl
Personal History	<input type="checkbox"/>	Weight	165 <input type="text"/> pounds	Body Fat	24 <input type="text"/> %
CVD	<input type="checkbox"/>	Age	48 <input type="text"/> Years	Systolic	140 <input type="text"/>
AF	<input type="checkbox"/>	Chol	215 <input type="text"/> mg/dl	Diastolic	80 <input type="text"/>
LVH	<input type="checkbox"/>	Tng	0 <input type="text"/> mg/dl	Fitness	Sedentary <input type="text"/>
Diabetes Type 1	<input type="checkbox"/>	LDL	0 <input type="text"/> mg/dl	Ethnicity	Caucasian <input type="text"/>
Diabetes Type 2	<input type="checkbox"/>	HDL	0 <input type="text"/> mg/dl	Gender	<input type="radio"/> Male <input checked="" type="radio"/> Female
Smoker	<input checked="" type="checkbox"/>				

Navigation: Record 1 of 1

Navigation:

FIG. 26

26/30


2700

Generate Reports

Current Test

File Number	First	Last	Test Time Stamp
123	Jill	Smith	12/8/98 10:16:01AM

<input checked="" type="checkbox"/> Full Program	Lifestyle Therapy	Lipid Drug	Blood Pressure Drug
<input checked="" type="checkbox"/> Cover Letter	None	None	None
<input checked="" type="checkbox"/> Summary	NCEP		
<input checked="" type="checkbox"/> Evaluation			
<input checked="" type="checkbox"/> Receipt			

 **Generate Report**

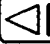
 **Previous**

FIG. 27

27/30

2802

	CORONARY RISK FACTORS	TEST RESULTS	IDEAL RANGE	GOALS
	GENDER	MALE		
IDEAL	PERSONAL HISTORY	NO	NONE	
IDEAL	FAMILY HISTORY	NO	NONE	
	CVD	NO		
	AF	NO		
	LVH	NO		
MODERATE	DIABETES (TYPE 1)	YES	NONE	
HIGH	SMOKER	YES	NO	
	HEIGHT	66 in		
	WEIGHT	155 lbs		
	AGE	44		
MODERATE	TOTAL CHOLESTEROL	211	< 200	
	TRIGLYCERIDES	200		
	HDL	N/A	45-65	
	LDL	N/A	65-135	
	GLUCOSE	N/A		
IDEAL	PERCENTAGE OF BODY FAT	N/A	18%	
IDEAL	BP SYSTOLIC	115	< 120	
LOW	BP DIASTOLIC	80	< 80	
LOW	FITNESS	MODERATE	HIGH	

2804

PERSONAL HEALTH CONSEQUENCES		
BODY MASS INDEX (BMI)	25	< 25
POUNDS OVERWEIGHT	0	0
CHOLESTEROL/HDL RATIO	N/A:1	< 3.5:1
CARDIAC RISK	5 YEARS - 5% 10 YEARS - 10%	
BIOLOGICAL AGE	47	< 44
STROKE RISK	10 YEARS - LOW RISK	

2806

EXTENDED HEALTH ASSESSMENT SUMMARY

CHRONOLOGICAL AGE	30	40	50	60	70
CARDIAC AGE	N/A	N/A	52	62	69

FIG. 28

28/30



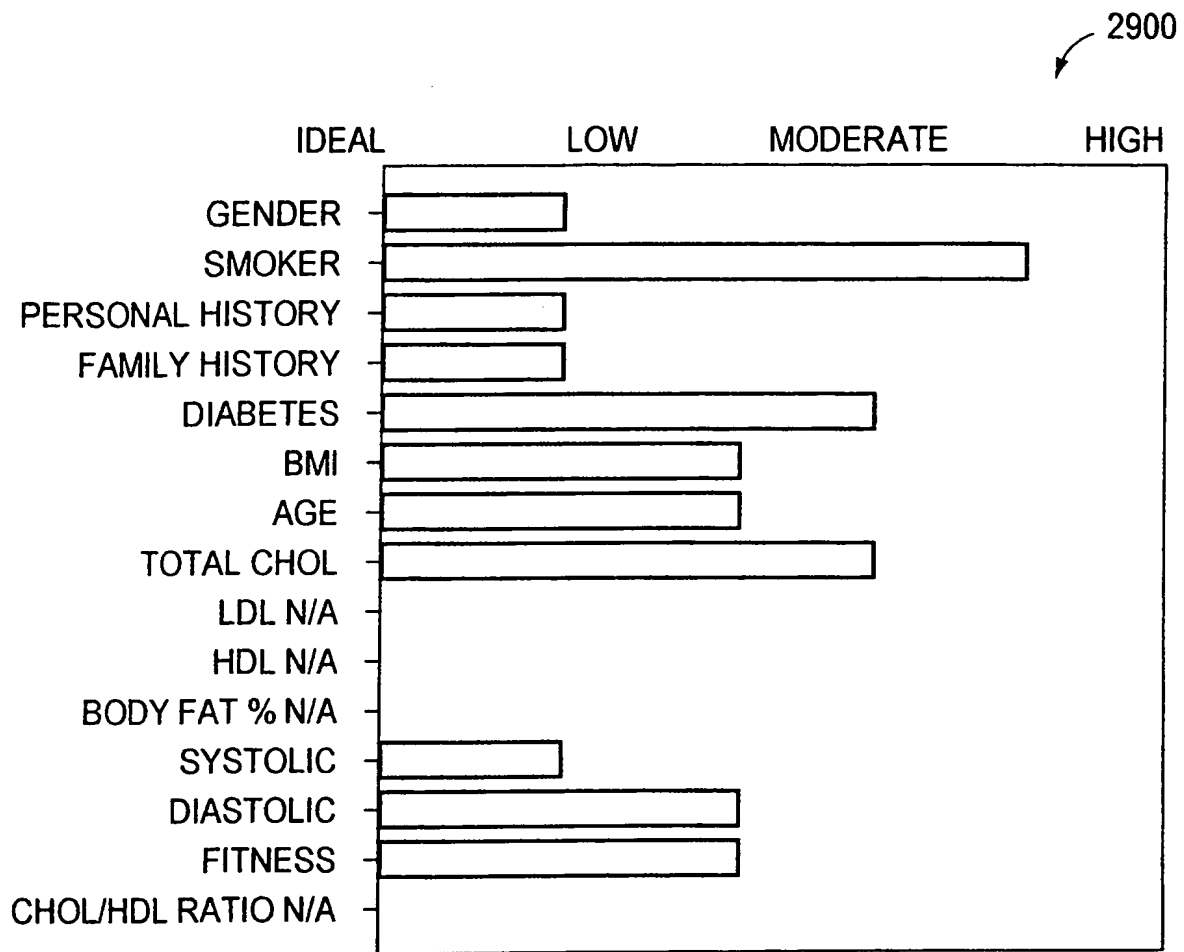


FIG. 29

**FIG. 30**

**FIG. 31**

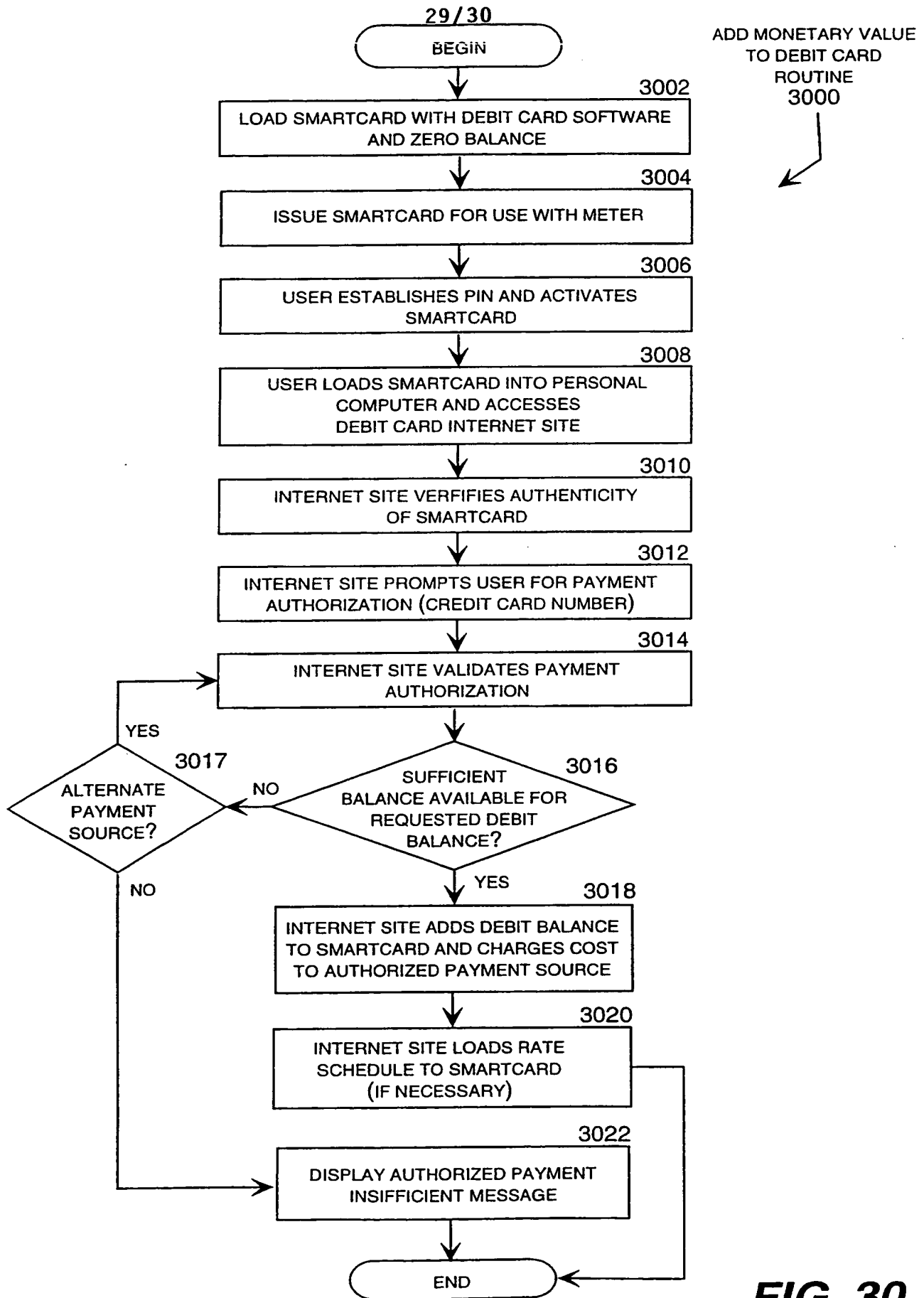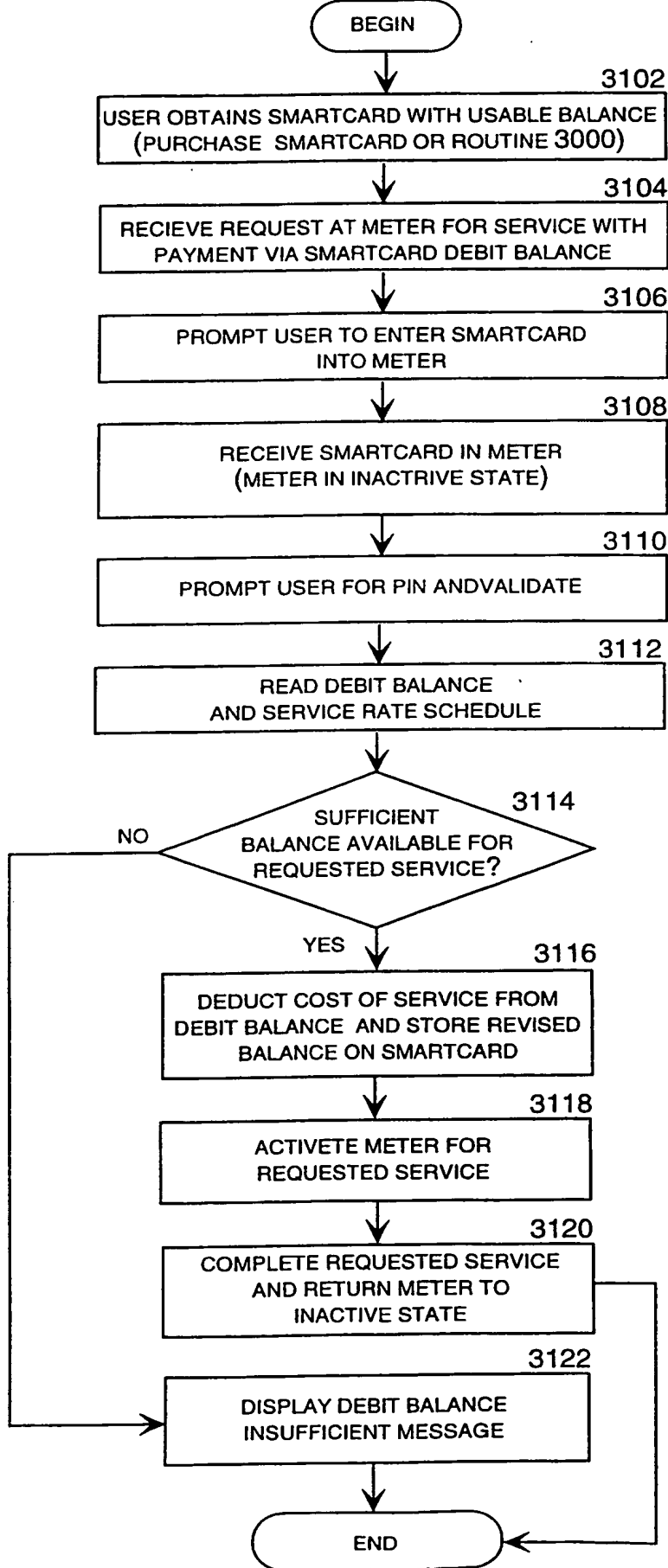

(81) **Designated States (national):** AE, AL, AM, AT, AU, AZ, BA, BB, BG, BR, BY, CA, CH, CN, CR, CU, CZ, DE, DK, DM, EE, ES, FI, GB, GD, GE, GH, GM, HR, HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MA, MD, MG, MK, MN, MW, MX, NO, NZ, PL, PT, RO, RU, SD, SE, SG, SI, SK, SL, TJ, TM, TR, TT, TZ, UA, UG, US, UZ, VN, YU, ZA, ZW.

(84) **Designated States (regional):** ARIPO patent (GH, GM, KE, LS, MW, SD, SL, SZ, TZ, UG, ZW), Eurasian patent (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM), European patent (AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU,

MC, NL, PT, SE), OAPI patent (BF, BJ, CF, CG, CI, CM, GA, GN, GW, ML, MR, NE, SN, TD, TG).

Published:

— With international search report.

(88) **Date of publication of the international search report:**
12 July 2001

For two-letter codes and other abbreviations, refer to the "Guidance Notes on Codes and Abbreviations" appearing at the beginning of each regular issue of the PCT Gazette.

INTERNATIONAL SEARCH REPORT

International Application No

PC, US 99/26521

A. CLASSIFICATION OF SUBJECT MATTER
IPC 7 G06F19/00

According to International Patent Classification (IPC) or to both national classification and IPC

B. FIELDS SEARCHED

Minimum documentation searched (classification system followed by classification symbols)

IPC 7 G06F G01N

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the international search (name of data base and, where practical, search terms used)

WPI Data, PAJ, INSPEC, MEDLINE, COMPENDEX, EPO-Internal

C. DOCUMENTS CONSIDERED TO BE RELEVANT

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
X	US 4 975 647 A (DALKE DAVID M ET AL) 4 December 1990 (1990-12-04)	1-4
Y	column 1, line 30 -column 3, line 10 column 3, line 39 -column 5, line 10 column 5, line 34 -column 7, line 6 ---	17-25
A	US 5 053 199 A (JACKSON CHARLES A ET AL) 1 October 1991 (1991-10-01) column 2, line 2 - line 2 column 2, line 39 -column 3, line 31 ---	1-4
A	US 5 780 304 A (DAFFERN GEORGE MICHAEL ET AL) 14 July 1998 (1998-07-14) column 3, line 65 -column 4, line 60 column 5, line 30 -column 6, line 35 --- -/--	1-4

☒ Further documents are listed in the continuation of box C.☒ Patent family members are listed in annex.

* Special categories of cited documents:

"A" document defining the general state of the art which is not considered to be of particular relevance

"E" earlier document but published on or after the international filing date

"L" document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)

"O" document referring to an oral disclosure, use, exhibition or other means

"P" document published prior to the international filing date but later than the priority date claimed

"T" later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention

"X" document of particular relevance: the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone

"Y" document of particular relevance: the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art.

"&" document member of the same patent family

Date of the actual completion of the international search

20 December 2000

Date of mailing of the international search report

04.01.01

Name and mailing address of the ISA

European Patent Office, P.B. 5818 Patentlaan 2
NL - 2280 HV Rijswijk
Tel. (+31-70) 340-2040, Tx. 31 651 epo nl,
Fax: (+31-70) 340-3016

Authorized officer

Barba, M

INTERNATIONAL SEARCH REPORT

International Application No

PC, JS 99/26521

C.(Continuation) DOCUMENTS CONSIDERED TO BE RELEVANT.

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to Claim No.
A	D.A. ADAMS ET AL: "Point of Care Technology: The i-STAT System for Bedside Blood Analysis" JOURNAL OF PEDIATRIC NURSING, vol. 10, no. 3, June 1995 (1995-06), pages 194-198, XP002139626 USA	1-4
X	page 194, left-hand column, line 1 - line 18 page 194, right-hand column, line 14 -page 195, right-hand column, line 14 page 196, left-hand column, line 33 -right-hand column, line 25 ---	5-16
A	N.L. ALPERT: "i-STAT Point of Care Testing System" CLINICAL INSTRUMENTS SYSTEMS, vol. 13, no. 4, 1994, pages 1-7, XP002139627 USA	1-4
X	page 1, left-hand column, line 12 -page 3, right-hand column, line 22 page 3, right-hand column, line 30 -page 4, right-hand column, line 4 ---	5-16
X	WO 94 19684 A (DIAMETRICS MEDICAL INC) 1 September 1994 (1994-09-01)	5-16
Y	page 1, line 23 -page 2, line 2 page 5, line 24 -page 7, line 32 page 8, line 25 -page 9, line 5 page 10, line 15 -page 13, line 34 page 14, line 7 - line 29 page 16, line 3 - line 33 ---	17-25
Y	US 5 822 715 A (WORTHINGTON DAVID R L ET AL) 13 October 1998 (1998-10-13) column 4, line 1 - line 60 column 5, line 30 -column 6, line 24 column 9, line 21 -column 10, line 22 column 12, line 13 -column 14, line 48 column 17, line 12 -column 18, line 12 ---	26-33
Y	WO 98 24358 A (SANDERS MATTHEW H ;FENSON EITAN (US); ENACT HEALTH MANAGEMENT SYST) 11 June 1998 (1998-06-11) page 3, line 3 - line 5 page 3, line 10 -page 5, line 2 page 8, line 16 -page 9, line 26 page 10, line 7 -page 11, line 20 page 12, line 13 - line 30 page 13, line 29 -page 15, line 26 ---	26-33
	--- -/--	

INTERNATIONAL SEARCH REPORT

International Application No

PC, JS 99/26521

C.(Continuation) DOCUMENTS CONSIDERED TO BE RELEVANT

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
A	GRIMSON J ET AL: "A CORBA-BASED INTEGRATION OF DISTRIBUTED ELECTRONIC HEALTHCARE RECORDS USING THE SYNAPSES APPROACH" IEEE TRANSACTIONS ON INFORMATION TECHNOLOGY IN BIOMEDICINE,US,IEEE SERVICE CENTER, PISCATAWAY, NJ, vol. 2, no. 3, September 1998 (1998-09), pages 124-138, XP000831158 ISSN: 1089-7771 page 131, left-hand column, line 26 -page 134, left-hand column, line 13 ---	34-44
A	WO 96 41288 A (E SYSTEMS INC) 19 December 1996 (1996-12-19) page 4, line 1 -page 7, line 10 page 9, line 18 -page 15, line 19 ---	34-44
A	WO 98 15910 A (SCHULTZ JOSEPH PAUL ;SCHULTZ MYRON G (US)) 16 April 1998 (1998-04-16) page 2, line 28 -page 5, line 5 page 11, line 6 - line 27 page 15, line 1 - line 30 page 16, line 19 -page 23, line 8 page 23, line 28 -page 24, line 3 page 33, line 14 -page 37, line 29 page 44, line 1 -page 46, line 2 ---	34-44
A	US 5 258 906 A (KROLL MARK W ET AL) 2 November 1993 (1993-11-02) column 2, line 1 - line 41 column 3, line 44 -column 5, line 36 column 6, line 1 - line 37 ---	45-49
A	US 4 706 275 A (KAMIL ZVI) 10 November 1987 (1987-11-10) the whole document ---	45-49
A	"LFST Signs Development and Production Agreement with Boehringer Mannheim" LIFESTREAM TECHNOLOGY, INC. INTERNET, [Online] 8 October 1997 (1997-10-08), pages 1-2, XP002156015 Retrieved from the Internet: <URL:http://lifestreamtech.com> [retrieved on 2000-11-20] the whole document --- -/--	5-25

INTERNATIONAL SEARCH REPORT

International Application No

PCT, JS 99/26521

C.(Continuation) DOCUMENTS CONSIDERED TO BE RELEVANT

Category	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
A	"FDA Clears New Hand-Held Instrument to Screen and Monitor Cholesterol" LIFESTREAM TECHNOLOGIES, INC. INTERNET, [Online] 5 October 1998 (1998-10-05), pages 1-2, XP002156016 Retrieved from the Internet: <URL:http://lifestreamtech.com> [retrieved on 2000-11-20] the whole document	5-25
A	POLEMI D: "Trusted third party services for health care in Europe" FUTURE GENERATION COMPUTER SYSTEMS, JUNE 1998, ELSEVIER, NETHERLANDS, vol. 14, no. 1-2, pages 51-59, XP002156017 ISSN: 0167-739X page 52, left-hand column, line 20 -page 53, right-hand column, line 13 page 55, left-hand column, line 23 -page 57, left-hand column, line 26 page 57, right-hand column, line 1 - line 20 -----	34-44

INTERNATIONAL SEARCH REPORT

International application No.
PCT/US 99/26521

Box I Observations where certain claims were found unsearchable (Continuation of item 1 of first sheet)

This International Search Report has not been established in respect of certain claims under Article 17(2)(a) for the following reasons:

1. ☐ Claims Nos.:
because they relate to subject matter not required to be searched by this Authority, namely:
2. ☐ Claims Nos.:
because they relate to parts of the International Application that do not comply with the prescribed requirements to such an extent that no meaningful International Search can be carried out, specifically:
3. ☐ Claims Nos.:
because they are dependent claims and are not drafted in accordance with the second and third sentences of Rule 6.4(a).

Box II Observations where unity of invention is lacking (Continuation of item 2 of first sheet)

This International Searching Authority found multiple inventions in this international application, as follows:

see additional sheet

1. ☒ As all required additional search fees were timely paid by the applicant, this International Search Report covers all searchable claims.
2. ☐ As all searchable claims could be searched without effort justifying an additional fee, this Authority did not invite payment of any additional fee.
3. ☐ As only some of the required additional search fees were timely paid by the applicant, this International Search Report covers only those claims for which fees were paid, specifically claims Nos.:
4. ☐ No required additional search fees were timely paid by the applicant. Consequently, this International Search Report is restricted to the invention first mentioned in the claims; it is covered by claims Nos.:

Remark on Protest

- ☐ The additional search fees were accompanied by the applicant's protest.
- ☒ No protest accompanied the payment of additional search fees.

FURTHER INFORMATION CONTINUED FROM PCT/ISA/ 210

This International Searching Authority found multiple (groups of) inventions in this international application, as follows:

1. Claims: 1-4

Test medium for use with a health monitoring device

2. Claims: 5-25

Health monitoring device

3. Claims: 26-33

System for remotely producing health reports

4. Claims: 34-44

Secure medical records maintenance system

5. Claims: 45-49

System to activate a health monitoring device for performance of a diagnostic clinical test based on a monetary balance associated with a patient

INTERNATIONAL SEARCH REPORT

Information on patent family members

International Application No

PC., JS 99/26521

Patent document cited in search report	Publication date	Patent family member(s)	Publication date
US 4975647 A	04-12-1990	DK 92989 A	02-05-1990
		EP 0366854 A	09-05-1990
		JP 2151763 A	11-06-1990
		PT 89844 A	31-05-1990
US 5053199 A	01-10-1991	AT 152844 T	15-05-1997
		AU 5189690 A	26-09-1990
		CA 2010400 A,C	21-08-1990
		DE 69030657 D	12-06-1997
		DE 69030657 T	11-09-1997
		EP 0460074 A	11-12-1991
		ES 2100881 T	01-07-1997
		JP 2897845 B	31-05-1999
		JP 4506721 T	19-11-1992
		KR 154314 B	01-12-1998
		WO 9010236 A	07-09-1990
US 5780304 A	14-07-1998	AU 709992 B	09-09-1999
		AU 3723595 A	27-03-1996
		CA 2199494 A	14-03-1996
		EP 0779984 A	25-06-1997
		JP 10505676 T	02-06-1998
		NO 971008 A	25-04-1997
		WO 9607908 A	14-03-1996
WO 9419684 A	01-09-1994	AU 6245694 A	14-09-1994
US 5822715 A	13-10-1998	US 5956501 A	21-09-1999
WO 9824358 A	11-06-1998	AU 5461298 A	29-06-1998
		AU 8165498 A	04-01-1999
		WO 9859487 A	30-12-1998
WO 9641288 A	19-12-1996	AU 6252496 A	30-12-1996
WO 9815910 A	16-04-1998	AU 4606597 A	05-05-1998
US 5258906 A	02-11-1993	NONE	
US 4706275 A	10-11-1987	NONE	

THIS PAGE BLANK (USPTO)